



Everyone in care will be kept **comfortable**,
safe and **secure** when protected by Alerta!



Alerta Medical are a world leader in the design, manufacture and supply of innovative and trusted durable medical equipment, delivering solutions to proactively maximise user safety and comfort.

Distributor Commitment

Alerta Medical supply exclusively through a network of distributor partners across the UK and in over 40 countries worldwide. The Alerta brand has become renowned for innovation, quality and reliability, which combined with our industry leading distributor support, provides a unique advantage to our distributors through enhanced brand image and increased profitability.

Delivering Frontline Solutions

Alerta Medical understand care providers need dependable and robust equipment to fulfil the rigorous requirements of care and meet high expectations. Alerta Medical was established through a belief that prevention is better than cure and everyone in care should enjoy the benefits of a proactively safer care environment.

Alerta Medical Mission

Proactively improving the comfort and quality of life for people in care by making cutting-edge medical equipment affordable to all through strategic distributor partnerships. Everyone in care will be kept comfortable, safe and secure when protected by Alerta.

Quality Accreditation

Our medical equipment is manufactured in ISO13485 approved production facilities with thorough quality control procedures, providing distributors with total assurance, strongly backed by our product expertise, delivery capability and customer support, underpinned with a clear vision for the future.

Alerta Medical Vision

To be the partner of choice for companies operating in the healthcare market, by building transparent relationships, unrelenting innovation and the pursuit of excellence.

Proven Experience

An independent and family owned company, Alerta Medical builds on a legacy of three generations, dedicated to improving clinical outcomes and user experience at the frontline of healthcare.

We look forward to serving you!



Certificate Number:
21061-QMS-001
ISO9001



Head Office

Alerta Medical Ltd
4 Symington Place
Riverside Business Park
Irvine KA11 5DE
United Kingdom



FAST PERSONAL SERVICE

Speed is everything. People buy from people. With Alerta you get both; fast-paced personal service is in the veins of the Alerta team. From an email response to a product despatch, we do it right and we do it right now.



EASY TO WORK WITH

As a global designer and manufacturer, we provide in depth product knowledge whenever you need it. Pragmatic fast response is in our DNA.

Quite simply, we are easy to work with.



ALWAYS HAVE STOCK

When you need it, and you need it big, and you need it now; Alerta are always right there for you, we overstock so you don't have to. Our growing distribution centres enable us to hold vast stocks you can always rely on.



UNBEATABLE QUALITY

Clinically proven over many years, our products are designed, tested and enhanced to consistently perform in the healthcare environment **Enhancing your brand image as a distributor. It's called: Reliability**



MANUFACTURER PRICES

Price is also important, as the first step in the supply chain, we offer an attractive pricing structure. **Increasing your profitability as a distributor.**



Contents

Beds & Accessories

Lomond Community Bed	
Lomond Standard Bed	
Lomond Low Bed	
Lomond Floor2 Bed	
Lomond Bariatric Bed	
Lomond Rotate-Stand-Up Bed	
Bed Accessories	
Bed Safety & Bumpers	
Patient Wedges	
Side Safety Wedge Set	
Slide Sheets	
Overbed Tables	

Laundry & Washables

Laundry Carts	
Laundry Bags	
Bed Pads	
Chair Pads	
Bibs	
Pillows & Duvets	
SensaStretch	

Walking Aids

Aluminium Indoor Rollator	
Aluminium Rollator with Solid Back	
Lightweight Walker with Seat & Bag	
Lightweight Walker with Basket, Tray & Bag	
Ultra-Lite Carbon Rollators	
Aluminium Foldable Rollator	
Aluminium 2-in-1 Wheelchair Rollator	
NEW Powered Rollator	
NEW Strato-Cruise Mobility Scooter	

Toileting, Showering, Moving & Handling

MoveAssist 200 Transfer Aid	
Powerlifter Stand-Aid 200	
Powerlifter Hoist Mini 135	
Powerlifter Hoist Midi 150	
Powerlifter Hoist Maxi 175	
Powerlifter Hoist Maxi 175-E	
Standing Turner	
Patient Support Belts	
Wheelchairs	
Aqua Shower Commode Chair	
Aqua Shower Commode Chair, Self-Propelled	
Aqua Shower Commode Chair, Bariatric	
Aqua Chair Tilt-in-Space	
Aqua Tilt-in-Space Hi-Lo Shower Commode Chair	
Showering & Toileting	
Evacuation	

Foam, Gel & Static Systems

Product Overview	
SensaFlex Underlay	

05-14	Foam Mattress SensaFlex 500	42
06	Foam Mattress SensaFlex 1000	43
06	Foam Mattress SensaFlex 3000	43
07	Foam Mattress SensaFlex 4000	44
07	Foam Mattress SensaFlex 4500	45
08	Foam Cushions	46
08	Gel Cushions	47
09	Static Air - Cushion	48
10-11	Static Air - Bootees	48
12	Deluxe Heel Protector	48
13	Heel Protector	48

Fall Prevention

49-61	Product Overview	50
15-19	Deluxe Alertamat	51
16	Wood Alertamat	52
16	Alertamat Crash Mat	52
17	Floor+ Alertamat	53
17	Bed Alertamat	54
17	Chair Alertamat	55
18	Detect Motion PIR Sensor	56
19	Wireless Window/Door Sensor	56
	Wall Point Receiver	57
20-27	Alerta Nurse Call Button	57
21	WetSense Mats	58
22	Accessories & Plug Identification	59
22	SafeHip Hip Protection	60-61

Alternating Air Systems

62-87	Product Overview	63
	Air Mattress Bubble 2 Overlay	64
	Air Mattress Pearl Overlay	65
	Air Mattress Emerald2 Overlay	66-67
	Air Mattress Emerald Auto Overlay	68-69
	Air Mattress Sapphire2 Replacement	70-71
	Air Mattress Ruby2 Replacement	72-73
	Air Mattress Ruby Auto Replacement	74-75
	Air Mattress Ruby Turn Lateral Turning	76-77
	Air Mattress Bariatric2 Replacement	78-79
	SensaFlo Hybrid Mattress	80-81
	SensaFlo Hybrid Cushion	82-83
	Partner Air Cushion	84
	Partner Plus Air Cushion	84
	Mobile Battery Air Cushion	85
	Product Selector Guide	86-87

Suction Machines & Nebulisers

88-89	Suction Machines, Mains & Battery	88
	Nebulisers	89

Diagnostics

90	Thermometers, Non-Contact	90
	Thermometers, Contact	90
	Blood Pressure Monitor	90
	Oximeter	90





Beds & Accessories

Care Beds



20CM EXTENSION AVAILABLE

The **Alerta Lomond Community Bed** provides the ultimate balance between user care, functionality, economy and quality. It has all the functionality of the Lomond Standard including bed rails, and benefits from Trendelenburg and reverse Trendelenburg tilting feature. Alerta Lomond Community bed boasts a 30cm to 70cm Operating Height Ratio providing the perfect harmony of working levels for the carer and user heights for the patient. It can be dismantled for storage or transporting and is available in an Oak wood finish.


ALT-LB-CO LOMOND COMMUNITY BED, OAK

KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 3" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling

TECHNICAL SPECIFICATION

- > Maximum Safe Working Load: 200kg
- > Height Adjustment: 30 - 70cm
- > Outer Dimensions: 102 x 209cm
- > Mattress Platform Dimensions: 90 x 200cm

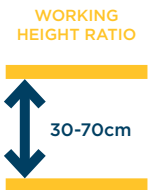

ALSO AVAILABLE IN OAK

20CM EXTENSION AVAILABLE

The **Alerta Lomond Standard Bed** provides the perfect balance between budget, functionality and quality, and comes with the added benefit of Trendelenburg and reverse Trendelenburg tilting feature.

ALT-LB-SW LOMOND STANDARD BED, WALNUT


ALT-LB-SO LOMOND STANDARD BED, OAK

KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 5" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting

TECHNICAL SPECIFICATION

- > Max Safe Working Load: 200kg
- > Height Adjustment: 40 - 80cm
- > Outer Dimensions: 102 x 213cm
- > Mattress Platform Dimensions: 90 x 200cm

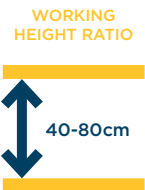

All Lomond Beds have a 3 year warranty on the bed frame and 2 years on the motors and electrics. The beds can be dismantled for storage or transporting and are available in Oak and Walnut wood finish.

alerta
lomond low

ALSO
AVAILABLE
IN WALNUT



20CM EXTENSION AVAILABLE

The **Alerta Lomond Low Bed** provides the ultimate balance between user care, functionality and quality. It has all the features of the Lomond Standard and comes with the added benefit of Trendelenburg and reverse Trendelenburg tilting feature. While being a very low bed, the Lomond Low also has the unique ability of its range extending from 20cm to 80cm high.

ALT-LB-LW LOMOND LOW BED, WALNUT

ALT-LB-LO LOMOND LOW BED, OAK

KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 3" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting

TECHNICAL SPECIFICATION

- > Max Safe Working Load: 250kg
- > Height Adjustment: 20 - 80cm
- > Outer Dimensions: 102 x 213cm
- > Mattress Platform Dimensions: 90 x 200cm

WORKING
HEIGHT RATIO



alerta
lomond floor2

ALSO
AVAILABLE
IN OAK



20CM EXTENSION AVAILABLE

The **Alerta Lomond Floor2 Bed** is a unique ultra-low bed for pioneering user safety and care. It has all the features of the Lomond Low and comes with the added benefit of its range going down to just 9cm from the ground, **and auto regression backrest**. This facilitates the most specific and complex care requirements, while also reducing impact if the user rolls out of bed.

ALT-LB-FW LOMOND FLOOR BED, WALNUT

ALT-LB-FO LOMOND FLOOR BED, OAK

KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Auto regression backrest
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 3" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting

TECHNICAL SPECIFICATION

- > Max Safe Working Load: 200kg
- > Height Adjustment: 9 - 70cm
- > Outer Dimensions: 102 x 213cm
- > Mattress Platform Dimensions: 90 x 200cm

WORKING
HEIGHT RATIO



All Lomond Beds have a 3 year warranty on the bed frame and 2 years on the motors and electrics. The beds can be dismantled for storage or transporting and are available in Oak and Walnut wood finish.



ALSO
AVAILABLE
IN OAK



20CM EXTENSION AVAILABLE

The **Alerta Lomond Bariatric 120cm Wide Bed** provides the ultimate balance between user care, functionality and quality. It has all the features of the Lomond Standard and has the added ability of its range extending from 20cm to 80cm high, and a 120cm wide platform.


- ALT-LB-BW LOMOND BARIATRIC BED, WALNUT
- ALT-LB-BO LOMOND BARIATRIC BED, OAK

KEY FEATURES AND BENEFITS

- > 10 function lockable handset
- > Trendelenburg and reverse Trendelenburg tilting feature
- > Knee brake adjustment
- > Click stop lower leg raising elevation
- > Electrically adjustable backrest
- > Four robust braking 3" castors
- > Four-section mattress support with mattress grips
- > Extra high side rails available
- > 20cm bed length extension kit available
- > Tool-free assembling and dismantling
- > Transport system for storage and transporting

TECHNICAL SPECIFICATION

- > Max Safe Working Load: 250kg
- > Height Adjustment: 20 - 80cm
- > Outer Dimensions: 132 x 213cm
- > Mattress Platform Dimensions: 120 x 200cm

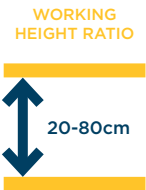

SCAN TO WATCH
PRODUCT VIDEO



The **Alerta Lomond Rotate-Stand-Up Bed** is a powered rotation bed to support users with unassisted transfers between lying in bed and standing, providing enhanced freedom, safety and care for the user.

The user can independently control the bed with the handset to transition from lying in bed to the seated position, and then rotating and standing up from the bed. The hand rails provide additional safety and assurance for the user when rotating and transferring in or out of the bed.

The **Alerta Lomond Rotate Bed** also features a battery back-up ensuring the bed never gets stuck in a sitting position in the event of a power cut. The bed can perform approx. 5 rotations on battery power.


- ALT-LB-RO LOMOND ROTATE-STAND-UP BED, OAK

KEY FEATURES AND BENEFITS

- > 12 function lockable handset
- > Knee brake adjustment
- > Electrically adjustable backrest
- > Four robust braking 4" castors
- > Four-section mattress support with mattress grips
- > Tool-free assembling and dismantling
- > Memory foam articulated mattress included
- > Hand rails to assist onto and off the bed in chair position
- > One button function from lying to sitting
- > Battery back-up in case of mains failure

TECHNICAL SPECIFICATION

- > Max Safe Working Load: 250kg
- > Height Adjustment: 40 - 80cm
- > Outer Dimensions: 102 x 213cm
- > Mattress Width: 82cm

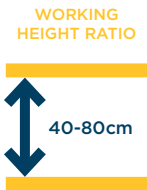

Bed Accessories

ALT-LB-EK90

**LOMOND BED 20CM
EXTENSION KIT, 90CM WIDE**



ALT-66

**ALERTA FOAM MATTRESS 20CM
EXTENSION INSERT, 90 X 14CM**



ALT-LB-LSRW

**LOMOND WALNUT LONG
SIDE RAILS, SET OF 4**



ALT-LB-EK120

**LOMOND BED 20CM
EXTENSION KIT, 120CM WIDE**

ALT-66B

**ALERTA FOAM MATTRESS 20CM
EXTENSION INSERT, 120 X 16CM**

ALT-LB-LSRO

**LOMOND OAK LONG
SIDE RAILS, SET OF 4**

ALT-LB-HSRW

**LOMOND WALNUT
SINGLE HIGH SIDE RAIL**

ALT-LB-HSRO

**LOMOND OAK
SINGLE HIGH SIDE RAIL**

ALT-LB-WP

**LOMOND BED
WALL PROTECTOR**



ALT-LB-AL

**LOMOND BED
ASSISTING LEVER**

ALT-LB-LP

LOMOND BED LIFTING POLE
> Comes with strap and handle

ALT-LB-MT

**LOMOND BED
MOVING TROLLEY**



**WATCH
MOVING TROLLEY
PRODUCT VIDEO**

Bed Safety & Bumpers

Alerta Medical bed protection is made with PU (Polyurethane Fabric) instead of PVC due to risks of Phthalates with PVC. PU is safer and higher quality, stretchable, breathable and serves a longer lifetime.



ALT-112 2 Bar Rail Bumpers with Netting



The **Alerta Bumpers with Net Inserts** are designed to zip onto side rails, providing protection for users in bed.

Net inserts are centred between the bed rails to improve visibility and comfort for the user.

Suitable for profiling beds with 2 side rails.

Also available in 20cm extended - code: ALT-115

ALT-113 3 Bar Rail Bumpers with Netting



The **Alerta Bumpers with Net Inserts** are designed to zip onto side rails, providing protection for users in bed.

Net inserts are centred between the bed rails to improve visibility and comfort for the user.

Suitable for profiling beds with 2 side rails and an additional 3rd high side rail attached.

ALT-114 Full Length Bed Rail Bumpers



The **Alerta Full Length Bed Rail Bumpers** are designed to velcro onto standard height wooden side rails, providing protection for users in bed.

Suitable for profiling beds with 2 side rails.

BREATHABLE | IMPERMEABLE | STRETCH | ANTIBACTERIAL | FR

ALT-118 2 Bar Rail Anti-Entrapment Bumpers with Netting



Alerta Anti-Entrapment Bumpers benefit from being connected with a mesh base which sits on the mattress platform of the bed underneath the mattress, forming a safety net to help prevent limb entrapment.


The **Alerta Anti-Entrapment Bumpers with Net Inserts** are designed to zip onto side rails, providing protection for users in bed. Net inserts are centred between the bed rails to improve visibility and comfort for the user.

Suitable for profiling beds with 2 side rails.



ALT-103 - Alerta Folding Crash Mat is designed to be placed next to the bed to minimise injury if the user should fall from bed. The mat incorporates triple fold design for easy storage and carry handles for ease of transport.

Alerta Folding Crash Mat is quick and easy to position and benefits from a PU wipe clean surface. When unfolded the mat is 190cm long x 60cm wide x 6cm high.

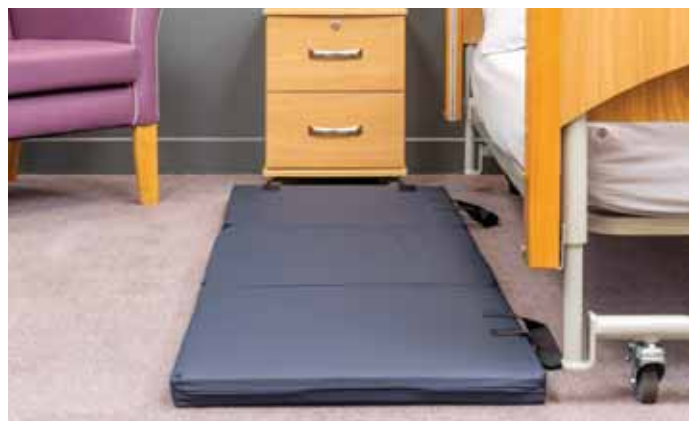

ALT-117 - Alerta Head & Foot Board Protector is designed to protect those in care from hard head and foot boards on their care bed. Made from a wipe clean fabric, each bumper is tapered to allow for tracking side rails and is easy to fit. Use in conjunction with side bumpers to create a comprehensive bed safety system.

Not compatible with Alerta Lomond Floor2 Bed.

KEY FEATURES & BENEFITS

- > Safety padding
- > 105cm wide
- > 58cm high
- > Wipe clean fabric
- > Tapered for neat side rail use
- > Use with Alerta side bumpers for all round protection

BREATHABLE | IMPERMEABLE | STRETCH | ANTIBACTERIAL | FR

Patient Wedges



ALT-PW-M1
ALERTA MEMORY FOAM
LIMB SUPPORT CUSHION

Limb support cushion with open-celled memory foam inner, used for ankle and knee spacing and positioning providing gentle cushioning and creating space between the legs.

25 x 19 x 13-15cm



ALT-PW-M2
ALERTA MEMORY FOAM
MULTI SUPPORT CUSHION

Multi-purpose support cushion with open-celled memory foam inner, used for all needed positioning areas, particularly the back.

49 x 32 x 8cm



ALT-PW-M3
ALERTA MEMORY FOAM
UNIVERSAL SUPPORT CUSHION

Universal support cushion with open-celled memory foam inner, used for all needed positioning areas, particularly knees, ankles, arms, neck and lower back.

49 x 19.5 x 10cm



ALT-PW-S
ALERTA PATIENT FOAM
WEDGE, SMALL (SLOPED)

Patient positioning wedge, small. Normally for use in prone position, to facilitate head, shoulder and arm movement, but also for use for support or elevation of the lower leg, or other positioning as required.


30 x 30 x 7.5-2.5cm



ALT-PW-M
ALERTA PATIENT FOAM
WEDGE, MEDIUM (SLOPED)

Patient positioning wedge, medium. Normally for use in prone position, to facilitate head, shoulder and arm movement, but also for use for support or elevation of the lower leg, or other positioning as required.

60 x 60 x 17.5-5cm

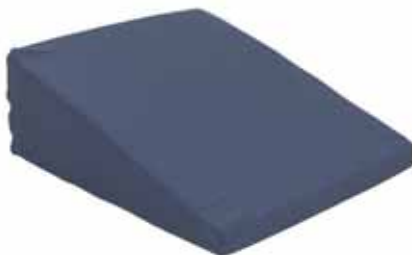

ALT-PW-L
ALERTA PATIENT FOAM
WEDGE, LARGE (SLOPED)

Patient positioning wedge, large. Normally for use in prone position, to facilitate head, shoulder and arm movement, but also for use for support or elevation of the lower leg, or other positioning as required.

60 x 60 x 25-5cm

ALT-116 - Alerta Safety Side Wedge Set is designed to protect those in care and help keep them comfortably positioned in bed. Made from high density foam with a wipe clean PU cover, it has a draw sheet and two built-in straps that secure the wedge set to the bed. 2 sets can be used down the length of the bed for extra protection if required.



KEY FEATURES & BENEFITS

- > Safety wedge
- > Draw sheet included
- > Wipe clean fabric
- > Straps for securing to the bed
- > 99cm wide x 94cm long x 20cm high



Packaged for retail

ALT-SS - Alerta Slide Sheets boast superior slide performance properties and are designed to facilitate nurses and carers with safe assisted repositioning and transfer of patients on a support surface.

Alerta Slide Sheets are tubular, reusable, and available in 5 convenient sizes. They are manufactured from a very strong highly durable quality fabric which creates an extremely low friction surface providing superior sliding performance during patient repositioning and transfers whilst protecting against skin shear. Alerta Slide Sheets are not intended for use to lift patients.

Fitted with a brightly coloured hanging loop allowing the sheet to be conveniently stored near the user.

KEY FEATURES AND BENEFITS

- > Brightly coloured hanging loop
- > Superior slide performance
- > Wash at a maximum of 80°C
- > Do NOT use fabric softener or chlorine
- > Slip hazard; Do not leave on the floor
- > Tubular and reusable

PRODUCT CODES

ALT-SS-S	SMALL 72 X 70CM
ALT-SS-M	MEDIUM 122 X 70CM
ALT-SS-L	LARGE 122 X 100CM
ALT-SS-W	WIDE 145 X 100CM
ALT-SS-XL	X-LARGE 200 X 140CM

First measurement is open side

Overbed Tables



Walnut without castors Oak with castors

The **Alerta Overbed Table** provides easy access to everyday conveniences for users seated in a bed, chair or wheelchair.

It is finished in oak or walnut and comes with or without castors. The table has a laminate top with two adjusting points which allows both height and tilt adjustment and comes with raised edges to prevent items slipping off.


TECHNICAL SPECIFICATION

- > Max weight capacity: 10kg
- > Height adjustment range: 64-97cm
- > Table top dimensions: 61 x 40cm
- > Max Tilt angle: 34°

PRODUCT CODES

ALT-751O	ALERTA OAK OVERBED TABLE, NO CASTORS
ALT-751W	ALERTA WALNUT OVERBED TABLE, NO CASTORS
ALT-751OC	ALERTA OAK OVERBED TABLE, WITH CASTORS
ALT-751WC	ALERTA WALNUT OVERBED TABLE, WITH CASTORS

KEY FEATURES AND BENEFITS

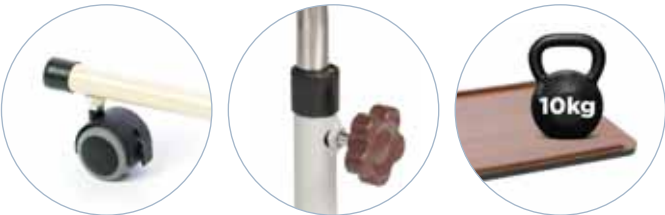

Castor option Height & tilt adjustment Won't topple over even with 10kg on top!



Walnut Oak

The **Alerta Gas Lift Overbed Table** features an easy lift handle to adjust the table height with minimal effort. It is finished in oak or walnut and comes with brakes on the two rear castors. The table has a laminate top and comes with either a straight leg base or a wheelchair accessible base design.

TECHNICAL SPECIFICATION

- > Straight Leg or Wheelchair access base
- > Braking castors
- > Max weight capacity: 15kg
- > Height adjustment range: 80-105cm
- > Table top dimensions: 80 x 40cm
- > Wheelchair Leg Base internal leg width clearance: 64 x 38cm
- > Table base height including castors: 10.5cm

PRODUCT CODES

ALT-GLS-O	ALERTA OAK GAS LIFT OVERBED TABLE, STRAIGHT LEG BASE
ALT-GLS-W	ALERTA WALNUT GAS LIFT OVERBED TABLE, STRAIGHT LEG BASE
ALT-GLW-O	ALERTA OAK GAS LIFT OVERBED TABLE, WHEELCHAIR LEG BASE
ALT-GLW-W	ALERTA WALNUT GAS LIFT OVERBED TABLE, WHEELCHAIR LEG BASE

KEY FEATURES AND BENEFITS



Rear castors Easy lift height adjustment Won't topple over even with 15kg on top!



Laundry
& Washables

Laundry Carts & Bags



The **Alerta Laundry Cart** range is specifically designed to simplify the life of the housekeeper in Care Home and Hospital environments. Crafted from robust steel construction, the carts glide easily over Care Home and Hospital floor surfaces on the easy-glide braked castors, further enhanced with plastic bumper protection.

The Single, Double, Triple and Quadruple Bag Cart design optimises, speeds up and streamlines the collection and segregation of soiled linen on an industrial scale, delivering labour efficiencies and improving hygiene standards.

Alerta Laundry Cart is supplied with colour coded lids for ease of sorting. Bags sold separately.

Carts can be ordered as single, double, triple or quadruple.

- | | |
|------------------|---|
| ALT-LYC/1 | ALERTA LAUNDRY CART, SINGLE BAG
(white lid supplied) |
| ALT-LYC/2 | ALERTA LAUNDRY CART, DOUBLE BAG
(white and red lids supplied) |
| ALT-LYC/3 | ALERTA LAUNDRY CART, TRIPLE BAG
(white, red and blue lids supplied) |
| ALT-LYC/4 | ALERTA LAUNDRY CART, QUADRUPLE BAG
(white, red, blue and green lids supplied) |
| ALT-LCL/B | ALERTA LAUNDRY CART SPARE LID BLUE |
| ALT-LCL/G | ALERTA LAUNDRY CART SPARE LID GREEN |
| ALT-LCL/R | ALERTA LAUNDRY CART SPARE LID RED |
| ALT-LCL/W | ALERTA LAUNDRY CART SPARE LID WHITE |

ALERTA LINEN BAG

- > Tough polyester construction
- > Heavy duty draw cord with toggle
- > Machine washable at 60-70°C
- > Machine dryable (low heat only)
- > Available in 5 colours
- > Linen Bag Size: 85x45cm

- | | |
|------------------|---------------|
| ALT-LYB/W | WHITE |
| ALT-LYB/B | BLUE |
| ALT-LYB/Y | YELLOW |
| ALT-LYB/G | GREEN |
| ALT-LYB/R | RED |



Washables

ALERTA WASHABLE BED PAD

- > 300 wash guarantee
- > Tough durable construction
- > Fully waterproof (integral barrier)
- > Generously sized tuck-in flaps
- > Machine washable
- > Machine dryable (low heat only)

DBP/2	2.1 LITRE ABSORBENCY, 72X90CM
DBP/3	3.2 LITRE ABSORBENCY, 86X90CM
DBP/D	4 LITRE ABSORBENCY, 86X140CM, DOUBLE




Packaged for retail

ALERTA WASHABLE CHAIR PAD

- > 300 wash guarantee
- > Tough, durable construction
- > Fully waterproof (integral barrier)
- > Generously sized (50cm x 60cm)
- > Machine washable
- > Machine dryable (low heat only)
- > 1.0 litre absorbency capacity
- > Available in 5 colours

DSP/BR	BROWN	DSP/PB	PETROL BLUE
DSP/GR	DARK GREEN	DSP/TA	TARTAN
DSP/LC	WINE		

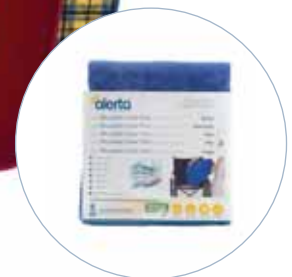

Packaged for retail

ALERTA BIB

- > 300 wash guarantee
- > Tough durable construction
- > Fully waterproof (integral barrier)
- > Machine washable
- > Machine dryable (low heat only)
- > Generous 45cm width
- > Length - 90cm
- > 2 stud fastening positions at neck




Packaged for retail



ALT-BIB-PS
TARTAN PINK
- 90CM



ALT-BIB-BL
TARTAN BLUE
- 90CM

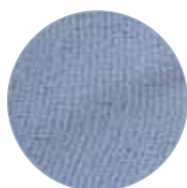

ALT-BIB-PB
PLAIN BLUE
- 90CM



ALT-BIB-PG
PLAIN GREEN
- 90CM



ALT-BIB-PW
PLAIN WINE
- 90CM

Pillows & Duvets



Flame Retardant
BS 7175



100% Quilted
Cotton



100% Polyester
Ball Fibre



Breathable



Washable

- ALT-LLP LONG-LIFE FR LUXURY WASHABLE PILLOW
- ALT-LLDS LONG-LIFE FR LUXURY WASHABLE DUVET - SINGLE
- ALT-LLDD LONG-LIFE FR LUXURY WASHABLE DUVET - DOUBLE

Alerta’s Long-Life Luxury Washable Pillow and Duvet feature a **100% quilted cotton outer**, offering a soft, smooth touch against the skin while allowing for natural breathability. The **100% polyester ball fibre filling** provides gentle yet supportive cushioning, offering excellent insulation and support.

In addition to their luxurious feel, they are both **flame retardant**, meeting rigorous safety standards for your peace of mind. The **washable design** ensures the pillow & duvet stay fresh and clean, even after multiple washes, maintaining their shape and support for the long term.

The 10.5 tog rating on the duvet provides the perfect balance of warmth and breathability, making it ideal for year-round use.

KEY FEATURES & BENEFITS

- > 100% quilted cotton outer
- > 100% polyester ball fibre filling
- > Washable
- > Breathable
- > Duvet is available in 2 different sizes
- > Flame retardant for added safety:
Meets safety standards BS 7175

TECHNICAL SPECIFICATION

- > Pillow size: 50 x 70cm
- > Duvet size (single): 140 x 200cm
- > Duvet size (double): 200 x 200cm

Available individually & bulk boxed

- > ALT-LLP: Boxed in 14’s
- > ALT-LLDS: Boxed in 12’s
- > ALT-LLDD: Boxed in 8’s



SensaStretch Soft-Knit Bed Sheets



ALT-SKS - SensaStretch Soft-Knit Bed Sheets are designed primarily for use with air alternating and foam pressure relief mattresses. The **SensaStretch** sheet is very different to a conventional bedsheet, it fits easily over all air alternating and foam mattress systems, both overlay and deepcell, and is a distinguished, soft, multi-way stretch, **Superior Soft-Knit** sheet.

Alerta SensaStretch Soft-Knit multi-way stretch sheet uniquely preserves and enhances the benefits of pressure relief afforded by alternating air and foam pressure relieving mattress systems, over against standard sheets which do not stretch and consequently reduce the pressure relief, thereby defeating the purpose.

KEY FEATURES & BENEFITS

SensaStretch Soft-Knit

- > Multi-way stretch, no elastic
- > Superior Soft-Knit Luxury
- > White & Cream
- > 90cm Width

User Benefits

- > Multi-way stretch
- > Improved pressure relief
- > Anti-wrinkle design
- > Superior Soft-Knit comfort

Caregiver Benefits

- > 1 person fit
- > Faster fitting time
- > 1 sheet fits all mattresses
- > No ironing required

Laundry Benefits

- > No ironing required
- > Less laundry time



Preserves & enhances effectiveness
of pressure relief mattresses

ALT-SKS/C

**SENSASTRETCH
SOFT-KNIT BED
SHEETS, CREAM**

ALT-SKS/W

**SENSASTRETCH
SOFT-KNIT BED
SHEETS, WHITE**



**Multi-Way Stretch
SensaStretch Soft-Knit**



Walking
Aids

Walking Aids



ALT-R010/B Blue



ALT-R010/G Green



ALT-R010/R Red

The **Alerta Four-Wheel Aluminium Indoor Rollator** is an excellent personal home-help friend, equipped with a useful tray, bag, pull up break & push down hand-brake lock, easy swivel wheels, and features a robust, foldable frame and adjustable handle height. The lightweight narrow construction design offers easy transportation, and the sturdy aluminium frame gives stability and support when in use.

TECHNICAL SPECIFICATION

- > Handle height: 89-93.5cm
- > Max user weight: 120kg
- > Product weight: 5.24kg
- > Product dimensions: 57w x 65l
- > Tray size: 38 x 29cm
- > Basket size: 37w x 27l x 18d cm
- > Handle width (widest point): 41cm
- > Handle width (narrowest point): 28cm
- > Wheel size: 15cm
- > Space between front wheels: 44.5cm
- > Space between rear wheels: 48.5

KEY FEATURES & BENEFITS

- > Available in three colours (blue, red, green)
- > Robust lightweight frame
- > Easy-fold for storage and transport
- > Handle lock safety feature
- > Storage basket
- > Contoured transport tray
- > Adjustable height handles
- > Robust 6" wheels for indoor use





ALT-R002

FOUR-WHEEL ALUMINIUM ROLLATOR WITH SOLID BACK

The **Alerta Four-Wheel Aluminium Rollator with Solid Back** is a feature rich walking aid to assist users both around the house and when out and about. It boasts a strong lightweight aluminium frame equipped with a padded seat for that essential rest, a bag and sturdy backrest, puncture proof tyres and robust build for the daily routines.

TECHNICAL SPECIFICATION

- > Handle height: 82 - 91cm
- > Seat height: 52cm
- > Max user weight: 136kg
- > Product weight: 6.5kg

KEY FEATURES & BENEFITS

- > Seat with sturdy backrest for users to rest
- > Robust lightweight aluminium frame
- > Quick fold design for easy transport and manoeuvrability
- > Handle lock safety feature
- > Spacious storage bag
- > Robust wheels for indoor and outdoor use
- > Detachable and foldable backrest
- > Hinged seat
- > Adjustable height handles



ALT-R004

THREE-WHEEL ALUMINIUM WALKER WITH SEAT & BAG

The **Alerta Three-Wheel Aluminium Walker** is equipped with a seat, backrest and storage pockets, and features a triangular foldable frame and adjustable handle height. It has robust 8" wheels and a lightweight frame for easy transport and manoeuvrability.

TECHNICAL SPECIFICATION

- > Handle height: 82-99cm
- > Seat height: 50cm
- > Max user weight: 135kg
- > Product weight: 6.25kg

KEY FEATURES & BENEFITS

- > Seat with backrest for users to rest
- > Robust aluminium frame
- > Aluminium and foldable for manoeuvrability and transport
- > Handle lock safety feature
- > Storage pockets
- > Robust 8" wheels for indoor and outdoor use



ALT-R014

THREE-WHEEL LIGHTWEIGHT WALKER WITH BASKET, TRAY & BAG

The **Alerta Three-Wheel Lightweight Walker** is equipped with a basket, tray and bag, and features a triangular foldable frame and adjustable handle height. It has robust 8" wheels and a lightweight frame for easy transport and manoeuvrability.

TECHNICAL SPECIFICATION

- > Handle height: 82-97cm
- > Max user weight: 115kg
- > Product weight: 6.5kg

KEY FEATURES & BENEFITS

- > Robust lightweight frame
- > Lightweight and foldable for manoeuvrability and transport
- > Handle lock safety feature
- > Spacious basket and storage bag
- > Includes tray with raised edges
- > Robust 8" wheels for indoor and outdoor use

Walking Aids



ALT-R016 Satin Black



ALT-R018 Satin Red – Low Seat

The **Alerta Four-Wheel Carbon Ultra-Lite Rollator** is one of the world's lightest rollators. It is equipped with a seat, backrest and bag, and features a robust foldable frame and adjustable handle height. It has 4 robust wheels for added stability and an ultra-lightweight carbon fibre frame for very easy transport and manoeuvrability.

TECHNICAL SPECIFICATION

- > Max patient weight: 200kg
- > Max bag load: 5kg
- > Seat height:
 - > ALT-R016: 64cm
 - > ALT-R018: 56cm
- > Seat width: 46cm
- > External folded width: 22cm
- > External opened width: 62cm
- > Total height folded: 82cm
- > Handle height:
 - > ALT-R016: 78-98cm
 - > ALT-R018: 70-90cm
- > Wheel diameter: 20cm
- > Product weight: 5.3kg

KEY FEATURES & BENEFITS

- > Ergonomic adjustable height handles
- > Ultra-lightweight robust carbon frame
- > Available in two colours
- > Handle lock safety feature
- > Spacious storage bag
- > Robust wheels for indoor and outdoor use
- > Easy quick-fold design
- > Reflectors on the front & back
- > Foot tilt mechanism

Walking Aids



SCAN TO WATCH
PRODUCT VIDEO



ALT-R006

FOUR-WHEEL ALUMINIUM DUAL-FOLDING ROLLATOR

The **Alerta Four-Wheel Aluminium Dual-Folding Rollator** is equipped with a seat, bag, and independently adjustable backrest. It boasts a unique patented robust quick-fold lightweight frame for easy transport and manoeuvrability, and also benefits from quick release adjustable height handles.

TECHNICAL SPECIFICATION

- > Handle height: 79-96cm
- > Seat height: 48cm
- > Max user weight: 136kg
- > Product weight: 7.35kg

KEY FEATURES & BENEFITS

- > Seat with independently adjustable backrest for users to rest
- > Robust lightweight aluminium frame
- > Unique patented quick-fold design for easy transport and manoeuvrability
- > Handle lock safety feature
- > Spacious storage bag
- > Robust wheels for indoor and outdoor use

PATENTED FOLD DESIGN
EASIEST TO FOLD ROLLATOR IN THE WORLD



ALT-R008

ALERTA FOUR-WHEEL ALUMINIUM 2-IN-1 WHEELCHAIR ROLLATOR

The **Alerta Four-Wheel Aluminium 2-in-1 Wheelchair Rollator** is equipped with a soft extra-comfort padded seat, sturdy comfort backrest, quick release adjustable height ergonomic shaped handles, pull up break & push down hand-brake lock, swing-away detachable footrests, fixed padded armrests for wheelchair position, and easy one-click pull rapid fold for storage and transport.

TECHNICAL SPECIFICATION

- > Seat height: 55cm
- > Handle height: 86-100cm (6 positions)
- > Max user weight: 136kg
- > Product weight: 9.9kg

KEY FEATURES & BENEFITS

- > Adjustable height handles
- > Adjustable back wheels
- > Spacious storage bag
- > Easy quick-fold design
- > Convert rollator to wheelchair
- > Walking stick holder
- > Seat with adjustable back rest
- > Detachable & adjustable footrests
- > Handle lock safety feature
- > Robust wheels for indoor and outdoor use

Powered Rollator



NEW



ALT-PR001/O Alerta Powered Rollator - Orange

ALT-PR001/N Alerta Powered Rollator - Navy

Crafted from lightweight aluminium, the **Alerta Powered Rollator** is a versatile **4-in-1 mobility solution** designed to enhance freedom, confidence, independence, and overall wellbeing helping you rediscover your natural get up and go.

Walk, sit, assist, or be assisted—the choice is always yours.

This futuristic powered rollator combines advanced **SmartAssist™** technology with everyday comfort features, including dual walking and seating controls, orthopaedically angled footrests, biometrically sensitive push handles, lux-padded ergonomic armrests, an integrated LED light, and a generous storage bag for daily essentials and shopping. All of this is housed within a streamlined, easy-fold design.

The 4-in-1 design allows the Powered Rollator to seamlessly operate as a Rollator, Powered **SmartAssist™** Rehabilitation Walker, Powered **SmartAssist™** Rehabilitation Wheelchair and Electric Wheelchair.

KEY FEATURES & BENEFITS

- > Available in 2 colours; Orange & Navy
- > Reversible backrest
- > 4-in-1 wheelchair rollator
- > 6 speed settings on front controller
- > 6 speed settings on rear controller
- > Lap belt
- > Comes with cupholder
- > Reflectors on the handles
- > Adjustable height handles
- > Handle lock safety feature
- > Spacious storage bag
- > Comfort wheels for indoor and outdoor use
- > Removeable footrests
- > Ergonomic lux-padded armrests
- > Footrest store
- > LED light – comes with USB & USB-C charging ports

TECHNICAL SPECIFICATION

- > Max user weight: 136kg
- > Product weight: 19.8kg
- > Max speed: 3.5mph
- > 36V 7.8Ah lithium battery
- > 10" front wheels
- > 8" rear wheels
- > Adjustable height handles: 86 – 106cm
- > Dimensions (unfolded): 100 x 66 x 89-110cm
- > Dimensions (folded): 71 x 34 x 89-110cm
- > Seat height: 54cm
- > Seat size: 45 x 36cm



ALT-SR-CG Alerta Strato-Cruise Mobility Scooter – Grey

ALT-SR-CR Alerta Strato-Cruise Mobility Scooter – Red

The **Alerta Strato-Cruise** is the tomorrow of mobility scooters, boasting many extra safety features; providing you with 25 miles of cruising freedom with its long lasting 37V battery, giving you cruises in the park, trips to coffee with friends, and cruise all your favourite trips around town with the confidence that you'll cruise back home effortlessly and safely with its LED headlight.

Strato-Cruise incorporates an innovative 2-step quick release folding mechanism making it convenient to travel with, in the boot of your car and easy to fold for storage.

The **Strato-Cruise** effortlessly handles common obstacles like driveway dropped kerbs, door thresholds, uneven pavements, and even the grass in your local park thanks to its powerful 250W brushless motor and large wheels (9" front & 10" rear). Rediscover the joy of a stable, comfortable cruise on almost any surface.

Strato-Cruise acts as your guardian, with auto-deceleration on hills and corners and a smart seat sensor that ensures you're safely seated before cruising, plus voice alerts on steep slopes giving you extra confidence and making every journey secure and serene. Get going from the moment you sit down, the intuitive 'press-to-go' lever and large, clear display puts you in effortless command - tailor your ride with variable speed modes, utilize the reverse gear for tight spots, and enjoy the flexibility of a manual push mode.

The fully adjustable and liftable armrests provide perfect support and easy access, ensuring you start and end every cruise in complete ease.



SCAN TO WATCH
PRODUCT VIDEO



KEY FEATURES



Built-in LED headlight



Display screen



3 speed settings



Anti-roll wheels



Adjustable armrests



Large basket



Voice prompts



Reflectors on back and sides



Bumper



Pneumatic Shock Absorbing Wheels

KEY FEATURES & BENEFITS

- > Available in two colours (grey & red)
- > Detachable and rechargeable battery
- > Drive system: "PAi" Smart Safety Riding System
- > Auto slows downhill & on turns
- > Dual brake system: Electromagnetic brake & electronic brake
- > Gyroscope
- > Basket with handlebar for storage
- > Hook under the handlebar for storage
- > Dual folding
- > Quick release folding mechanism
- > Built-in LED headlight
- > Reflectors on back and sides
- > Fold back armrests
- > Horn
- > 3 speed settings
- > Reverse function
- > Walk function (allows you to walk alongside it without being on it)

TECHNICAL SPECIFICATION

- > Max speed: 6mph (10km/h)
- > Range per charge: 25m (40km)
- > 37V 12Ah lithium battery
- > Motor spec: 250W brushless motor
- > Max climbing angle: 8°
- > Scooter weight: 28.1kg (exc battery)
- > Max user weight: 140kg
- > Max obstacle crossing height: 5cm (2")
- > Front wheels: 23cm (9")
- > Rear wheels: 25.4cm (10")
- > Ground clearance: 7cm (3")
- > Seat size: 40.6 x 35.6cm (16 x 14")
- > Seat height (from ground): 56cm (22")
- > Seat height (from base): 40.5cm (16")
- > Backrest height: 31cm (12.2")
- > Overall length: 90cm (35.4")
- > Overall width: 55cm (21.6")
- > Overall height: 86.5cm (34")
- > Overall height folded (w/o basket): 39.5cm (15.5")
- > Overall height folded (with basket): 61cm (24")
- > Basket size (LxWxH): 23 x 33 x 22cm (9 x 13 x 8.5")
- > Handlebar width: 40cm (15.7")
- > Back of seat to handlebar: 56cm (22")
- > Handlebar height: 69cm (27")
- > Armrest length: 24.5cm (9.6")



Toileting, Showering,
Moving & Handling



ALERTA MOVEASSIST 200

ALT-ST200 - Alerta MoveAssist - a transfer assist unit which makes transfer functions of all types quick and easy, requiring minimal caregiver assistance. **Alerta MoveAssist** promotes user engagement, enhancing mobility.

Each unit is equipped with a crossbar where users can grasp and pull their self up into a standing position using their own strength. The padded split seats swing out allowing the user to stand up and board the unit. The caregiver then flips the split seats back into place to form a secure and comfortable seat for the user to sit on for the remainder of the transfer.

A user who qualifies to use the **Alerta MoveAssist** must have enough leg and lower body strength to stand up and remain in the standing/ sitting position. Adequate arm strength is required if the user must use the crossbar. Users who satisfy these criteria and have difficulty walking or present a potential for falling down will find the **Alerta MoveAssist** a useful and safe transfer device. For users who lack these requirements, a standard hoist is recommended.

TECHNICAL SPECIFICATION

- > Adjustable legs for ease of manoeuvring
- > Ideal for bed-chair-wheelchair-commode transfers
- > Perfect for room to room transfers
- > MoveAssist can be wheeled directly up to a toilet, offering enhanced dignity
- > Braked rear castors
- > Max user weight 200Kgs
- > Sling hooks for extra security
- > Leg height clearance: 10cm
- > Overall Length: 90cm



ALT-PL-CHB
T-Motion Fixed Charging
Base Available

ALERTA STAND 200 POWERLIFTER

ALT-PLS200 - Alerta Stand 200 Powerlifter - a lifting hoist for assisting users with standing from a seated position. Suitable for users up to 200kg. Alerta Stand lifter is designed to be used in hospital, nursing and care home environments.

It provides security, functionality and outstanding reliability, easy to operate and easy to move. The Alerta Stand lifter also gives users a very restful and calming lifting experience.

The wide range of features make the Alerta Stand lifter a highly versatile and cost-effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Braked rear castors
- > Disassembles into two parts
- > Adjustable arm with 5 positions
- > Adjustable legs for ease of manoeuvring
- > Max user weight: 200kg

Patient Transfer Aids



* 1 year warranty on batteries



- Safe working load
135kg
- Maximum lift height
1350mm
- Minimum lift height
330mm
- Overall width
530mm
- Overall length
1130mm



ALT-PL-CHB
T-Motion Fixed Charging
Base Available

ALERTA POWERLIFTER MINI 135

ALT-PL135 - Alerta Powerlifter Mini 135 – a transfer hoist for handling and transferring users. Suitable for users up to 135 kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



* 1 year warranty on batteries



- Safe working load
150kg
- Maximum lift height
1470mm
- Minimum lift height
340mm
- Overall width
540mm
- Overall length
1230mm



ALT-PL-CHB
T-Motion Fixed Charging
Base Available

ALERTA POWERLIFTER MIDI 150

ALT-PL150 - Alerta Powerlifter Midi 150 – a transfer hoist for handling and transferring users. Suitable for users up to 150 kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



* 1 year warranty on batteries



Safe working load
175kg



Maximum lift height
1590mm




Minimum lift height
450mm



Overall width
630mm



Overall length
1360mm

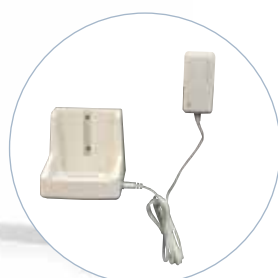

ALT-PL-CHB
T-Motion Fixed Charging
Base Available

ALERTA POWERLIFTER MAXI 175

ALT-PL175 - Alerta Powerlifter Maxi 175 – a transfer hoist for handling and transferring users. Suitable for users up to 175 kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



* 1 year warranty on batteries



Safe working load
175kg



Maximum lift height
1590mm



Minimum lift height
450mm



Overall width
630mm



Overall length
1360mm



ALT-PL-CHB
T-Motion Fixed Charging
Base Available

ALERTA POWERLIFTER MAXI 175-E

ALT-PL175-E - Alerta Powerlifter Maxi 175-E is a transfer hoist with electric legsread, used for handling and transferring users. Suitable for users up to 175kg. **Alerta Powerlifter** hoists are designed to be used in hospital, nursing and care home environments.

They provide security, functionality and outstanding reliability, easy to operate, easy to move and incorporate anti-entrapment techniques. **Alerta Powerlifter** hoists also give users a very restful and calming lift experience with a thermoplastic rubber coated, anti-swing spreader bar.

The wide range of features make the **Alerta Powerlifter** hoists a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

TECHNICAL SPECIFICATION

- > TiMotion control system with emergency stop button
- > Sound/light alarm battery level
- > Operate hoist from handset or top of control box
- > Spreader bar with 360° rotation and anti-swing
- > Braked rear castors
- > Disassembles into two parts



ALT-ST002 - Alerta Standing Turner - a transfer device to assist carer/attendant with transfer of the user between seat, bed, toilet, commode, wheelchair by standing and turning. Alerta Standing Turner promotes user engagement and enhances mobility.

It benefits from Alerta's GlideTurn easy swivel base, height adjustable handle, height adjustable padded knee pads, sturdy steel frame and tool free assembly.

To use this device, the user requires enough arm, leg and body strength to stand up and remain in the standing position. For users who lack these attributes, an Alerta Standaid hoist or MoveAssist is recommended.

Do not use this device to transport a user. When patient is on board DO NOT tilt or lever the Standing Turner, the whole footplate should remain in contact with the floor at all times.

TECHNICAL SPECIFICATION

- > Base width: 41cm
- > Handle width: 47cm
- > Handle Height Adjust: 89cm to 126cm
- > Product weight: 17.5kg
- > 2 x 7cm rear castors for transporting between users
- > Safe Working Load: 160kg



The **Alerta Deluxe Support Belt** with leg straps is designed to assist with the safe transfer of individuals from one position to another - such as bed, wheelchair, toilet, dining transfers, or walking with the user - without the need to hold the individual's clothing or body, enabling freedom of movement and a greater sense of independence.

Alerta Deluxe Support Belt is NOT for use as a device for patient lifting.

ALT-216 ALERTA DELUXE SUPPORT BELT, MEDIUM

ALT-217 ALERTA DELUXE SUPPORT BELT, LARGE

KEY FEATURES & BENEFITS

- > Multiple transfer type uses
- > Gait training
- > Promotes user independence - no physical touching
- > Non-slip inside surface
- > Detachable leg straps included
- > Adjustable Velcro & buckle fastening
- > 10 firm handhold positions

WASHING INSTRUCTIONS

- > Material: 100% Polyester
- > Do not dry clean
- > Machine washable at up to 30°C (velcro must be closed)

SIZES

ALT-216 Medium 90x12cm
(to suit waist sizes: 62cm (24") to 88cm (34"))

ALT-217 Large 115x12cm
(to suit waist sizes: 82cm (32") to 110cm (43"))

Wheelchairs



ALT-1100

ALERTA CAR TRANSIT WHEELCHAIR, CRASH TESTED

The **Alerta Car Transit Wheelchair** is crash tested and comes with drop down armrests, swing away detachable footrests, arm pads, lap belt and rear brakes, with 30cm rear wheels and 20cm front wheels, all with puncture proof tyres. Seat is 40cm deep x 43cm wide.

TECHNICAL SPECIFICATION

- > 30cm rear wheels
- > 20cm front wheels
- > Seat size: 40cm deep x 43cm wide
- > Max user weight: 115kg
- > Product weight: 15.5kgs

KEY FEATURES & BENEFITS

- > Crash tested
- > Flip up table-entry arms
- > Padded armrests
- > Lap belt
- > Swing away detachable footrests
- > Puncture proof tyres
- > Braked rear wheels
- > Heel loops



ALT-1300

ALERTA SELF-PROPELLED WHEELCHAIR, CRASH TESTED

The **Alerta Self-Propelled Wheelchair** is crash tested and comes with self-propelling rear wheels, detachable arms, padded armrests, lap belt, swing away and detachable footrests, heel loops and rear brakes, with 60cm rear wheels and 20cm front wheels, all with puncture proof tyres. Seat is 40cm deep x 45cm wide.

TECHNICAL SPECIFICATION

- > 60cm rear self-propelled wheels
- > 20cm front wheels
- > Seat size: 40cm deep x 45cm wide
- > Max user weight: 120kg
- > Product weight: 17kg

KEY FEATURES & BENEFITS

- > Self-propelled rear wheels
- > Wave profile self-propel rim
- > Flip up table-entry arms
- > Padded armrests
- > Lap belt
- > Swing away detachable footrests
- > Puncture proof tyres
- > Rear wheel brakes
- > Heel loops



ALT-1100/WHD

ALERTA CAR TRANSIT WIDE HEAVY-DUTY ALUMINIUM WHEELCHAIR

The **Alerta Car Transit Wide Heavy-Duty Aluminium Wheelchair** comes with folding back, lift up / swing away armrests, table-entry armrests, swing away detachable footrests, 30cm braked rear wheels and 20cm front wheels, all with puncture proof tyres. The wide heavy-duty wheelchair has anti-tip bars at the back beside the rear wheels to help prevent tipping backwards when crossing an obstacle.

TECHNICAL SPECIFICATION

- > 30cm braked rear wheels
- > 20cm front wheels
- > Seat size: 46cm deep x 60cm wide
- > Max user weight: 182kg
- > Product weight: 14.8kg

KEY FEATURES & BENEFITS

- > Swing away / Lift up armrests
- > Table entry padded arms
- > Handle lock safety feature
- > Folding backrest
- > Hand brakes
- > Rear anti-tip bars
- > Lap belt
- > Swing away detachable footrests
- > Puncture proof tyres
- > Heel loops
- > Carry handles once folded
- > Double crossbar for extra strength
- > Reinforced extra thick upholstery
- > Rear seat pocket



SHOWER CHAIR | COMMODE CHAIR | TRANSFER CHAIR | TOILETING CHAIR

ALT-SC150

ALERTA AQUA SHOWER COMMODE CHAIR

ALT-SC150/SP

ALERTA AQUA SHOWER COMMODE CHAIR, SELF-PROPELLED

ALT-SC250

ALERTA AQUA SHOWER COMMODE CHAIR, BARIATRIC



The **Alerta Aqua Chair** combines 4 chairs in 1, performing seamlessly as a shower chair, a commode chair, a toileting chair and a transfer chair.

While it can be used for showering, transferring and as a mobile commode, the pot can be removed and the chair can be wheeled to the toilet, making this a very easy to use toileting chair.

The Alerta Aqua is stylish, effortlessly manoeuvrable and manufactured from easy to clean plastic coated steel.

A selection of soft padded seat sections are also available as accessories if desired.

KEY FEATURES & BENEFITS

- > Available with and without self-propelled wheels
- > Five position footrest
- > 4 individually braked wheels
- > Weight capacity:
 - > ALT-SC150: 150kg
 - > ALT-SC150/SP: 150kg
 - > ALT-SC250: 240kg
- > Removable arms
- > 2 year warranty
- > Seat Size:
 - > ALT-SC150: 45 x 43cm
 - > ALT-SC150/SP: 45 x 43cm
 - > ALT-SC250: 45 x 65.7cm





ALT-SC300

ALERTA AQUA COMMUNITY, TILT-IN-SPACE SHOWER COMMODE CHAIR

ALT-SC350

ALERTA AQUA TILT-IN-SPACE HI-LO SHOWER COMMODE CHAIR

The **Alerta Aqua Tilt-in-Space Shower Commode Chairs** are for use in Care Homes, Hospitals and homecare environments. They are designed specifically for assisted hygiene routines, showering and toileting.

The transfer friendly design offers safe solutions for residents with differing mobility levels, covering needs ranging from a bedside commode to a mobile, multifunction shower and commode chair. When used as a toilet chair, it can be placed over a conventional toilet or used in a freestanding position. The Aqua Chair Tilt-In-Space has an adjustment controlled by gas spring for the seat angle for active and passive seating and is designed for repeated use.

KEY FEATURES & BENEFITS

- > Easy, simple, safe to operate Tilt-in-Space
- > Removable swing-away footrests with angle-adjustable footplate
- > Tool-less seat height adjustment
- > Adjustable and detachable footrests
- > Multi-angle directional and height adjustable head support
- > 5" casters with wheel & swivel locking mechanism
- > Commode pan with lid included
- > Moulded seat and backrest for high comfort
- > All upholstery is waterproof, fire-proof and non-toxic

ALT-SC300

- > Adjustable height footrests
- > 1 Year warranty

ALT-SC350

- > Rear anti-tip bars
- > Width & angle adjustable pivoting upper body side supports with quick release

ALT-SC300 TECHNICAL SPEC

- > Overall height: 102cm
- > Seat size: 40 x 46cm
- > Seat height: 59.5 - 64.5cm
- > Backrest height: 49cm
- > Max user weight: 175kg

ALT-SC350 TECHNICAL SPEC

- > Adjustable height: 133-159cm
- > Seat size: 46 x 42.5cm
- > Seat height: 52.5 - 62.5cm
- > Max seat angle forwards: 5°
- > Max seat angle backwards: 36.5°
- > Max user weight: 160kg

Showering & Toileting

ALT-BE002

ALERTA MOBILE SHOWER CHAIR



This mobile shower chair has 2 rear braked castors. 136kg weight capacity.

Boxed in 2's

ALT-BE004

ALERTA MOBILE SHOWER CHAIR, ADJUSTABLE HEIGHT



This mobile shower chair has adjustable height and 2 rear braked castors. 136kg weight capacity.

Seat Height: 44.5-54.5cm

Boxed individually

ALT-BE005

ALERTA STATIONARY SHOWER CHAIR, ADJUSTABLE HEIGHT



This stationary shower chair has adjustable height. 136kg weight capacity.

Seat Height: 41.5-53.3cm

Boxed in 4's



The **Alerta Aluminium Shower Stool & Aluminium Shower Stool with Back** enjoys a simple but cosmetically attractive design, with a lightweight aluminium frame. It is the perfect companion to aid support and transfer whilst showering, benefitting from adjustable height for ease of use, integral handles, and the added safety features of moulded drainage holes and nonslip rubber feet.

Constructed from corrosion proof materials, extra stability is afforded with angled leg design and suction style foot tips. It also benefits from a left and right-hand side built in shower head holder.

Max user weight: 130kg

Boxed individually



ALT-BE010 ALERTA ALUMINIUM SHOWER STOOL

ALT-BE011 ALERTA ALUMINIUM SHOWER STOOL WITH BACK



The **Alerta Raised Toilet Seats** increase the toilet seat height for easier toilet access in nursing homes and home care.

Max user weight: 200kg

Boxed in 12's

ALT-BE006 ALERTA RAISED TOILET SEAT, 2 INCHES

ALT-BE007 ALERTA RAISED TOILET SEAT, 4 INCHES

ALT-BE008 ALERTA RAISED TOILET SEAT, 6 INCHES

ALT-BE009**ALERTA STACKING COMMODE, ADJUSTABLE HEIGHT**

The **Alerta Stacking Commode** is manufactured from powder coated tubular steel and comes with a comfortable padded upholstery on the seat and back, which is heat sealed allowing it to be cleaned without the risk of water ingress. The toilet seat is oval shaped with a clip on/clip off design. A rubber buffer sits in the seat to avoid seat abrasion against the frame. The commode holds a 5 litre round potty.

Max user weight: 180kg

Seat Height: 45-60cm

Boxed in 3's

ALT-BE001**ALERTA TOILET SEAT AID, ADJUSTABLE HEIGHT**

The **Alerta Toilet Seat Aid** is plastic coated and has a large seating area. It features cut-outs at the front and rear for easy cleaning access. The toileting seat simply clicks into position on the frame and rotates when required.

Max user weight: 165kg

Seat Height: 48-62cm.

Boxed in 2's

ALT-BE003**ALERTA PORTABLE TOILET SURROUND, ADJUSTABLE HEIGHT & WIDTH**

The **Alerta Portable Toilet Seat Surround** with height extension legs is designed to make using the toilet easier for users in nursing homes and home care.

Height Range: 68-80.7cm

Outer Product Width: 55-60cm

Boxed individually

ALT-1200**ALERTA COMMODE & TRANSFER CHAIR**


The **Alerta Commode & Transfer Chair** comes with drop down armrests, swing away detachable footrests, arm pads, four braked castors and wipeable PVC seat and back. Seat is 53cm deep x 46cm wide.

- > Stable, robust powder coated steel frame
- > Drop down padded armrests
- > Swing away detachable footrests
- > Full length arm pads
- > Push handle for optimal positioning
- > 4 solid steel 5" swivel braked castors
- > Fold-up footrests
- > Wipeable PVC seat and back
- > Max user weight: 125kg

**ALT-1400****ALERTA ADJUSTABLE HEIGHT COMMODE & TRANSFER CHAIR**

The **Alerta Adjustable Height Commode & Transfer Chair** comes with drop down armrests, swing away detachable footrests, arm pads, four braked castors and wipeable PVC seat and back. Seat is 53cm deep x 46cm wide.

- > Stable, robust powder coated steel frame
- > Drop down padded armrests
- > Swing away detachable footrests
- > Full length arm pads
- > Push handle for optimal positioning
- > 4 solid steel 5" swivel braked castors
- > Fold-up footrests
- > Wipeable PVC seat and back
- > Height ratio: 52.5 - 65cm
- > Max user weight: 125kg

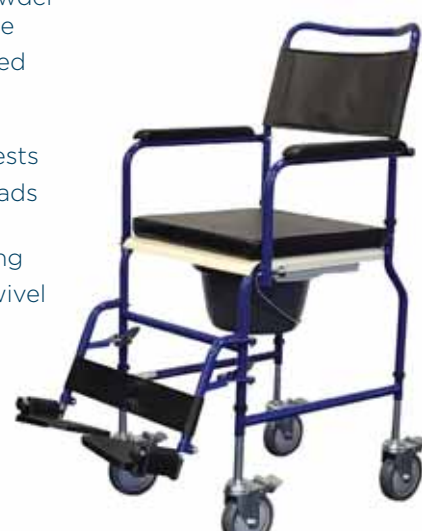



Evacuation



ALT-E001 EVACUATION SHEET

The **Alerta Evacuation Sheet** fixes permanently to the underside of a flexible foam mattress by means of strong elastics. It has safety belts to hold the user on the mattress and pulling straps at each end.

The combination of slippery sheet and flexible mattress will allow controlled and easy access through doorways and downstairs.

Max user weight: 350kg

Body strap length: 160cm

Head/Foot evac total strap length: 150cm



Wall mounted
zip storage bag

ALT-E002 EVACUATION SLEDGE

The **Alerta Evacuation Sledge** is used for transporting a user to safety during an evacuation. The user is placed on the sledge, secured with the chest and leg straps and pulled to safety using the straps at each end.

The flexible sledge and slippery underneath will allow controlled and easy access through doorways and downstairs.

Max user weight: 350kg

Body strap length: 106cm



ALT-E004 BARIATRIC EVACUATION SLEDGE

The **Alerta Bariatric Evacuation Sledge** is used for transporting a user to safety during an evacuation. The user is placed on the sledge, secured with the chest and leg straps and pulled to safety using the straps at each end.

The flexible sledge and slippery underneath will allow controlled and easy access through doorways and downstairs.

Max user weight: 350kg

Width: 100cm

Body strap length: 160cm

Head/Foot evac total strap length: 212cm



ALT-E003 AIR MATTRESS EVACUATION SHEET

The **Alerta Air Mattress Evacuation Sheet** is intended to be placed permanently underneath an alternating air mattress ready to be used for transporting a user to safety during an evacuation. Alerta's Air Mattress Evacuation Sheet is suitable for both Overlay and Replacement air mattress systems.

The Air Mattress Evacuation Sheet enables the whole mattress to be turned into an evacuation sledge in the event of an evacuation. The flexible and slippery underneath will allow controlled and easy access through doorways and downstairs.

Fits Overlay's and Deep cell Replacement air mattresses.



200cm long x 90cm wide





Foam, Gel &
Static Systems

Product Overview

All foam products are supplied vacuum roll packed

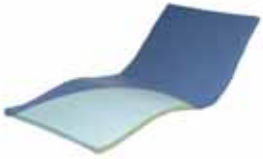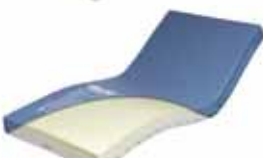	 Alerta[®] sensaflex underlay	Foam Mattress Underlay	
	 Alerta[®] sensaflex500	Foam Mattress	
	 Alerta[®] sensaflex1000	Foam Mattress	
	 Alerta[®] sensaflex3000	Memory Foam Mattress	
	 Alerta[®] sensaflex4000	Gel-Topped Memory Foam Mattress	
	 Alerta[®] sensaflex4500	Heel Slope Mattress Memory & Gel Foam	
	 Alerta[®] sensaflex200	Foam Cushion	
	 Alerta[®] sensaflex250	Memory Foam Cushion	
	 Alerta[®] sensaflex300	Memory Foam Cushion	
	 Alerta[®] sensaflex350	Memory Foam Bariatric Cushion	

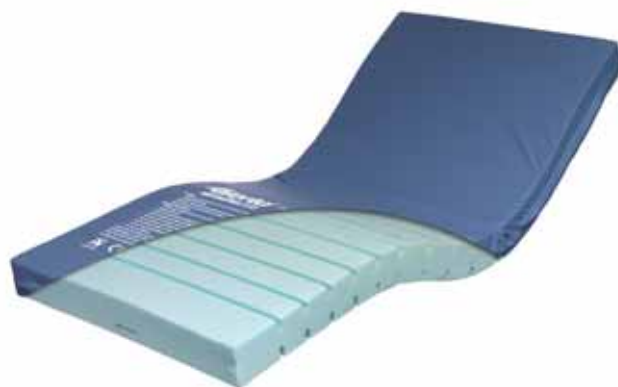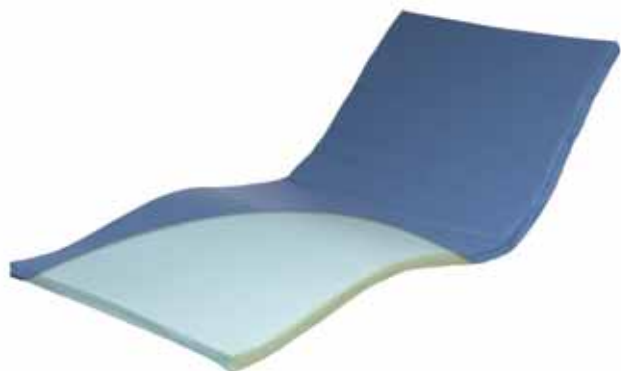

ALT-UND – Alerta 2 inch Foam Mattress Underlay

– a foam mattress underlay for use underneath an air mattress system in the event of power failure or during the use of CPR. The underlay provides effective comfort and care for users in hospital, nursing and care home environments.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

TECHNICAL SPECIFICATION

Mattress Specification

- > Mattress dimensions: 198 x 90 x 5 cm (78 x 35.5 x 2")
- > Multi-stretch and vapour permeable PU cover
- > Waterproof and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5

ALT-500 – Alerta Sensaflex 500 – a medium risk profiling foam mattress built with castellated foam. The sophisticated mattress design provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

TECHNICAL SPECIFICATION

Mattress Specification

- > Mattress dimensions: 198 x 90 x 12 cm (78 x 35.5 x 4.8")
- > Max load capacity: 120kg (265lb)
- > Min load capacity: 50kg (110lb)
- > Multi-stretch and vapour permeable PU cover
- > Waterproof and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5



To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

- > Multi-stretch and vapour permeable PU cover
- > Waterproof and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5



ALT-1000 – Alerta Sensaflex 1000 – a high risk profiling foam mattress built with castellated foam. The sophisticated mattress design provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

TECHNICAL SPECIFICATION

Mattress Specification

- > Mattress dimensions:
198 x 90 x 14 cm
(78 x 35.5 x 5.5")
- > Max load capacity: 180kg (397lb)
- > Min load capacity: 40kg (88lb)



ALT-3000 – Alerta Sensaflex 3000 – a high risk profiling memory foam mattress. The sophisticated mattress design provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

TECHNICAL SPECIFICATION

ALT-3000 Mattress Specification

- > Mattress dimensions:
198 x 90 x 14 cm
(78 x 35.5 x 5.5")
- > Max load capacity: 180kg (397lb)
- > Min load capacity: 40kg (88lb)

ALT-3000/4 Mattress Specification

- > Mattress dimensions:
198 x 120 x 16 cm
(78 x 47 x 6.3")
- > Max load capacity: 200kg (441lb)
- > Min load capacity: 40kg (88lb)

**BARIATRIC
120CM
WIDE**



LOW
MINIMUM
USER
WEIGHT



ALT-4000 - Alerta Sensaflex 4000 – a very high risk profiling memory foam mattress built with a high density base, foam centre and a castellated gel memory foam top. The mattress is designed with a section of v-cut memory foam which creates additional air movement for enhanced pressure relieving properties, providing effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

TECHNICAL SPECIFICATION

Mattress Specification


- > Mattress dimensions:
198 x 90 x 16 cm
(78 x 35.5 x 6.3")
- > Max load capacity: 200kg (441lb)
- > Min load capacity: 30kg (66lb)
- > Multi-stretch and vapour permeable PU cover
- > Waterproof and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5



All foam mattresses (except 2" underlay) and cushions are shipped vacuum rolled for ease of transport and minimising storage



Heel Slope Mattress

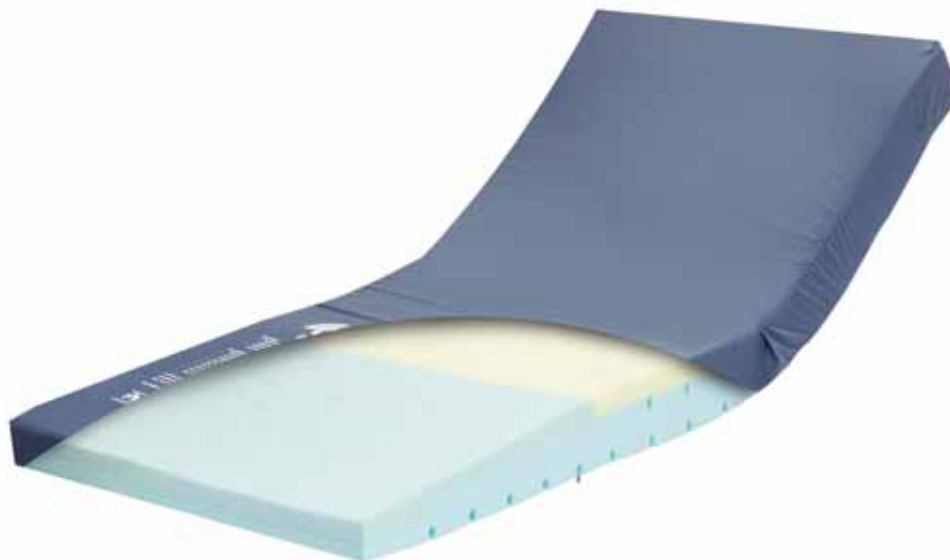

ALT-4500 - Alerta Sensaflex 4500 - is a very high risk profiling memory foam mattress built with a high density base, memory foam top, and gel-infused memory foam sloping heel section. It is designed with a section of v-cut base foam which creates additional air movement for enhanced pressure relieving properties, providing effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

- > Heel slope: assists pressure redistribution by alleviating pressure from the high risk heel area.
- > SensaGel® heel zone: gel-infused memory foam to prevent pressure sores on the heel.

To launder the cover in a washing machine, set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

- > Multi-stretch and vapour permeable PU cover
- > Antimicrobial cover with white underside
- > Waterproof and machine washable cover
- > BS7177: Crib 5



TECHNICAL SPECIFICATION

Mattress Specification

- > Mattress Dimensions:
198 x 90 x 8-14cm
(78 x 35.5 x 3.1-5.5")
- > Max Load: 180kg (397lb)
- > Min Load: 40kg (88lb)



To launder the cover in a washing machine set the main wash to 65°C for no less than 10 minutes or 95°C for no less than 3 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

alerta[®] **sensaflex200**



ALT-200 – Alerta Sensaflex 200 – a medium risk foam cushion built with castellated foam. The foam cushion provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

TECHNICAL SPECIFICATION

- > Cushion dimensions:
43 x 43 x 10 cm (17 x 17 x 4")
- > Max Load: 150kg (331lb)
- > Min Load: 40kg (88lb)



- > Multi-stretch and vapour permeable PU cover
- > Waterproof and machine washable cover
- > Antimicrobial cover with white underside
- > BS7177: Crib 5

alerta[®] **sensaflex250**



ALT-250 – Alerta Sensaflex 250 – a high risk memory foam cushion which provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

TECHNICAL SPECIFICATION

- > Cushion dimensions:
43 x 43 x 10 cm (17 x 17 x 4")
- > Max Load: 150kg (331lb)
- > Min Load: 40kg (88lb)



alerta[®] **sensaflex300**



ALT-300 – Alerta Sensaflex 300 – a very high risk memory foam cushion built with a high density base, foam centre and a castellated gel memory foam top. The cushion provides effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments.

TECHNICAL SPECIFICATION

- > Cushion dimensions:
43 x 43 x 10 cm (17 x 17 x 4")
- > Max Load: 150kg (331lb)
- > Min Load: 40kg (88lb)



alerta[®] **sensaflex350**




ALT-350 – Alerta Sensaflex 350 – a high risk bariatric memory foam cushion which provides effective comfort, care and pressure redistribution for bariatric users in hospital, nursing and care home environments.

TECHNICAL SPECIFICATION

- > Cushion Dimensions:
45 x 56 x 12 cm (17.7 x 22 x 5")
- > Max Load: 280kg (617lb)
- > Min Load: 50kg (110lb)




Gel Cushions

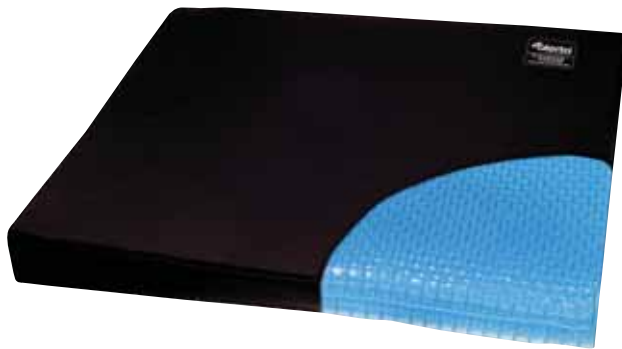

ALT-204 - Alerta GelCube Cushion – a patented dual honeycomb air-gel cushion to provide advanced comfort, care and pressure redistribution for users in hospital, nursing and care home environments. It is built with an intelligent dual honeycomb structure to provide optimal blood circulation and skin health, and balanced pressure distribution. The combination of flexible air-gel cells provides a non-slip surface for superior comfort over an extended period of use.

The bottom of the cover is coated with silicon for non-slip and the top of the cover is made with polyester fabric.


- > Max Load: 200kg
- > Cushion Dimensions: 46 x 42 x 5 cm (18 x 16.5 x 2")

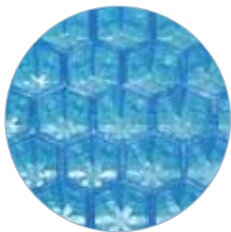

Patented dual honeycomb structure



Intelligent design for balanced pressure distribution

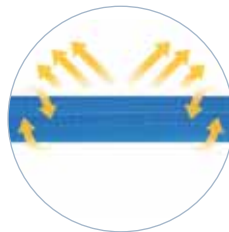

Provides enhanced blood circulation and skin health



Very easy to clean



Packaged for retail



Cover provided





Packaged for retail



ALT-201 - Alerta Gel Cushion – a gel cushion containing polyurethane SensaGel which acts as a layer of artificial fat, to provide effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments. It can slide and move laterally to effectively reduce shear forces and is guaranteed not to harden over time, neither leak nor dry if accidentally punctured.

- > PU and fire resistant classified M1 cover
- > Waterproof and vapour permeable cover
- > Max Load: 100kg
- > Cushion Dimensions: 43 x 43 x 2.5 cm (17 x 17 x 1")

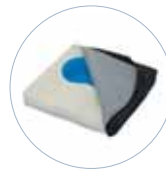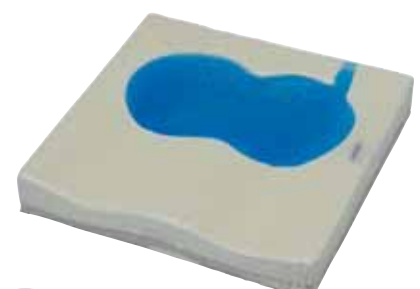

Cover provided



Packaged for retail



ALT-202 - Alerta Gel-Visco Cushion – a gel-visco cushion built with a contoured moulded viscoelastic foam base, connected to a polyurethane SensaGel insert with two densities, to provide effective comfort, care and pressure redistribution for users in hospital, nursing and care home environments. The viscoelastic SensaGel provides maximum relief to bony areas and the shape ensures a good posture while seated and prevents users from sliding forward.

- > PU and fire resistant classified M1 cover
- > Waterproof and vapour permeable cover
- > Max Load: 155kg
- > Cushion Dimensions: 43 x 43 x 8cm (17 x 17 x 3")



ALT-203 - Alerta Air Cushion with Pump

– an aid in the prevention of pressure ulceration for all users including those assessed as high risk, and for the treatment of existing stage 1 and 2 pressure ulcers. Clinical supervision is advised if skin damage is severe. Once cleaned this product can be reused by another user.

Cushion Dimensions: 46 x 46 x 3 cm (18 x 18 x 1.2")



ALT-210 - Alerta Air Heel Bootees with Strap – provides effective protection for users by offloading pressure from the heel to minimise the risk of pressure ulceration and injury to the foot, including those assessed as very high risk. If existing pressure related tissue injury is severe then clinical supervision is recommended.



ALT-211 Alerta Deluxe Heel Bootee - is designed to support and elevate the heel, providing relief to individuals who experience pain, discomfort or pressure on the heels and can be worn by the user in many positions.

The heel can be a very difficult area to protect against pressure ulcers. Alerta's Deluxe Heel Bootee provides superior deluxe comfort and pressure relief, through genuine heel offloading, reducing the possibility of a pressure ulcer developing and helping existing pressure ulcers to heal.

The cool brushed nylon fabric limits perspiration, providing a comfortable surface and eliminates friction while the superior deluxe memory foam internal padding reduces pressure risk.

Retail Packaged.



ALT-212 Soft Heel Protectors - will help reduce the risk of heel pressure ulcers by eliminating friction and reducing pressure. The cool natural cotton limits perspiration while the padding provides a comfortable surface and raises the foot slightly to reduce contact.

Alerta Heel Protector fits most sizes; simply fasten the straps to provide a comfortable fit.

Fabric Content: Cotton surface.

Siliconised polyester hollow fibre fill.

Retail Packaged.



Fall Prevention

Product Overview

		Premium floor mat, Wired and Wireless	
		Premium floor mat, Wired	
		Crash Mat, Wireless	
		Floor mat, Wired and Wireless	
		Bed mat system, Wired and Wireless	
		Chair mat system, Wired and Wireless	
		PIR motion sensor, Wired and Wireless	
		Wireless Window and Door Sensor	
		Nurse call button, Wired and Wireless	
		WetSense Chair Mat	
		WetSense Bed Mat	



AVAILABLE IN 5 COLOURS



The **Deluxe Alertamat** is a robust long-life luxury carpet Alertamat; the world's first and only genuine carpeted Alertamat. Available in 5 colours to suit an existing colour scheme, the **Deluxe Alertamat** provides a more refined touch for the care environment.

The **Deluxe Alertamat** is most effective when placed at the bedside or doorway; activating the alarm when a user leaves their safe zone by stepping on the mat. This alerts the caregiver to ascertain the user's safety. Selecting a colour that blends with the existing flooring makes it an ideal solution for situations where a dementia user may otherwise try to avoid the mat.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.

- > Genuine non-slip base
- > Luxury carpet finish
- > Choice of 5 colours
- > Long lasting welded seals and strong wiring
- > Large size avoids chance of overstepping
- > Dementia friendly
- > Dimensions: 60 x 90 cm (23.5 x 35.5")



WIRED DELUXE ALERTAMAT

DAM - Wired Deluxe Alertamat - can be connected to a nurse call system with a 3 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (sold separately).

- > Plug-and-play system
- > 0, 2 or 4 second delay to reduce false alarms
- > System includes:
1x Deluxe Alertamat and 1x 3 metre cable

The Deluxe Alertamat comes with a unique built-in "Y" cable, providing a spare socket for other Alerta devices to be connected - using only one nurse call point.



Connect to nurse call system or stand-alone alarm monitor

Breakaway cable to reduce trip hazard and prevent damage



WIRELESS DELUXE ALERTAMAT

W-DAM-T - Wireless Deluxe Alertamat - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > 95m wireless range from transmitter to receiver
- > 0, 2 or 4 second delay to reduce false alarms
- > System includes: 1x Deluxe Alertamat and 1x Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

Alerta Wireless Transmitter

- > Battery powered
- > Includes 2x 1.5v AAA batteries



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





WAM-FHOAK - Wood Alertamat - a robust long-life luxury Alertamat with a wood safety floor finish; the world's first and only genuine wood finish Alertamat. The Wood Alertamat provides a more refined touch for the care environment.

The **Wood Alertamat** is most effective when placed at the bedside or doorway; activating the alarm when a user leaves their safe zone by stepping on the mat. This alerts the caregiver to ascertain the user's safety. The light colour makes it an ideal solution for situations where a dementia user may try to avoid that mat.

It can be connected to a nurse call system with a 3 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (sold separately).

The **Wood Alertamat** comes with a unique built-in "Y" cable, providing a spare socket for other Alerta devices to be connected - using only one nurse call point.

KEY FEATURES AND BENEFITS

- > Plug-and-play system
- > Genuine non-slip base
- > Luxury wood finish
- > Long lasting welded seals and strong wiring
- > Hygienic, antibacterial easy-clean material
- > Large size avoids chance of overstepping
- > Dementia friendly
- > Dimensions: 55 x 80 cm (21.5 x 31.5")
- > System includes: 1x Wood Alertamat and 1x 3 metre cable



Connect to nurse call system or stand-alone alarm monitor



Breakaway cable to reduce trip hazard and prevent damage



SCAN TO WATCH PRODUCT VIDEO

**WORLD'S FIRST
100% WIRELESS
ALARM CRASH MAT**

W-ACM-T - Alertamat Crash Mat with Wireless Nurse Call Alert - is designed to be placed next to the bed to minimise injury if the user should fall from bed and activates the nurse call system when a fall occurs.

The mat incorporates sensor pads and a Wireless Transmitter within the foam which connects to the Alerta Wall Point Receiver. It can be connected to a nurse call system for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Wall Point Receiver only. The volume can be deactivated for silent in-room operation when connected to a nurse call system.

This helps prevent a user remaining on the mat undetected for an unknown length of time after falling, by alerting the caregiver to provide timely assistance.

The Alertamat Crash Mat is quick and easy to position and benefits from a PU wipe clean surface. When unfolded the mat is 180cm long x 60cm wide x 6cm high. The mat folds in half for easy storage and transportation.

Patent Number: 2203815.2

KEY FEATURES AND BENEFITS

- > Get alerted when a user falls from bed
- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > 95m wireless range from transmitter to receiver
- > Genuine non-slip base
- > Dimensions: 180 x 60cm (71 x 23.5")
- > Thickness: 6cm (2.5")
- > System includes: 1x Alertamat Crash Mat and 1x Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

Alerta Wireless Transmitter

- > Battery powered
- > Includes 2x 1.5v AAA batteries



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





The **Floor+ Alertamat** is a robust long-life Alertamat and is most effective when placed at the bedside or doorway; activating the alarm when a user leaves their safe zone by stepping on the mat. This alerts the caregiver to ascertain the user's safety. The light colour makes it an ideal solution for situations where a dementia user may try to avoid the mat.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.



WATCH WIRELESS
PRODUCT VIDEO

- > Genuine non-slip base
- > Long lasting welded seals and strong wiring
- > Hygienic, antibacterial easy-clean material
- > Large size avoids chance of overstepping
- > Dementia friendly
- > Dimensions: 60 x 90 cm (23.5 x 35.5")



WIRED FLOOR+ ALERTAMAT

PAM - Wired Floor+ Alertamat - can be connected to a nurse call system with a 3 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (sold separately).

- > Plug-and-play system
- > 0, 2 or 4 second delay to reduce false alarms
- > System includes:
1x Floor+ Alertamat and 1x 3 metre cable

The Floor+ Alertamat comes with a unique built-in "Y" cable, providing a spare socket for other Alerta devices to be connected - using only one nurse call point.



Connect to nurse call system or
stand-alone alarm monitor

Breakaway cable to
reduce trip hazard
and prevent damage



WIRELESS FLOOR+ ALERTAMAT

W-PAM-T - Wireless Floor+ Alertamat - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > 95m wireless range from transmitter to receiver
- > 0, 2 or 4 second delay to reduce false alarms
- > System includes: 1x Floor+ Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

Alerta Wireless Transmitter

- > Battery powered
- > Includes 2x 1.5v AAA batteries



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





The **Bed Alertamat** is most effective when concealed under the bed sheet on top of the mattress; activating the alarm when a user leaves their bed. This alerts the caregiver to ascertain the user's safety. The alarm is reset automatically when the user returns to the bed.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.

- > Hygienic, antibacterial easy-clean material
- > Dimensions: 25 x 76 cm (10 x 30")



WATCH WIRELESS
PRODUCT VIDEO



WIRED BED ALERTAMAT

BAM - Wired Bed Alertamat - can be connected to a nurse call system with a 2 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (supplied as part of the system).

- > Plug-and-play system
- > System includes: 1x Bed Alertamat, 1x Alerta Alarm Monitor and 1x 2 metre cable

Alerta Alarm Monitor

- > High, low and silent volume control
- > Multiple choice alarm tone
- > 0, 2 or 4 second delay to reduce false alarms
- > Low battery alert
- > Mains or battery powered
- > Takes 1x 9v battery




Battery operated



Mains operated



Alerta Alarm Monitor

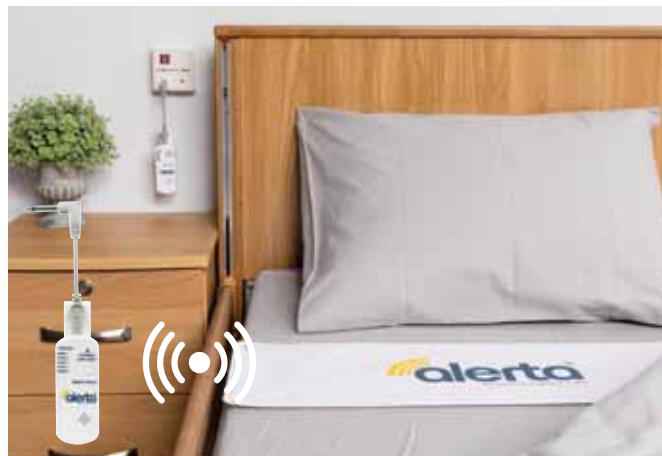

WIRELESS BED ALERTAMAT

W-BAM-T - Wireless Bed Alertamat - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > 95m wireless range from transmitter to receiver
- > System includes: 1x Bed Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

Alerta Wireless Transmitter

- > Battery powered
- > Includes 2x 1.5v AAA batteries
- > 0, 2 or 4 second delay to reduce false alarms



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries





The **Chair Alertamat** is most effective when concealed under the chair cushion; activating the alarm when a user leaves their chair. This alerts the caregiver to ascertain the user's safety. The alarm is reset automatically when the user returns to the chair.

Both **Wired** and **Wireless** systems can be connected to a nurse call system for centralised activation, or used on their own, independent of a nurse call system to give localised alerts.

- > Hygienic, antibacterial easy-clean material
- > Dimensions: 25 x 38 cm (10 x 15")



WATCH WIRELESS
PRODUCT VIDEO



WIRED CHAIR ALERTAMAT

CAM - Wired Chair Alertamat - can be connected to a nurse call system with a 2 metre cable provided for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Alarm Monitor (supplied as part of the system).

- > Plug-and-play system
- > System includes: 1x Chair Alertamat, 1x Alerta Alarm Monitor and 1x 2 metre cable

Alerta Alarm Monitor

- > High, low and silent volume control
- > Multiple choice alarm tone
- > 0, 2 or 4 second delay to reduce false alarms
- > Low battery alert
- > Mains or battery powered
- > Takes 1x 9v battery




Battery operated



Mains operated



Alerta Alarm Monitor

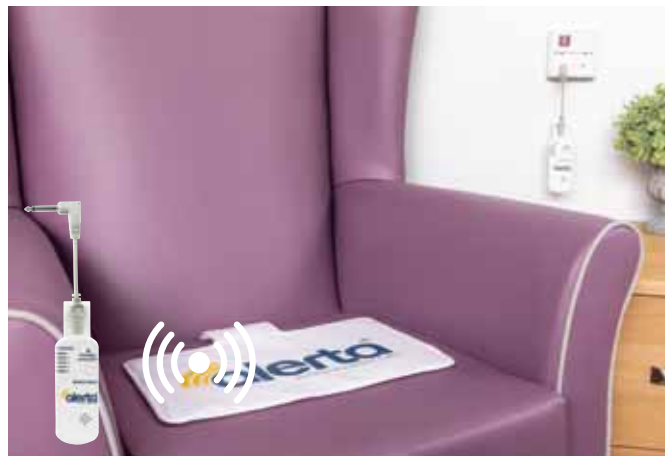

WIRELESS CHAIR ALERTAMAT

W-CAM-T - Wireless Chair Alertamat - comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > 95m wireless range from transmitter to receiver
- > No cables, extra safe and versatile
- > System includes: 1x Chair Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

Alerta Wireless Transmitter

- > Battery powered
- > Includes 2x 1.5v AAA batteries
- > 0, 2 or 4 second delay to reduce false alarms



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



The **Alerta Detect** uses infrared technology to activate the alarm when body movement is detected within the range of the infrared beam.

The **Alerta Detect** is most effective when mounted at the door or bedside; activating the alarm when a user leaves their safe zone and breaks the infrared beam. This alerts the caregiver to ascertain the user's safety. It is ideal for providing users with the freedom and safety of their own space while still proactively giving the caregiver advanced warning when a user leaves their safe zone.



WIRED ALERTA DETECT

ALT-DET - Wired Alerta Detect - connects to a nurse call system with a 2 metre cable provided for centralised activation.

Alerta Detect Motion Sensor

- > Plug-and-play system
- > Mains or battery powered
- > Low battery alert
- > In-built alarm
- > Adjustable alarm volume (silent, low, high)
- > Hold/snooze button
- > System includes: 1x Alerta Detect motion sensor and 1x 2 metre cable
- > Takes 3x 1.5v AA batteries



WATCH WIRELESS
PRODUCT VIDEO

WIRELESS ALERTA DETECT

ALT-DET - Wireless Alerta Detect - can be connected to a nurse call system through the Alerta Wall Point Receiver for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Wall Point Receiver only.

- > Wireless connect-and-play system
- > 95m wireless range from transmitter to receiver
- > No cables, extra versatile
- > System includes: 1x Alerta Detect motion sensor (can be used with Alerta Wall Point Receiver - sold separately)

Alerta Detect Motion Sensor

- > Mains or battery powered
- > Low battery alert
- > In-built alarm
- > Adjustable alarm volume (silent, low, high)
- > Hold/snooze button
- > Wireless signal indicator
- > Takes 3x 1.5v AA batteries



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



WIRELESS WINDOW AND DOOR SENSOR



W-WDS - Wireless Alerta Window and Door Sensor - can be mounted on a door & door frame, or window & window frame and then connected to a nurse call system through the Alerta Wall Point Receiver for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through the Alerta Wall Point Receiver only.


- > Wireless connect-and-play system
- > 95m wireless range from transmitter to receiver
- > No cables, extra versatile
- > System includes: 1x Alerta Window and Door Sensor (must be used with Alerta Wall Point Receiver - sold separately) and takes 2x 1.5 AAA batteries

WIRED ALERTA NURSE CALL BUTTON

PCL - Wired Alerta Nurse Call Button - provides the user with reassurance that help is at hand. With a large easy to press button it activates the alarm when a user presses the button. This alerts the caregiver to ascertain the user's need.

It can be plugged directly into the nurse call system or via the spare 'Y' cable that is provided with the wired Deluxe Alertamat, Wood Alertamat and Floor+ Alertamat.

- > 2 metre cable

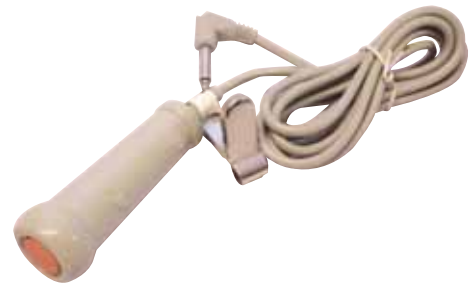

ALERTA WALL POINT RECEIVER

W-WPR - Alerta Wall Point Receiver - can be connected to a nurse call system for centralised activation, or used on its own, independent of a nurse call system to give localised alerts through this device only. The volume can be deactivated for silent in-room operation when connected to a nurse call system.

The **Alerta Wall Point Receiver** can be wirelessly paired with up to 5 Alerta wireless products and is compatible with the Deluxe Alertamat, Floor+ Alertamat, Bed Alertamat, Chair Alertamat, Alerta Detect and the Alerta Button, Alertmat Crash mat, WetSense Chair Alertamat, and WetSense Bed Alertamat.

It uses radio frequency technology to connect, with a wireless range of up to 95 metres. Once paired with Alerta wireless devices it only transmits to those devices and will not interfere with other Alerta Wall Point Receivers or devices in use nearby.

- > Wireless connect-and-play system
- > No cables, extra versatile
- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



WIRELESS ALERTA NURSE CALL BUTTON

W-NCB - Alerta Wireless Nurse Call Button is supplied with a wall mountable holder that can be mounted to a wall or bedside unit for holding the Wireless Nurse Call Button meaning it is always in one place if the user is bed bound, also making it easy for those affected by dementia to remember its location.

The **Alerta Wireless Nurse Call Button** is also supplied with a safety lanyard for use over the neck, and a wrist strap enabling it to be wrist worn. The button can be mounted in the wall holder with the wrist strap still attached.

- > Can be strapped round the wrist
- > Provided with a lanyard to put over the patients neck
- > Wall mountable
- > Can be worn when showering
- > Wristband material: Silicon
- > Wireless connect-and-play system
- > No cables, extra versatile
- > 95m wireless range from transmitter to receiver
- > Includes 1x CR2450 battery
- > System includes:
 - 1x Alerta Button (must be used with Alerta Wall Point Receiver - sold separately)



The **Alerta Wireless WetSense Bed and Chair Alertamats** are designed to inform caregivers when enuresis needs to be addressed. **Alerta WetSense** mats can be used under a sheet on a bed, or on a chair to wirelessly activate the alarm/nursecall system when wetness is detected. This alerts the caregiver to take appropriate action and ensure a comfortable resolution to the individual's needs.

The **Wireless WetSense Alertamats** can be connected to a nursecall system for centralised activation, or used on their own, independent of a nursecall system to give localised alerts.



WIRELESS WETSENSE CHAIR ALERTAMAT

W-CAM-WET-T - Wireless WetSense Chair Alertamat
- comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > Dimensions: 25 x 38 cm (10 x 15")
- > 95m wireless range from transmitter to receiver
- > System includes: 1x WetSense Chair Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

Alerta Wireless Transmitter

- > Battery powered
- > Includes 2x 1.5v AAA batteries
- > 0, 2 or 4 second delay to reduce false alarms



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



WIRELESS WETSENSE BED ALERTAMAT

W-BAM-WET-T - Wireless WetSense Bed Alertamat
- comes with an Alerta Wireless Transmitter which connects to the Alerta Wall Point Receiver.

- > Wireless connect-and-play system
- > No cables, extra safe and versatile
- > Dimensions: 25 x 76 cm (10 x 30")
- > 95m wireless range from transmitter to receiver
- > System includes: 1x WetSense Bed Alertamat and 1x Alerta Wireless Transmitter (must be used with Alerta Wall Point Receiver - sold separately)

Alerta Wireless Transmitter

- > Battery powered
- > Includes 2x 1.5v AAA batteries
- > 0, 2 or 4 second delay to reduce false alarms



Alerta Wall Point Receiver

- > High, low and silent volume control
- > Low battery alert
- > Mains or battery powered
- > Takes 2x 1.5v AA batteries



These products are designed as a monitor alerting the staff to user enuresis. They do not monitor a user leaving the bed or chair.

These mats have a 3 month warranty due to the nature of the substance they are expected to be in regular contact with. The WetSense Alertamats can be replaced on their own and the Alerta Wireless Transmitter and Alerta Wall Point Receiver be reused with a replacement mat.

Accessories & Plug Identification

W-FBCT

ALERTA WIRELESS TRANSMITTER

The Alerta Wireless Transmitter is supplied with the following wireless systems: Deluxe Alertamat, Floor+ Alertamat, Bed Alertamat, Chair Alertamat, Alertmat Crash mat, WetSense Chair Alertamat, and WetSense Bed Alertamat.

It uses radio frequency technology to connect, with a wireless range of up to 95 metres. Once paired with the Alerta Wall Point Receiver it will only transmit to that Wall Point Receiver and will not interfere with other Alerta Wall Point Receivers or devices in use nearby. Includes 2x 1.5v AAA batteries.

It incorporates an adjustable timer setting for alarm delay if required.

Low battery light and beep every 3 seconds.



AUAM

ALERTA ALARM MONITOR

The Alerta Alarm Monitor can be used with all wired Alertamats. It includes 1x 9v battery or can be mains operated.



MA

MAINS ADAPTOR

Works with the Alerta Alarm Monitor.



EXT

EXTENSION LEAD

Please specify the required length at point of order.



DPA

DOUBLE PLUG ADAPTOR

Allows two devices to be connected to one nurse call point.



Identify your plug type

Most nurse call systems take one of the following plug types:

SINGLE RING PLUG



DOUBLE RING PLUG



WHITE PLUG, 6 GROOVE



CLEAR PLUG, 6 GROOVE



CLEAR PLUG, 4 GROOVE



Please ensure that your clients know the plug type they require at point of order to ensure that their Alertamats will function correctly when they receive them.

Plug types can sometimes be difficult to identify. If your client is unsure or if you need assistance, please contact Alerta Medical experts on +44 (0)3452 088 097 or sales@alertamedical.com

SCAN HERE FOR
ALERTA MEDICAL
PLUG TYPE
SELECTOR





HIP PROTECTION

CLINICALLY PROVEN

SEAMLESS

WASHABLE

Hip protection helps prevent hip fractures

People who have suffered a fall – or who are at risk of falling – can benefit from effective hip protection. What starts out as an innocent fall may lead to an operation, a long stay in hospital and an increased demand for professional care. In addition, 30% of all people who suffer a hip fracture die within 2 years (1).

A growing number of elderly people indicates a future rise in the occurrence of hip fractures. Worldwide, it is estimated that the number of hip fractures will increase from 1.66 million into 2.6 million in 2025 and 6.26 million in 2050 (2).

Hip protection should be seen as an extra security that contributes to the user's own safety on the same level as seatbelt, lifejacket, or motorcycle helmet. Hip protection reduces the risk of serious injury, saves money and can also improve quality of life by promoting a sense of security, thus allowing the user to remain active, social and reduce the fear of falling.

SAFEHIP[®] AirX[®] is a 100% textile hip protector with the protector shield based on **AirX[®]** spacer fabric technology, which is a breathable, skin friendly and soft textile material. The **AirX[®]** technology is your guarantee for optimal compliance.

SAFEHIP[®] AirX[®] hip protectors come in three models – a female, male and unisex model in white. The pants are flexible and comfortable to wear, not too warm, and may be used with incontinence pads as a fixation pant. They have been designed to perform like regular underwear, with a fly opening on the male model and lace-look detail on the female. When the gender of the patient is not known in advance, the Unisex brief is suitable for both.

All SAFEHIP[®] AirX[®] pants are very durable and can be washed in industrial washing machines at up to 95°C/200°F and tumble dried at highest temperatures.

Unique features:

- > Patented horse-shoe shaped shield designed to disperse energy from a fall away from the hip joint
- > The shields adapt to the body perfectly, so the hip protector is not noticeable under clothes
- > Excellent breathability (test winner in EMPA test no. 448819)
- > Prevents moisture build-up and feels comfortable against skin
- > Minimises skin irritation such as reddening and itching
- > Pants are very comfortable to wear - no irritating seams
- > High percentage of elastane ensures a perfect fit, which also makes the pants suitable for incontinent users




SAFEHIP [®]	Cm	Inch	Find your size
S	75-95	30-38	
M	90-110	36-44	
L	100-120	40-48	
XL	110-140	44-56	

Sizes overlap. To find the user's correct size, measure the hips at the widest point. If the hip measurement lies between two sizes, choose the smaller size if the patient has slim thighs, and the larger size if thighs are wide or an incontinence pad will be worn.

Protector

100% AirX[®] soft textile

Seamless Briefs
58% Cotton
36% Polyamide
6% Elastane

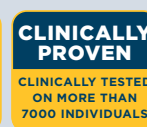

SAFEHIP AIR-X MALE



Male version with fly

ALT-SAM-S	Air-X Male - S
ALT-SAM-M	Air-X Male - M
ALT-SAM-L	Air-X Male - L
ALT-SAM-XL	Air-X Male - XL

SAFEHIP AIR-X FEMALE



Female version with flower detail at front

ALT-SAF-S	Air-X Female - S
ALT-SAF-M	Air-X Female - M
ALT-SAF-L	Air-X Female - L
ALT-SAF-XL	Air-X Female - XL

SAFEHIP AIR-X UNISEX



Neutral design for both sexes

ALT-SAU-S	Air-X Unisex - S
ALT-SAU-M	Air-X Unisex - M
ALT-SAU-L	Air-X Unisex - L
ALT-SAU-XL	Air-X Unisex - XL

SAFEHIP AIR-X OPEN



Open version

ALT-SAO-S	Air-X Open - S
ALT-SAO-M	Air-X Open - M
ALT-SAO-L	Air-X Open - L
ALT-SAO-XL	Air-X Open - XL

1) IOF - EFFE Jay magaziner, PhD, MSHyg, Eleanor M. Simonsick, PhD, T. Michael Kashner, PhD, J Richard Hebel, PhD, and John E. Konzona, MD Hvidovre Symposium Open Bridge 08/2000
2) Cooper et al. Osteoporosis Int. 1992; 2:285-9 & Melton L.J. Bone 1996; 18:121S-5S
3) Hindsø K, Lauritzen JB. Intervention study with hip protectors [abstract]. Osteoporosis International, 1998; 8:119.




Alternating Air Systems

Product Overview

		Overlay pad	
		Overlay system	
		Overlay system	
		Overlay system	
		Replacement system	
		Full Replacement system	
		Full Replacement system	
		Full Replacement system	
		Full Replacement system	
		Full Replacement system	
		Cushion system	
		Bariatric Cushion System	
		Battery cushion system	



Extension flaps

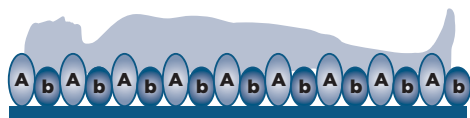

ALT-9001

ALERTA BUBBLE 2



with tuck-in extension flaps



Alerta Bubble 2 – a light weight and compact overlay alternating pressure relieving mattress system for effective prevention and treatment of users at low risk of developing a pressure ulcer in hospital, nursing and care home environments. **Alerta Bubble 2** has been designed to give care providers the basic essentials for pressure care and the pump is a basic entry level pump.

With simple to use settings and functionality, the **Alerta Bubble 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Bubble 2** has been manufactured to comply with the most stringent quality and in-use guidelines.


*If putting in washing machine, ensure the air inlet tubes are watertight sealed.

Mattress Specification

- > Mattress Dimensions:
200 x 90 x 7 cm (78.5 x 35.5 x 3")
- > Mattress Weight: 2.5kg (5.5lb)
- > Cell Structure: 130 Bubble cells
- > Material: Medical grade PVC
- > Max Load: 135kg (297.5lb)
- > Machine washable*

Pump Specification

- > Pump Dimensions:
24 x 11 x 8 cm (9.4 x 4.3 x 3.1")
- > Pump Weight: 0.9kg (1.7lb)
- > Cycle Time: 12 minutes
- > Power Input:
AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets

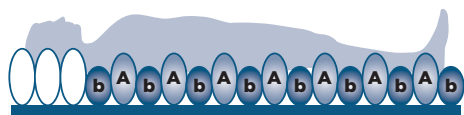

ALT-9004

ALERTA PEARL OVERLAY SYSTEM



ALT-9004 – Alerta Pearl – an overlay alternating pressure relieving mattress system with 5" air cells for effective prevention and treatment of users at medium to high risk of developing a pressure ulcer in hospital, nursing and care home environments. **Alerta Pearl** has been designed to give care providers the essentials for pressure care and the pump is an entry level pump.

With simple to use settings and functionality, the **Alerta Pearl** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Pearl** has been manufactured to comply with the most stringent quality and in-use guidelines.

Key Features

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode

Mattress Specification

- > Mattress Dimensions: 200 x 90 x 12.7 cm (78.5 x 35.5 x 5")
- > Mattress Weight: 5.5kg (12.1lb)
- > CPR: Quick release toggle
- > Max Load: 160kg (353lb)

Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification


- > Pump Dimensions: 28 x 14 x 9 cm (11 x 5.5 x 3.5")
- > Pump Weight: 1.4kg (3lb)
- > Cycle Time: 12 minutes
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure and power failure

Visual and audible warning of air loss and power failure




Weight setting control knob

Static/alternating modes

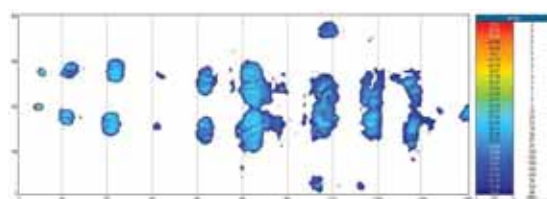

Clinical pressure mapping evaluation

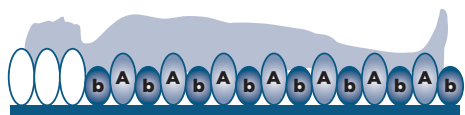

ALT-9005

ALERTA EMERALD 2 OVERLAY SYSTEM



ALT-9005 – Alerta Emerald 2 – an overlay alternating pressure relieving mattress system with 5" air cells for effective prevention and treatment of users at high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Emerald 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Emerald 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Alerta Emerald 2 incorporates cutting-edge **SensaCare TruFlo Technology®** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

Key Features

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode

Mattress Specification

- > Mattress Dimensions: 200 x 90 x 12.7 cm (78.5 x 35.5 x 5")
- > Mattress Weight: 5.5kg (12.1lb)
- > CPR: Quick release toggle
- > Max Load: 160kg (353lb)


Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

KEY FEATURES AND BENEFITS

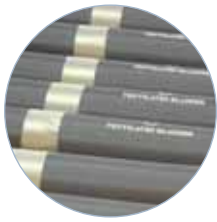

Nylon TPU coated cells




Independently removable cells



Covert sealed air cells



Branded cells

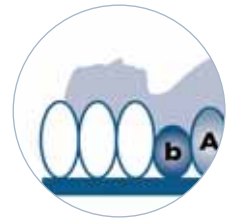

Pillow function




CPR release




Anchoring straps



Zipped storage bag

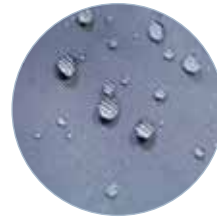

Multi-stretch, waterproof & vapour permeable PU cover

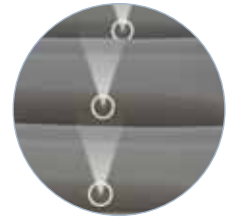

Micro-perforated cells



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode

CONTROL PANEL

Auto pressure sensing for consistent distribution of air


Visual and audible warning of air loss, power failure and alternate failure

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)




Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency

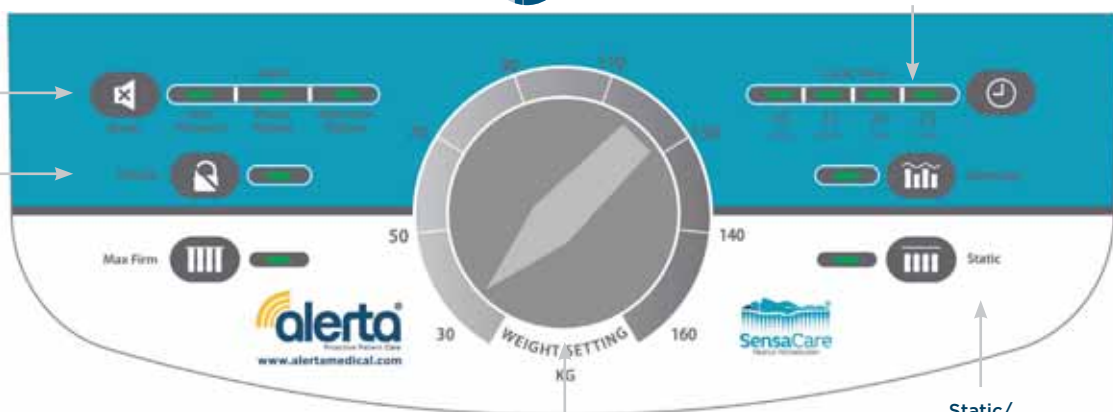


Weight setting control knob

Static/alternating modes

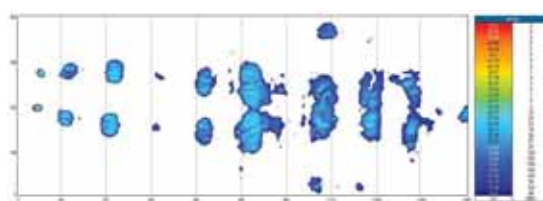

Clinical pressure mapping evaluation

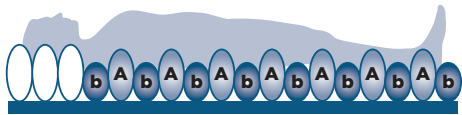

with remote

ALT-9006

ALERTA EMERALD AUTO OVERLAY SYSTEM



ALT-9006 - Alerta Emerald Auto – a fully **Auto Weight Sensing** overlay alternating pressure relieving mattress system with 5" air cells for effective prevention and treatment of users at high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Emerald Auto** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease. The pump has the latest cutting-edge **Auto Weight Sensing** technology, continuously adjusting interface pressures according to the user's weight.

The wide range of features make the **Alerta Emerald Auto** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Alerta Emerald Auto incorporates cutting-edge **SensaCare TruFlo Technology®** – constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

Key Features

- > Auto weight sensing
- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Robust pump connection
- > Transport mode

Mattress Specification

- > Mattress Dimensions: 200 x 90 x 12.7 cm (78.5 x 35.5 x 5")
- > Mattress Weight: 5.5kg (12.1lb)
- > CPR: Quick release toggle
- > Max Load: 160kg (353lb)


Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

KEY FEATURES AND BENEFITS

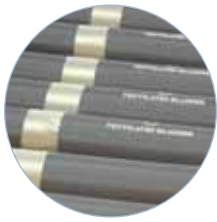

Nylon TPU coated cells




Independently removable cells



Covert sealed air cells



Branded cells

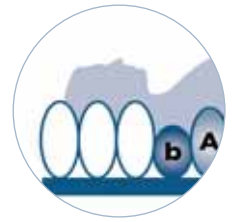

Pillow function




CPR release



Anchoring straps



Zipped storage bag

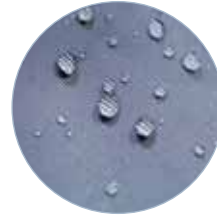

Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Auto weight



Remote control

CONTROL PANEL

Auto pressure sensing for consistent distribution of air



Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency

Auto weight control

Remote control communication


Visual and audible warning of air loss, power failure and alternate failure




Power on

Static/alternating modes

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)

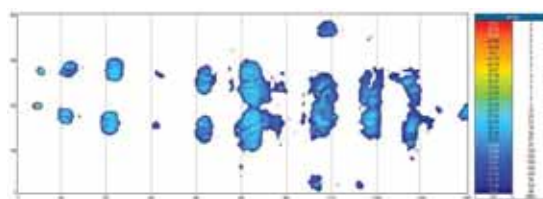

Clinical pressure mapping evaluation

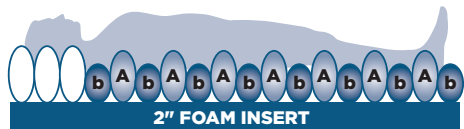

ALT-9007

ALERTA SAPPHIRE 2 REPLACEMENT SYSTEM



ALT-9007 – Alerta Sapphire 2 – a replacement alternating pressure relieving mattress system with 5" air cells and a self-inflating 2" sealed foam base forming an in-use height of 7" for effective prevention and treatment of users at high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Sapphire 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Sapphire 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Alerta Sapphire 2 incorporates cutting-edge **SensaCare TruFlo Technology®** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

Key Features

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel

Mattress Specification

- > Mattress Dimensions: 200 x 90 x 12.7 cm + 5cm foam base (78.5 x 35.5 x 5" + 2" foam base)
- > Mattress Weight: 5.5kg (12.1lb)
- > Cell Structure: Cell-on-foam
- > CPR: Quick release toggle
- > Max Load: 200kg (441lb)
- > Heel function (3 x heel cells with one way valves)


Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

KEY FEATURES AND BENEFITS

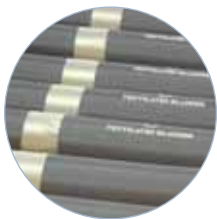

Nylon TPU coated cells



Independently removable cells



Covert sealed air cells



Branded cells



Pillow function



CPR release




Anchoring straps



Zipped storage bag



Multi-stretch, waterproof & vapour permeable PU cover

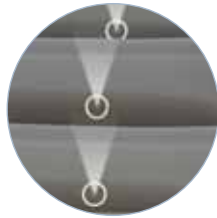

Micro-perforated cells



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Removable heel cells

CONTROL PANEL

Auto pressure sensing for consistent distribution of air




Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency


Visual and audible warning of air loss, power failure and alternate failure

Control panel lock to prevent unauthorised Or accidental tampering (hold for 2 sec)

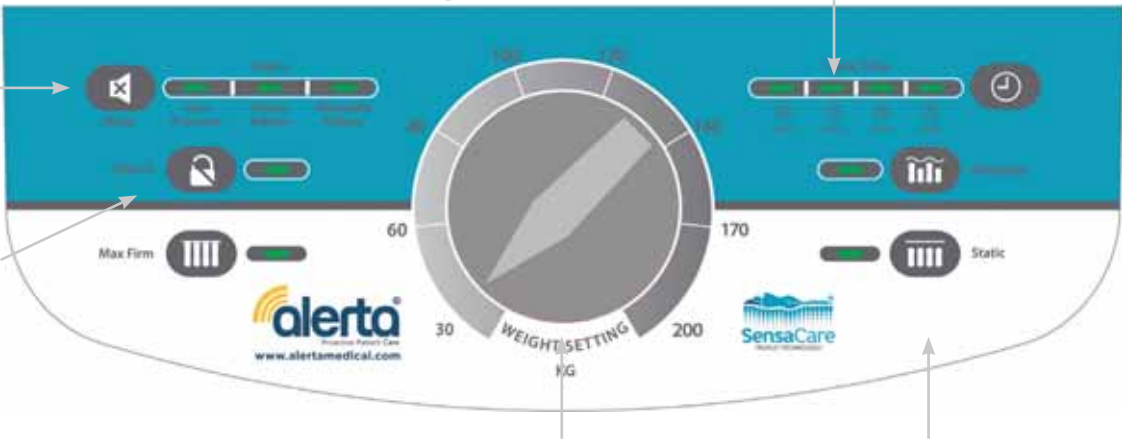


Weight setting control knob

Static/alternating modes

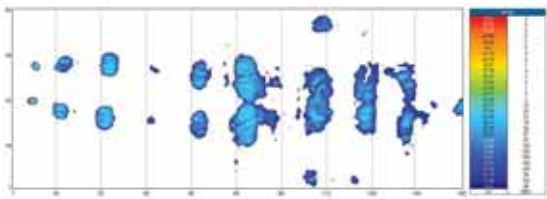

Clinical pressure mapping evaluation

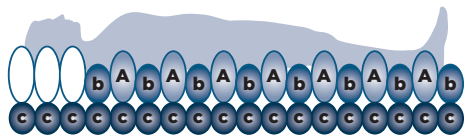

with remote

ALT-9008

ALERTA RUBY 2 REPLACEMENT SYSTEM



ALT-9008 – Alerta Ruby 2 – a replacement alternating pressure relieving mattress system with cell-on-cell air cells, a static base section and an in-use height of 8" for effective prevention and treatment of users at very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Ruby 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Ruby 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Alerta Ruby 2 incorporates cutting-edge **SensaCare TruFlo Technology®** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

Key Features

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Micro-perforated cells for air circulation
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel
- > Seat inflation feature for additional support when sitting upright
- > Digital user weight input

Mattress Specification

- > Mattress Dimensions: 200 x 90 x 20 cm (78.5 x 35.5 x 8")
- > Mattress Weight: 10kg (22lb)
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 200kg (441lb)
- > Heel function (3 x heel cells with one way valves)


Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

KEY FEATURES AND BENEFITS

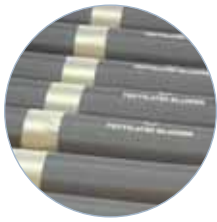

Nylon TPU coated cells



Independently removable cells



Removable heel cells



Branded cells



Pillow function




CPR release




Anchoring straps



Zipped storage bag

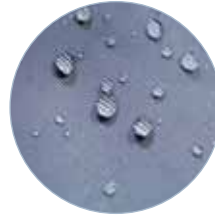

Multi-stretch, waterproof & vapour permeable PU cover

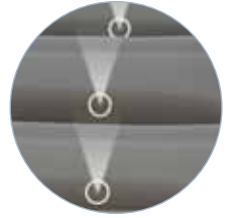

Micro-perforated cells



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Remote control

CONTROL PANEL



Digital LED display for easy use

Digital user weight input


Multiple cycle time control for efficiency



Ultra-Silent™ pump operation to avoid disturbance

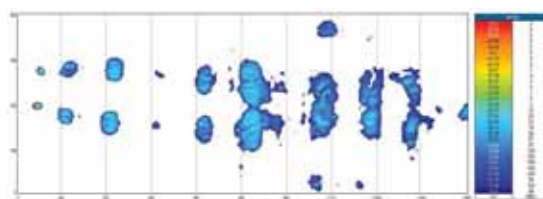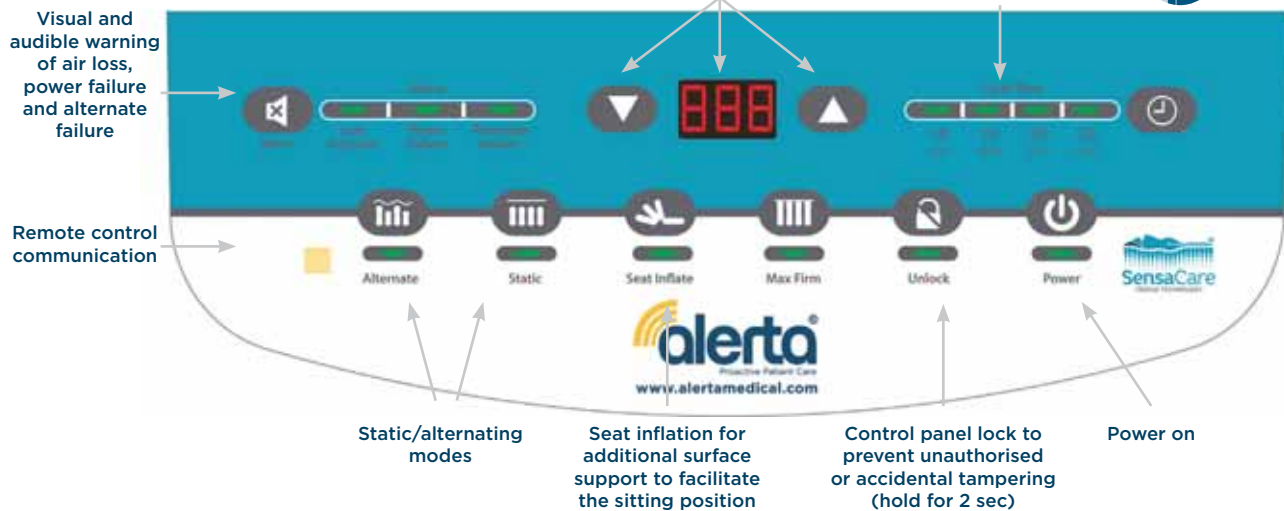

Clinical pressure mapping evaluation

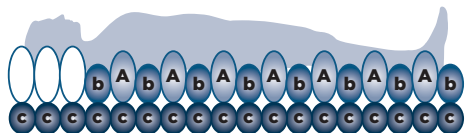

with remote

ALT-9009

ALERTA RUBY AUTO REPLACEMENT SYSTEM



ALT-9009 - Alerta Ruby Auto – a fully **Auto Weight Sensing** replacement alternating pressure relieving mattress system with cell-on-cell air cells, a static base section and an in-use height of 8" for effective prevention and treatment of users at very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Ruby Auto** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease. The pump has the latest cutting-edge **Auto Weight Sensing** technology, continuously adjusting interface pressures according to the user's weight.

The wide range of features make the **Alerta Ruby Auto** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Alerta Ruby Auto incorporates cutting-edge **SensaCare TruFlo Technology®** – constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

Key Features

- > Auto weight sensing
- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel
- > Seat inflation feature for additional support when sitting upright

Mattress Specification

- > Mattress Dimensions: 200 x 90 x 20 cm (78.5 x 35.5 x 8")
- > Mattress Weight: 10kg (22lb)
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 200kg (441lb)
- > Heel function (3 x heel cells with one way valves)


Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Operating System: SensaCare TruFlo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

KEY FEATURES AND BENEFITS

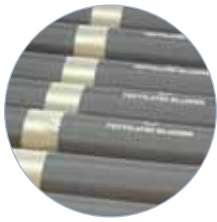

Nylon TPU coated cells



Independently removable cells



Removable heel cells



Branded cells



Pillow function




CPR release



Anchoring straps



Zipped storage bag

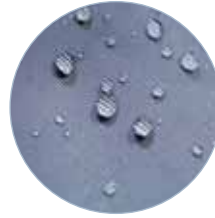

Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



Auto weight



Remote control

CONTROL PANEL

Digital LED display for easy use




Ultra-Silent™ pump operation to avoid disturbance

Multiple cycle time control for efficiency

Auto weight control

Remote control communication

Visual and audible warning of air loss, power failure and alternate failure

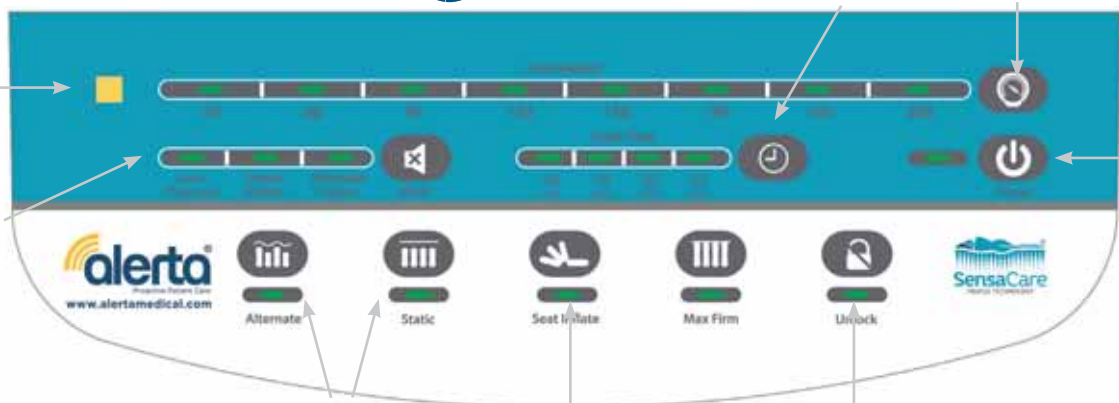


Power on


Static/alternating modes

Seat inflation for additional surface support to facilitate the sitting position

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)

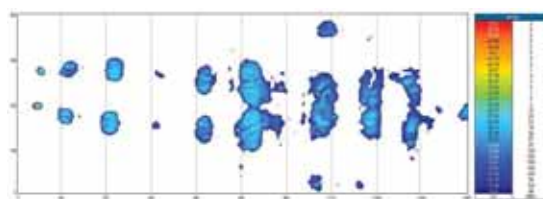

Clinical pressure mapping evaluation

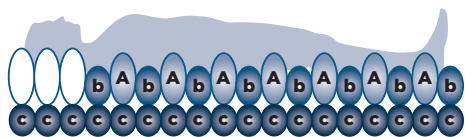

SCAN TO WATCH
PRODUCT VIDEO



ALT-9450

ALERTA RUBY TURN LATERAL TURNING REPLACEMENT ALTERNATING MATTRESS SYSTEM



ALT-9450 - Alerta Ruby Turn – a lateral turning full replacement alternating pressure relieving mattress system with cell-on-cell air cells, a static base section and an in-use height of 8" for effective prevention and treatment of users at very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Ruby Turn** uses automated lateral turning programs, maximising nursing time, allowing staff to focus on other critical duties with confidence. Built with 20 deep air cells in a cell-on-cell design, this mattress features a sturdy static base, supporting up to 230kgs. The programmable lateral turning system optimizes individual patient care by allowing caregivers to select left, right, or dual turning modes, featuring precise angle control adjustable to clinical requirements using 20°, 25° or 35°.

The wide range of features make the **Alerta Ruby Turn** a highly versatile and a cost-effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Alerta Ruby Turn incorporates cutting-edge **SensaCare Truflo Technology®** – constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with **Ultra-Silent™** pump operation.

Key Features

- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel
- > Seat inflation feature for additional support when sitting upright

Mattress Specification

- > Mattress Dimensions: 200 x 90 x 20cm (78.5 x 35.5 x 8")
- > Mattress Weight: 8.3kg (18.3lb)
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 230kg (500lb)
- > Heel function (3 x heel cells with one-way valves)


Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Operating System: SensaCare Truflo Technology®
- > Pump Dimensions: 37.8 x 17 x 26.6cm (14.8 x 6.7 x 10.5")
- > Pump Weight: 6.4kg (14lb)
- > Alternating cycle time: 10, 15, 20, 25mins
- > Lateral turning cycle time: 10, 20, 30, 60mins
- > Adjustable turning angles
- > Memory function
- > Air Output: 42-50 litres/minute
- > Power Input: AC220-240V, 50-60Hz 2A
- > Power Consumption (when turning): 18W
- > Power Consumption (daily usage): 10W
- > Alternating and static modes
- > Seat inflate
- > Ultra-Silent™ pump operation
- > Visual and audible alarm for low pressure, power and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

KEY FEATURES AND BENEFITS

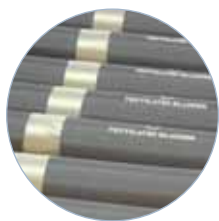


Nylon TPU coated cells



Independently removable cells



Removable heel cells

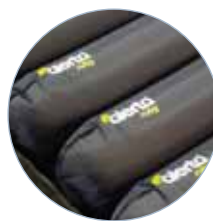

Branded cells



Pillow function



CPR release



Anchoring straps



Zipped storage bag



Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable




360 Degree zipper



Robust pump connection



Transport mode

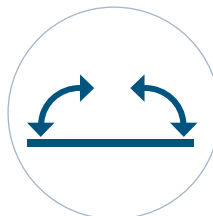

Lateral turning



Adjustable turn angles

CONTROL PANEL

Digital LED display for easy use




Ultra-Silent™ pump operation to avoid disturbance

Visual and audible warning of air loss, power failure and alternate failure

Adjustable cycle time control for efficiency

Adjustable rotation mode

Adjustable angle control

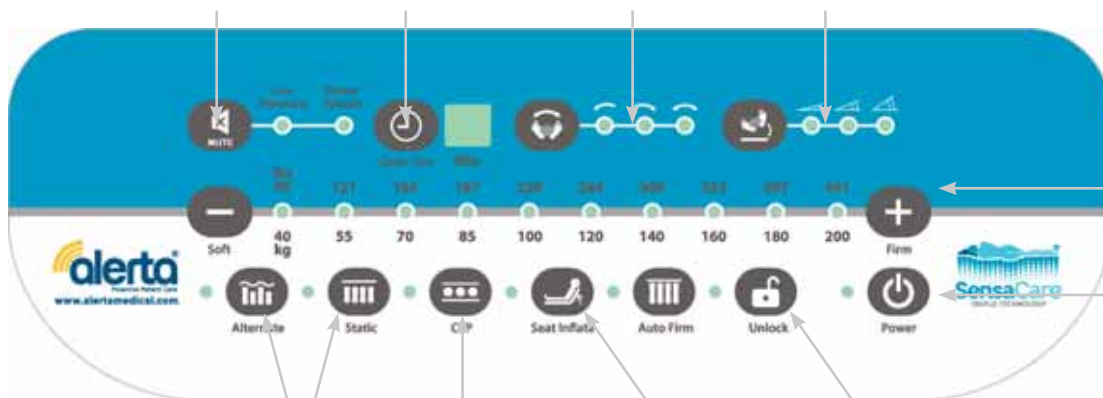

Auto weight control


Power on/off

Static/alternating modes


Press to set the system to static therapy mode at minimum pressure value

Seat inflation for additional surface support to facilitate the sitting position

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)

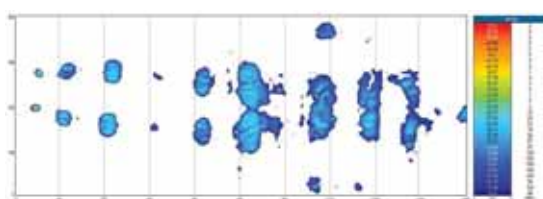

Clinical pressure mapping evaluation

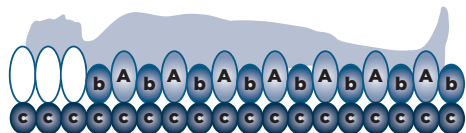

with remote

ALT-9300

ALERTA BARIATRIC 2 REPLACEMENT SYSTEM



ALT-9300 – Alerta Bariatric 2 – a replacement alternating pressure relieving mattress system with cell-on-cell air cells, a static base section and an in-use height of 10" for effective prevention and treatment of bariatric users up to very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Bariatric 2** enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the **Alerta Bariatric 2** a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Alerta Bariatric 2 incorporates cutting-edge **SensaCare Truflo Technology®** - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.

Key Features

- > 350kg weight capacity
- > Pillow function
- > CPR release
- > Anchoring straps for securing to bed
- > Robust pump connection
- > Transport mode
- > Individually removable heel cells to create zero pressure under the heel
- > Seat inflation feature for additional support when sitting upright
- > Digital user weight input

Mattress Specification

- > Mattress Dimensions: 200 x 120 x 25 cm (78.5 x 47 x 10")
- > Mattress Weight: 12.5kg (27.6lb)
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 350kg (772lb)
- > Heel function (3 x heel cells with one way valves)


Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Operating System: SensaCare Truflo Technology®
- > Pump Dimensions: 34.5 x 24 x 12.5 cm (13.5 x 9.4 x 4.9")
- > Pump Weight: 2.8kg (6.1lb)
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 12-13 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 9W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control panel lock feature to avoid accidental change of settings
- > Max firm mode

KEY FEATURES AND BENEFITS

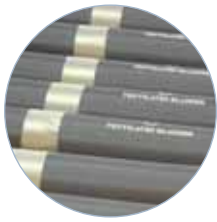

Nylon TPU coated cells



Independently removable cells



Covert sealed air cells



Removable heel cells



Pillow function



CPR release



Anchoring straps



Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable



360 Degree zipper



Robust pump connection



Transport mode



350KG Capacity



Remote control

CONTROL PANEL

Digital LED display for easy use

Digital user weight input

Multiple cycle time control for efficiency



Ultra-Silent™ pump operation to avoid disturbance

Visual and audible warning of air loss, power failure and alternate failure

Remote control communication




Static/alternating modes

Seat inflation for additional surface support to facilitate the sitting position

Control panel lock to prevent unauthorised or accidental tampering (hold for 2 sec)

Power on

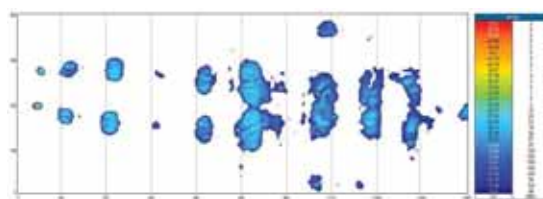


Clinical pressure mapping evaluation



BARIATRIC 120CM VERSION

Available in 90cm
and 120cm

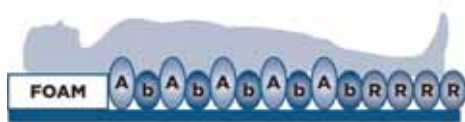

ALT-HYB2

ALERTA SENSAFLO
HYBRID MATTRESS ONLY

ALT-HYB2/4

BARIATRIC 4FT WIDE
ALERTA SENSAFLO HYBRID
MATTRESS ONLY

ALT-HYB2/05

ALERTA SENSAFLO HYBRID
MATTRESS WITH PUMP

ALT-HYB2/05/4

BARIATRIC 4FT WIDE
ALERTA SENSAFLO HYBRID
MATTRESS WITH PUMP



The **Alerta Sensaflo Hybrid** Mattress is the world's first 4 in 1 hybrid mattress, including a built-in fire evacuation system. It benefits from foam-filled air cells, SensaGel adaptive foot cells and an in-use height of 6" for effective prevention and treatment of users up to very high risk of developing a pressure ulcer in hospital, nursing, and care home environments.

When used without a pump the **Alerta Sensaflo Hybrid** performs the same as a high specification foam mattress and when used with a pump it performs the same as an alternating air system. There is no requirement to remove the user from bed to change from a foam system to an alternating air system.

The **Alerta Sensaflo Hybrid** has a built-in fire evacuation system with straps and handles, allowing the mattress to be quickly moved with the user in position during emergency situations. There are four adaptable gel-foam heel zone cells that can be individually removed and repositioned to create zero pressure under the heel.

By inflating the mattress and sealing off with the transport cap you create a balanced-air self-adjusting semi-active mattress. The air cells are contained within a high-density foam perimeter which creates a smooth transition from the side rail to the air cells.


The **Alerta Sensaflo Hybrid** incorporates cutting-edge SensaCare TruFlo Technology® – constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with Ultra-Silent™ pump operation.



SCAN TO WATCH
PRODUCT VIDEO

Ergonomically castellated
pillow function

In the centre of the mattress are 10
foam filled air cells which the user
interfaces with directly

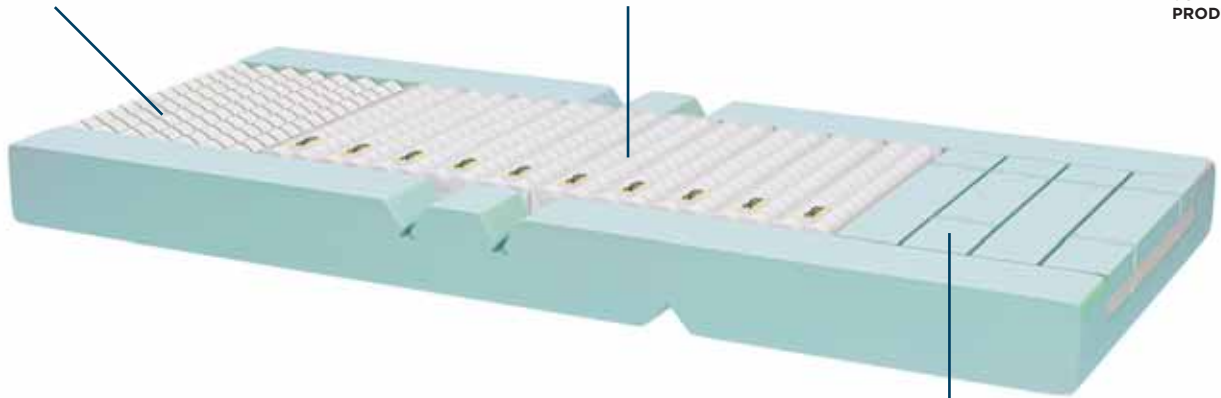

Four adaptable gel-foam heel zone
cells can be individually removed and
repositioned to create zero pressure
under the heel



Built-in fire evacuation system with
straps and handles, allow the mattress
to be quickly moved with the user in
position during emergency situations

Key Features

- > Pillow function
- > CPR release
- > Robust pump connection
- > Transport mode
- > Profiling side sections
- > Adaptable gel-foam heel zone
- > Built-in evacuation system

Mattress Specification

- > For users up to very high risk of developing a pressure ulcer
- > CPR quick release toggle
- > Max load: 250kg
- > Min load: 30kg
- > Mattress depth: 150mm (6")

Cover Specification

- > Multi-stretch and vapour permeable PU cover
- > Waterproof and machine washable cover
- > Antimicrobial cover with white underside
- > 360 degree zipper
- > BS7177: Crib 5
- > Mains cable management system

Pump Specification

- > Alternating and static modes
- > Visual and audible alarm for low pressure, power failure and alternate failure
- > Control Panel lock feature to avoid accidental change of settings
- > Fitted with hanging brackets
- > Max firm mode



ALT-HYB-C-SYS

ALERTA SENSAFLO HYBRID CUSHION WITH PUMP



ALT-HYB-C-SYS - Alerta Sensaflo Hybrid cushion system – an alternating foam-filled air cell pressure relieving cushion system for effective prevention and treatment of users up to very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Sensaflo Hybrid** cushion enables the care provider to quickly set up the cushion system and have it operating at optimal pressure for a specific user with ease. There is no requirement to remove the user from a chair to change from a foam cushion to an alternating air system.

Alerta Sensaflo Hybrid cushion can be used without a pump and boasts a unique SensaGel gel-infused memory foam topper for unrivalled pressure redistribution and improved circulation. When used with a pump it performs the same as an alternating air system. By inflating the cushion and sealing off with the transport cap you can create a balanced-air self-adjusting semi-active cushion.

The **Alerta Sensaflo Hybrid** incorporates cutting-edge **SensaCare TruFlo Technology®**; with constant pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure, and ultra-silent pump operation to promote an even ambience for natural comfort and stable blood circulation.

Cushion Specification


- > Cushion Dimensions:
45x45x12cm inflated including SensaGel insert
45x45x10cm uninflated
45x45x8cm uninflated & SensaGel insert removed
- > Cushion Weight: 1.5kg (3.3lb)
- > Cell Structure: Foam filled cell
- > Max Load: 160kg (353lb)

Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5

Pump Specification

- > Pump Dimensions:
27 x 15.5 x 10 cm (10.5 x 6.1 x 3.9")
- > Pump Weight: 1.6kg (3.5lb)
- > Cycle Time: 10 minutes
- > Power Input:
AC220-240V, 50-60Hz
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure

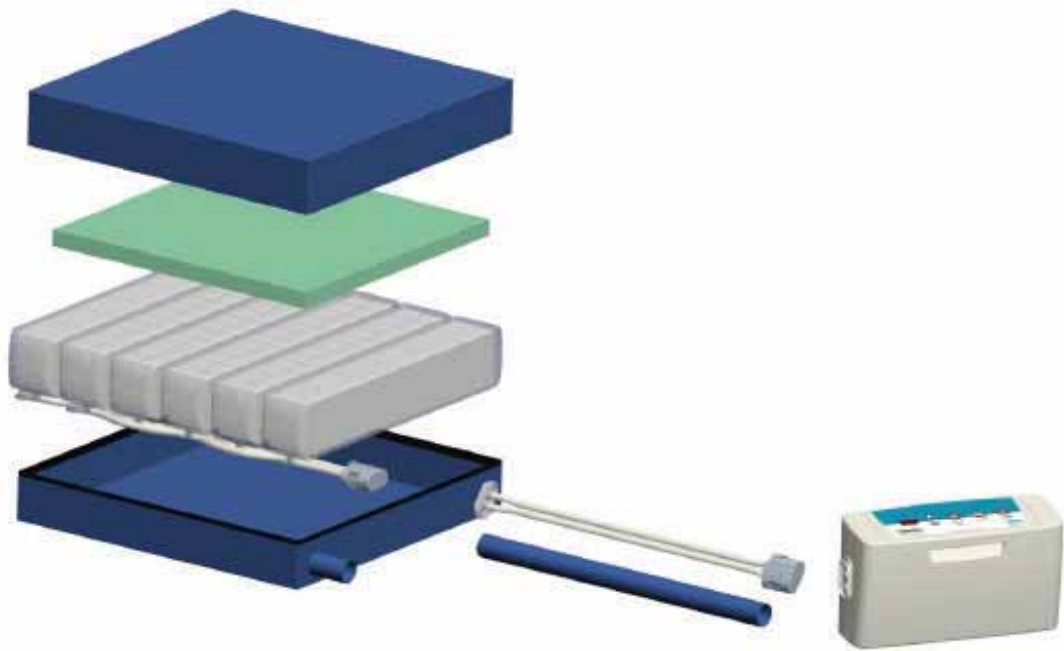


CONTROL PANEL

Digital LED display for easy use

Digital user weight input




Ultra-Silent™ pump operation to avoid disturbance

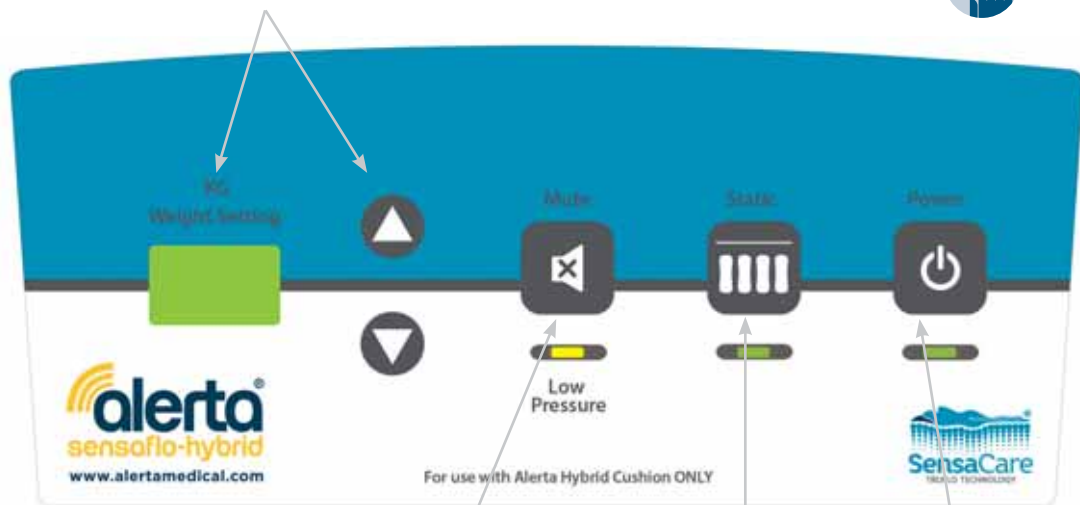

Visual and audio warning of air loss

Static mode (on when LED lit)

Power on

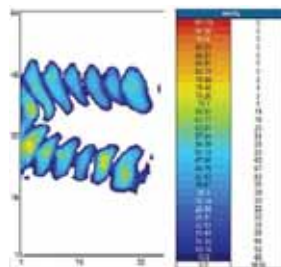

Clinical pressure mapping evaluation



ALT-213/02 ALERTA PARTNER CUSHION SYSTEM

ALT-213/02 - Alerta Partner – an alternating pressure relieving cushion system with a sealed 2" foam base for effective prevention and treatment of users up to high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the **Alerta Partner** enables the care provider to quickly set up the cushion system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Partner** is a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Cushion Specification

- > Cushion Dimensions: 45 x 50 x 10 cm (17.5 x 19.5 x 4")
- > Cushion Weight: 1.4kg (3.1lb)
- > Cell Structure: Cell-on-foam
- > Max Load: 150kg (331lb)

Cover Specification


- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover
- > BS7177: Crib 5

Pump Specification

- > Pump Dimensions: 27 x 13.5 x 9 cm (10.5 x 5 x 3.5")
- > Pump Weight: 1.6kg (3.5lb)
- > Cycle Time: 12 minutes
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets
- > Alternating and static modes
- > Visual and audible alarm for low pressure and power failure

CONTROL PANEL

Visual and audible warning of air loss and power failure

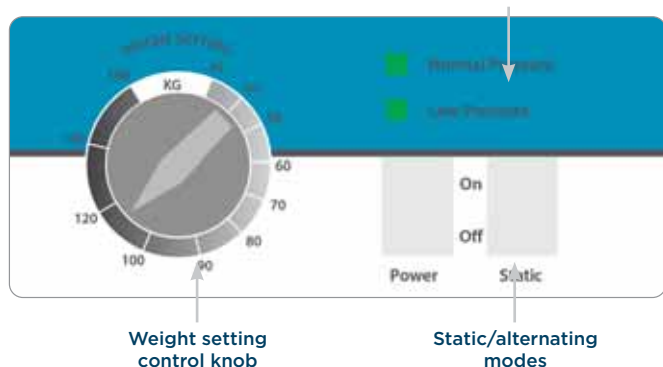

KEY FEATURES AND BENEFITS



Nylon TPU coated cells



Independently removable cells



Multi-stretch, waterproof & vapour permeable PU cover



ALT-214/02B - Alerta Partner Plus – an alternating pressure relieving cushion system with a sealed 2" foam base for effective prevention and treatment of users up to very high risk of developing a pressure ulcer in hospital, nursing and care home environments.

Cover Specification

- > Cushion Dimensions: 63 x 51 x 12cm (24.8 X 20 X 4. 7")
- > Max user weight: 250kg (551lb)





ALT-215/04

ALERTA MOBILE BATTERY CUSHION SYSTEM



ALT-215/04 – Alerta Mobile – a battery powered alternating pressure relieving cushion system for effective prevention and treatment of users up to high risk of developing a pressure ulcer in hospital, nursing and care home environments. The **Alerta Mobile** control unit is equipped with a long lasting Lithium Ion battery. In static mode the battery gives a remote operation time of more than 24 hours and in alternating mode the system can last at least one day.

With simple to use settings and functionality, the **Alerta Mobile** enables the care provider to quickly set up the cushion system and have it operating at optimal pressure for a specific user with ease.

The **Alerta Mobile** is a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Cushion Specification

- > Cushion Dimensions: 45 x 50 x 10 cm (17.5 x 19.5 x 4")
- > Cushion Weight: 2.4kg (5.3lb)
- > Cell Structure: Cell-on-foam
- > Max Load: 130kg (286.6lb)

Cover Specification



- > Multi-stretch PU
- > Vapour permeable
- > Waterproof
- > Machine washable
- > Antimicrobial
- > Enclosed zipper
- > White underside of cover
- > BS7177: Crib 5

Pump Specification


- > Pump Dimensions: 21 x 10 x 6 cm (8.2 x 4 x 2.5")
- > Pump Weight: 0.7Kg (1.54lb)
- > Cycle Time: 10, 15, 20 minutes
- > Power Input: AC220-240V, 50-60Hz
- > Battery: 12V Lithium 2200mAh
- > Enclosed pump case fitted with hanging straps
- > Alternating and static modes
- > Visual and audible alarm for low pressure and low battery power

CONTROL PANEL

KEY FEATURES AND BENEFITS

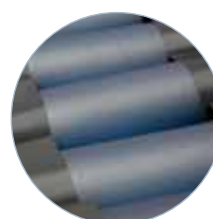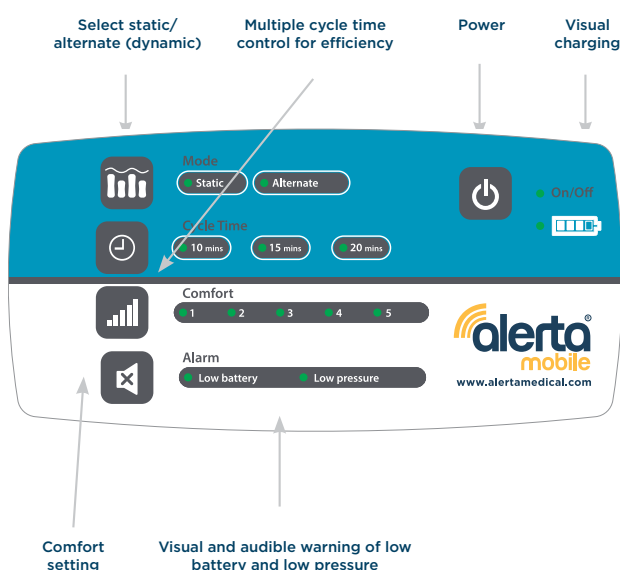

TPU coated cells

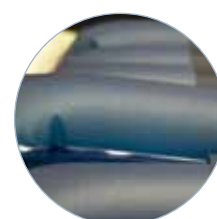

Independently removable cells



Multi-stretch, waterproof & vapour permeable PU cover



Antimicrobial, machine washable



360 Degree zipper

Product Selector

					
					
User Risk Level	High/Very High	Very High	Med/High	High	High
Max Load	250kg	160kg	160kg	160kg	160kg
Cell structure	Foam and air/foam	Foam and air/foam	Cell-on-cell	Cell-on-cell	Cell-on-cell
Mattress/cushion depth	15cm	12cm	12.7cm	12.7cm	12.7cm
Overlay system			✓	✓	✓
Replacement system	✓	✓			
Auto weight sensing					✓
Built-in evacuation system	✓				
Nylon TPU coated cells			✓	✓	✓
Independently removable cells	✓	✓	✓	✓	✓
Heel zone	✓				
Covert sealed air cells			✓	✓	✓
Pillow function	✓		✓	✓	✓
CPR	✓		✓	✓	✓
Anchoring straps			✓	✓	✓
Zipped storage bag			✓	✓	✓
Multi-stretch PU cover	✓	✓	✓	✓	✓
Waterproof & vapour permeable cover	✓	✓	✓	✓	✓
Antimicrobial & washable cover	✓	✓	✓	✓	✓
360 degree zipper	✓	✓	✓	✓	✓
Robust pump connection	✓	✓	✓	✓	✓
Transport mode	✓	✓	✓	✓	✓
Remote control					✓
Clinical pressure mapping evaluation	✓	✓	✓	✓	✓

alerta
sapphire 2

alerta
ruby 2

alerta
ruby-auto

alerta
ruby turn

alerta
bariatric 2

alerta
partner

alerta
partner plus

alerta
mobile



High	Very High	Very High	Very High	Very High	High	High	High
200kg	200kg	200kg	230kg	350kg	150kg	250kg	130kg
Cell-on-cell -on-foam	Cell-on-cell	Cell-on-cell	Cell-on-cell	Cell-on-cell	Cell-on-cell -on-foam	Cell-on-cell -on-foam	Cell-on-cell -on-foam
17.7cm	20cm	20cm	20cm	25cm	10cm	12cm	10cm
✓	✓	✓ ✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓
✓	✓	✓	✓	✓	✓	✓	✓

Suction Machines

AS1003A

3A PROFESSIONAL ASPIRATOR SINGLE PUMP, 1000ML BOTTLE



- > 16 litre/minute flow rate
- > Adjustable vacuum level
- > Machine only weighs 2.3kg
- > Carry handle

What's included?

- > Machine
- > 1000ml jar
- > Suction tubing
- > 1 suction catheter
- > Spare anti-bacterial filter
- > Mains power lead
- > Connector/flow regulator
- > Disposable jar liner

AS903A

3A PROFESSIONAL ASPIRATOR DOUBLE PUMP, 1000ML BOTTLE



- > 22 litre/minute flow rate
- > Adjustable vacuum level
- > Machine only weighs 3.2kg
- > Carry handle

What's included?

- > Machine
- > 1000ml jar
- > Suction tubing
- > 1 suction catheter
- > Spare anti-bacterial filter
- > Mains power lead
- > Connector/flow regulator
- > Disposable jar liner

AS457A

3A MINI RECHARGEABLE PORTABLE ASPIRATOR



- > 30 litre/minute flow rate
- > Adjustable vacuum level
- > Machine only weighs 4.6kg
- > Carry handle
- > Rechargeable for independent use
- > Can also be mains operated
- > 12V connection facility

What's included?

- > Machine
- > 1000ml jar
- > Suction tubing
- > 1 suction catheter
- > Spare anti-bacterial filter
- > Mains power lead
- > Connector/flow regulator
- > Disposable jar liner
- > Shoulder bag/carry case

Nebulisers

AP508A

**HAPPYNEB III
NEBULISER**



- > 14 litre/minute flow rate
- > .45ml/minute nebulisation rate
- > Machine weighs 2kg
- > Mains operated

What's included?

- > Nebuliser
- > Power cable
- > Tubing
- > 2 face masks
- > Chamber
- > Mouth piece
- > Nose piece
- > Filter

Accessories

3A1385

**ANTIBACTERIAL FILTER
FOR 3A PROFESSIONAL
ASPIRATORS**



3A561

**SILICON SUCTION
TUBING 1.3M LONG**



3A560

**MANUAL FLOW
REGULATOR
& CONNECTOR**



3A1055

**MULTI VOLTAGE
CHARGER UNIT FOR
AS457A BATTERY
PORTABLE ASPIRATOR**



3AC285

**REPLACEMENT
VASE 1000ML**



3A1687

**DISPOSABLE SUCTION
LINER**



3AC170

**FACE MASK KIT
FOR NEBULISER**

What's included?

- > Face mask
- > Air tube
- > Chamber



Diagnostics

ALT-DE100

**DIGITAL NON-CONTACT
FOREHEAD & EAR
THERMOMETER**

An easy-to-use functional thermometer for recording body temperature quickly and accurately. It holds a memory of the last 30 readings and automatically shuts down after 15 seconds of no use.



ALT-DE102

**DIGITAL BODY
THERMOMETER**

An easy-to-use functional thermometer for recording body temperature quickly and accurately. It holds a memory of the last reading and automatically shuts down after 10 minutes of no use.



ALT-DE200

UPPER ARM DIGITAL BLOOD PRESSURE MONITOR

An easy-to-use functional device for recording blood pressure quickly and accurately. It measures systolic pressure, diastolic pressure and pulse. With automatic inflation and deflation, this device is comfortable to wear and provides a clinically precise reading in seconds.

It has a memory capacity for 2 sets of 90 readings, provides an average of the latest 3 readings and automatically shuts down after no use.

A cuff is provided and fits upper arm circumferences of 22-32cm.




ALT-DE300

**FINGERTIP PULSE
OXIMETER**

An easy-to-use functional oximeter for recording blood oxygen saturation and pulse rate quickly and accurately by simply placing on the finger. This device is comfortable to use and provides a clinically precise reading in seconds. It automatically turns off after no use and has a low battery indicator.





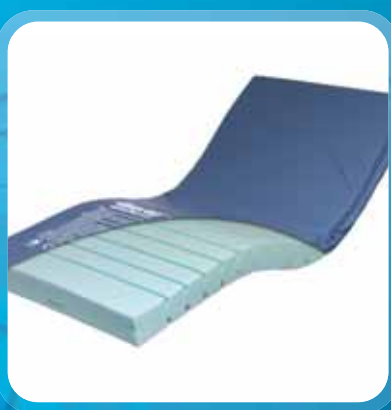

BEDS & ACCESSORIES | OVERBED TABLES | BED RAIL BUMPERS | PATIENT POSITIONING | EVACUATION
LAUNDRY & WASHABLES | PATIENT TRANSFERS | WHEELCHAIRS | WALKING AIDS | SHOWERING &
TOILETING | FOAM PRESSURE CARE | GEL PRESSURE CARE | FALL PREVENTION | HIP PROTECTION
ALTERNATING PRESSURE CARE | HYBRID PRESSURE CARE | SUCTION & NEBULISATION