



2025

Product Guide

26TH
Anniversary
Celebration
SINCE 1999



**your health,
we care.**

ISO13485 FDA CE MDR MDMA

A Leading Manufacturer Of Home Medical Equipment

XIAMEN SENYANG CO., LTD

ABOUT US

Xiamen Senyang Co., Ltd was found by Mrs. Susan Wang in 1999, with the intention to make more people healthier.

Senyang is a professional alternating pressure air mattress system for pressure ulcer prevention developer and manufacturer. After got successful with products for Alternating pressure mattress system in shortly time, Senyang took also compression therapy system and transfer lift into the portfolio.

Today Senyang presents itself as a leading manufacturer in home medical equipment.

10

Production lines

24,000

Square meters

200+

Employees

5,000

Air mattresses per day

FDA



ISO



CE





What can we offer to you?

Focus on optimizing the production of home care products for the elderly disabled. Senyang continuous development in R&D technology , and runs by a core team of highly qualified professionals with a deep understanding of the products and services. It helps Senyang to satisfy the needs of customers from developing and designing to manufacture requests. Senyang always try to convey our services with the best care, knowledge and a passion for excellence.

YOUR HEALTH, WE CARE.



CONTENTS

01

COMPANY PROFILE

04

MANUFACTURING ENVIRONMENT

05

CORPORATE PROMOTION

06

SYSTEM COMBINATIONS

07

PUMP FUNCTIONS OVERVIEW

08-40

PRODUCT INTRODUCTION

08-17	ALTERNATING PRESSURE MATTRESS WITH PUMP
18	ANTI-DECUBITUS FOAM MATTRESS
19	TRANSFER MATTRESS
20-24	PRESSURE SEAT CUSHION
25-30	AIR PRESSURE THERAPY SYSTEM
31-39	NURSING TRANSFER LIFT
40	MEDICAL AUXILIARY PRODUCT

MANUFACTURING ENVIRONMENT

Experienced

We have over 26 years of industry production experience and are committed to providing professional elderly care and rehabilitation products and solutions.

Excellent quality control

We have over 200 professional skilled workers, advanced production equipment, 24000 square meters of production area, and adopt strict product quality control processes.

Our Qualifications: ISO13485 : 2016.

Our principle: Quality, Innovation, Honesty, Development.

Continuous and efficient R & D

We have a professional R&D team that continuously develops products and technologies to ensure that we remain at the forefront of the industry.

Automated production material processing flow



Air mattress production equipment




Mattress performance test




Air pump test



3D printing equipment

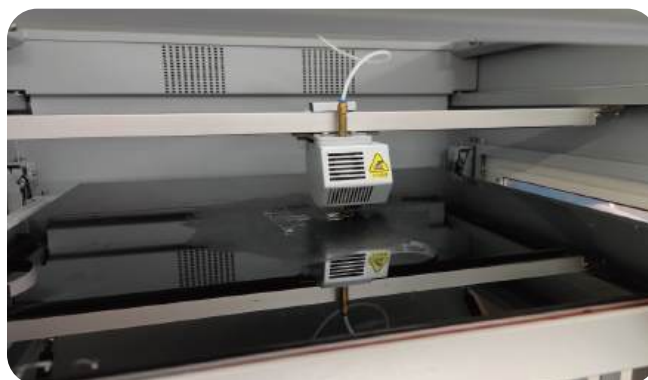


CORPORATE PROMOTION

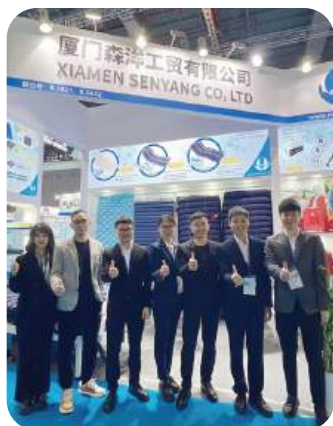

CMEF



REHACARE

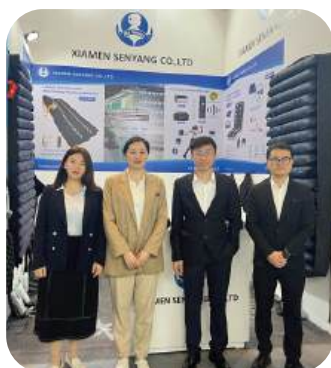

ARAB HEALTH



MEDICA



FIME



Canton fair

Customized services

Our company undertakes OEM & ODM business, and this catalog does not display products customized by all global brands.

We can meet customers' personalized customization requirements and even develop new products based on their needs.

Exhibition time

We maintain active contact and communication with customers through professional exhibitions held around the world every year, and pass on our manufacturing capabilities and professionalism to more people. We have also made many new friends, laying a solid foundation for more extensive cooperation in the future.

Welcome to visit Senyang at professional exhibitions around the world every year.



ALTERNATING PRESSURE MATTRESS WITH PUMP SYSTEM

STAGE I-IV FOR ANTI DECUBITUS & PRESSURE ULCER

As the color deepens, it means the risk of bedsores increases



SYSTEM COMBINATIONS	STAGE I	STAGE II	STAGE III	STAGE IV
	PVC	Nylon PVC/Nylon TPU	Nylon TPU	Nylon TPU
	Basic treatment	Medium treatment	Intensive care	Intensive care
	 B01 Mattress + TKS2012 Pump	 T01 Mattress + P01 Pump	 T05 Mattress + P08A Pump	 T04 Mattress + P02 Pump
	 B01 Mattress + TKS02 Pump	 T01 Mattress + P01S Pump	 T05 Mattress + P09 Pump	 T04 Mattress + P06A Pump
	 B01 Mattress + TKS03 Pump	 T01 Mattress + TKS03 Pump	 T06 Mattress + P06A Pump	 T08 Mattress + P06B Pump
	 S01 Mattress + P01 Pump	 T01 Mattress + P05 Pump	 T06 Mattress + P09 Pump	 T09 Mattress + P06C Pump
	 S01 Mattress + P01S Pump	 T01S Mattress + TKS03 Pump	 T07 Mattress + P08B Pump	 T11 Mattress + P11 Pump
	 S01 Mattress + TKS05 Pump	 T02 Mattress + P01S Pump	 T07 Mattress + P09 Pump	 T13 Mattress + P13 Pump
	 S01 Mattress + P05 Pump	 T02 Mattress + P05 Pump	More new products are under development UPCOMING	More new products are under development UPCOMING

AIR PUMP FOR ALTERNATING PRESSURE MATTRESS SYSTEM

FUNCTIONS OVERVIEW

Air pump model	P01 Series		TKS Series		P02		P05		
									
P06 Series	P08 Series		P09		P11		P13		
									
Model Function	P01 Series	TKS Series	P02	P05	P06 Series	P08 Series	P09	P11	P13
A-B cycle alternating mode	✓	✓	✓	✓	✓	✓	✓	✓	✓
A-B-C cycle alternating mode	—	—	✓	—	✓	✓	✓	✓	✓
Static mode	✓	—	✓	✓	✓	✓	✓	✓	✓
CLP mode	—	—	—	—	—	—	—	✓	✓
Max firm mode	—	—	✓	—	—	—	—	✓	✓
Soft hardness adjustment	✓	✓	✓	✓	✓	✓	✓	✓	✓
Cycle time adjustment	—	—	✓	—	✓	✓	✓	✓	✓
Low pressure alarm	—	—	✓	—	✓	✓	—	✓	✓
Seat inflating mode	—	—	—	—	✓	—	✓	—	—
Auto turn Mode	—	—	—	—	—	—	—	✓	✓
Trun left/right or no turn key	—	—	—	—	—	—	—	✓	✓
Head up or foot up key	—	—	—	—	—	—	—	—	✓
Quiet mode	—	—	✓	—	✓	—	—	✓	✓
Panel lock	—	—	✓	—	✓	—	✓	✓	✓
CPC quick connector	—	—	✓	—	✓	✓	✓	✓	✓




ALTERNATING PRESSURE MATTRESS WITH PUMP

ISO13485 FDA CE MDR MDMA

B01+TKS2012

Multiple optional colors

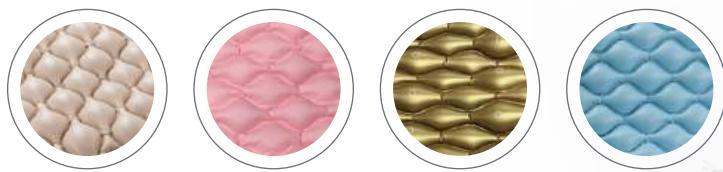

- ▼ Stage I basic treatment
- ▼ A-B lines air tubes alternating for anti bedsores/decubitus
- ▼ Individual bubble cells for maximum comfort and effectiveness
- ▼ Dual channel intake design, fast inflation speed
- ▼ Warranty: 1 year

Multiple combinations with air pumps:



MATTRESS B01

Size:	79"L*35"W*2.8"H
No. of bubbles:	130 air bubbles
End flaps:	400mm
Weight capacity:	135Kg
*Multiple color customization schemes	

Multiple combinations with air pumps:

Item	Dimension
B01 130 bubbles	79" x 35"
B02 96 bubbles	73" x 32"
B03 96 bubbles	73" x 28"
B05 130 bubbles	74" x 35"

AIR PUMP DATA


	TKS Series	P01 Series	P05
Power:	110V/60Hz or 220V/50Hz		
Air output:	≥6L/min		
Cycle time:	10~12mins		
Net Weight:	0.72~0.78Kg	1.15Kg	1.15Kg



ALTERNATING PRESSURE MATTRESS WITH PUMP

S01+P01S

ISO13485 FDA CE MDR MDMA

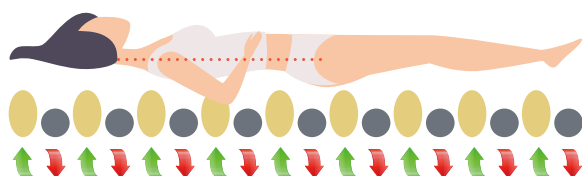

Aircells Cyclic Fluctuation, Alternating Pressure Release



Air cell A



Air cell B



Super quiet design



Biological histocompatibility



Water vapor permeation



Waterproof



Antibacterial (option)



Flame resistance (option)

- ▼ Stage I basic treatment
- ▼ A-B lines air tubes alternating for anti bedsores/decubitus
- ▼ Integrated molding
- ▼ Dual channel intake design, fast inflation speed
- ▼ Warranty: 1 year

MATTRESS S01

Size: 73"L*32"W*2.6"H after inflated

Air cell number: 22 cells

End flaps: 16"

Weight capacity: 135Kg

*Multiple color customization schemes

PUMP P01S

Power: 110V/60Hz or 220V/50Hz

Air output: ≥ 6L/min

Cycle time: 10~12 mins

Power consumption: 12VA

Fuse: 0.5A

Size: 275*120*100mm(L*W*H)

Net Weight: 1.15Kg

Multiple combinations with air pumps:



AIR PUMP DATA


	TKS Series	P01 Series	P05
Power:	110V/60Hz or 220V/50Hz		
Air output:	≥6L/min		
Cycle time:	10~12mins		
Net Weight:	0.72~0.78Kg	1.15Kg	1.15Kg

ALTERNATING PRESSURE MATTRESS WITH PUMP



T01+P08B

ISO13485 FDA CE MDR MDMA

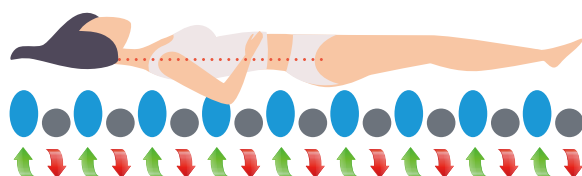

Aircells Cyclic Fluctuation, Alternating Pressure Release



Air cell A

Air cell B



Super quiet design



Biological histocompatibility



Water vapor permeation



Waterproof



Corrosion resistance (option)



Flame resistance (option)

- ▼ Stage II medium treatment
- ▼ A-B lines air tubes alternating for anti bedsore/decubitus
- ▼ CPR for quick air release in emergency
- ▼ 2-way stretch PU cover with snap fixture, fire retardant, waterproof (option)
- ▼ Warranty: 1 year

MATTRESS T01

Size: 80"L*32"~36"W*4.5"H(18cells)
5"H 16 cells and 4"H 24 cells for more options

Mattress material: Nylon PVC or Nylon TPU

Cover material: PU or thin cover

Net weight: 5.31Kg

Weight capacity: 145Kg

*Multiple color customization schemes

PUMP P08B

Power: 110V/60Hz or 220V/50Hz

Air output: 7L/min

Cycle time: 10~12 mins

Power consumption: 12VA

Fuse: 0.5A

Size: 270*130*90mm(L*W*H)

Net Weight: 1.29Kg

T01S

26 Cells



Size: 75"L*34"W*3.4"H(26cells)

Mattress material: Nylon PVC

Cover material: PU or thin cover

Bottom material: PVC tarpaulin

Net weight: 5.15Kg

Weight capacity: 145Kg

P01S



Static mode
Size: 275*120*100mm(L*W*H)
Net Weight: 1.15Kg

P05



Size: 275*100*100mm(L*W*H)
Net Weight: 1.1Kg

TKS03



Size: 215*110*95mm(L*W*H)
Net Weight: 0.73Kg

P08A




Static mode
Low pressure indicator light
Power failure alarm
Size: 270*130*95mm(L*W*H)
Net Weight: 1.29Kg

ALTERNATING PRESSURE MATTRESS WITH PUMP



T05+P09


ISO13485 FDA CE MDR MDMA

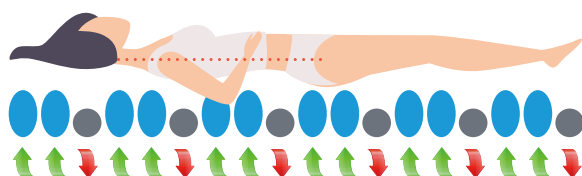

A/B/C cycle alternating mode

During the alternating process, there are always 2 sets of air cells in a fully inflated state, which means there are 12 air cells that can share and support the patient's weight.

Compared with the A/B alternating mode (where only 9 air columns remain inflated at the same time), the A/B/C cyclic alternating mode can better share the weight of the user. Users feel more comfortable and have better sleep quality.

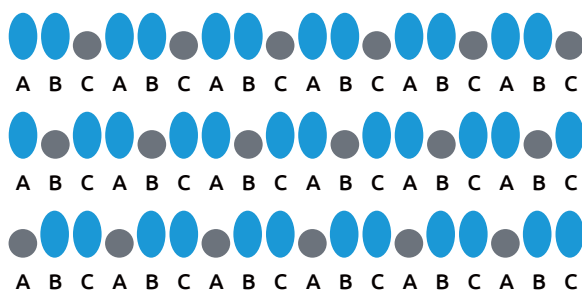

Schematic diagram of A/B/C cycle alternating mode

MATTRESS T05

Alternating mode:	A/B or A/B/C
Size:	80"L*32"~36"W*4.5"H(18cells) 4"H(24cells) and 5"H(16cells) for more options
Mattress material:	Nylon TPU
Cover material:	PU or thin cover
Net weight:	4.58Kg
Weight capacity:	145Kg

*CPR for quick air release; Pillow function.

PUMP P09

Power:	110V/60Hz or 220V/50Hz
Air output:	7L/min
Cycle time:	14mins (A-B lines alternation), 21mins (A-B-C lines alternation)
Size:	270*130*90mm(L*W*H)
Net Weight:	1.3Kg
Function:	Static/ pressure indication/ alarm/ adjustable cycle time



Nylon TPU material

TPU is more elastic and soft, providing a more comfortable sleeping experience. And it performs better in extreme temperatures, maintaining stability without becoming brittle. Tpu also offers better sealing and resistance to seepage.

P08A



Static mode
Low pressure indicator light
Power failure alarm
Size: 270*130*95mm(L*W*H)
Net Weight: 1.29Kg

More combinations can be made according to customer needs.

ALTERNATING PRESSURE MATTRESS WITH PUMP




T06+P06A

ISO13485 FDA CE MDR MDMA

In AB alternating mode, the alternating pressure mattress offers several benefits:

- 1 Pressure Relief:** Alternating inflation reduces pressure on specific body areas, preventing pressure ulcers.
- 2 Improved Circulation:** The pressure shift promotes blood flow, reducing discomfort from prolonged lying.
- 3 Enhanced Comfort:** Even pressure distribution adapts to the body, improving sleep quality.
- 4 Ideal for Long-Term Bedridden:** Provides support and comfort for patients or individuals who are bedridden for extended periods.

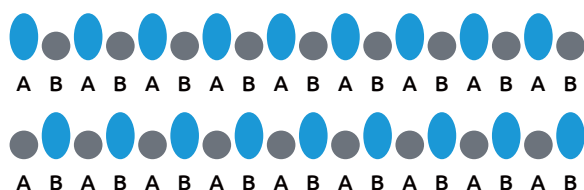

Schematic diagram of A/B cycle alternating mode



Bottom sponge pad

Adding a foam pad underneath a cyclic alternating air mattress can enhance overall comfort and support, reducing the sensation of movement caused by air flow. Additionally, the foam provides better insulation, helping to maintain warmth and preventing cold from affecting sleep.

MATTRESS T06

Alternating mode:	A/B
Size:	80"L*32"~36"W*4.5"H(18cells), 4"H(24cells) and 5"H(16 cells) for more options
Mattress material:	Nylon TPU
Cover material:	PU or Nylon PVC or pure TPU
Net weight:	6.95Kg
Weight capacity:	210Kg

*Multiple color customization schemes/ Bottom sponge pad

PUMP P06A

Power:	110V/60Hz or 220V/50Hz
Air output:	7/12 L/min
Cycle time:	6/9/12/25 mins
Size:	345*115*185mm(L*W*H)
Net Weight:	2.55Kg

*Knob adjusts the pressure

P06B



Pressure display
Size: 345*115*185mm(L*W*H)
Net Weight: 2.55Kg

P06C




Touch screen
Size: 345*115*185mm(L*W*H)
Net Weight: 2.55Kg

P09



Size: 270*130*90mm(L*W*H)
Net Weight: 1.3Kg

More combinations can be made according to customer needs.

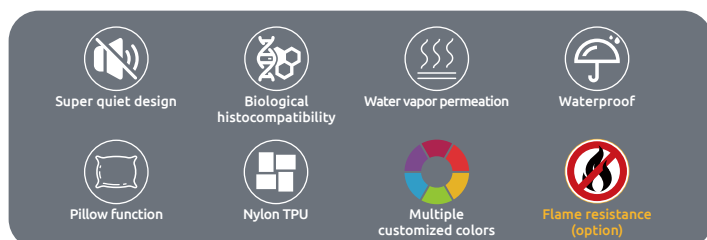


ALTERNATING PRESSURE MATTRESS WITH PUMP




T04+P02

ISO13485 FDA CE MDR MDMA

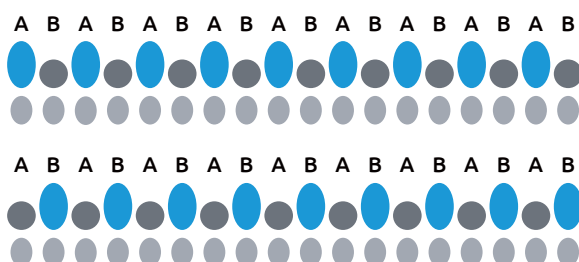

A/B cycle alternating mode

- ▼ Stage IV Intensive care
- ▼ Bottom layer can do static and upper layer A-B or A-B-C or A-B-C~S or A-B~S cycle alternating
- ▼ 2-way stretch PU cover with snap fixture, waterproof, fire retardant(option)
- ▼ CPR for quick air release in emergency
- ▼ Ventilated bladder
- ▼ Warranty: 1 year

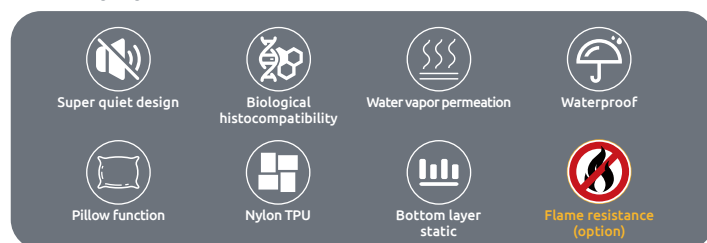

Static bottom air column provides extra support and comfort, and can enhance overall sleep quality.



The vents on a cyclic alternating air mattress are primarily used to regulate air pressure and ventilation, preventing overheating or moisture buildup.

MATTRESS T04

Size:	80"L*32"~36"W*8"H(18cells)
Bariatric size:	42"/48"(width)
Mattress material:	Nylon TPU or Pure TPU, Nylon PVC
Cover material:	PU
Net weight:	8.35Kg
Weight capacity:	210Kg
*Alternating mode: A-B or A-B-C or A-B-C~S or A-B~S; Pillow function	

PUMP P02

Power:	110V/60Hz or 220V/50Hz
Air output:	12L/min
Cycle time:	15mins
Fuse:	1A*2
Size:	258*85*250mm(L*W*H)
Net Weight:	2.5Kg
*3 or 4 air hose connection Function: 30mins static mode/ key lock/ alarm/ low pressure indicaton light	

P06A



P06B



P06C




More combinations can be made according to customer needs.

ALTERNATING PRESSURE MATTRESS WITH PUMP




T08+P06C

ISO13485 FDA CE MDR MDMA

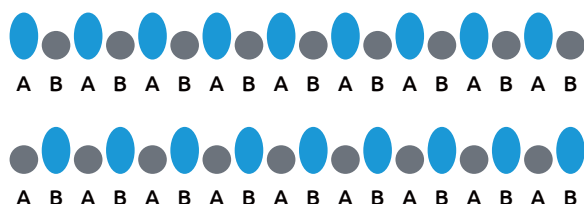

Schematic diagram of A/B cycle alternating mode

- ▼ Stage IV Intensive care
- ▼ 2-way stretch PU cover with zipper closure, waterproof, fire retardant(option)
- ▼ CPR for quick air release in emergency
- ▼ Ventilated bladder
- ▼ Warranty: 1 year

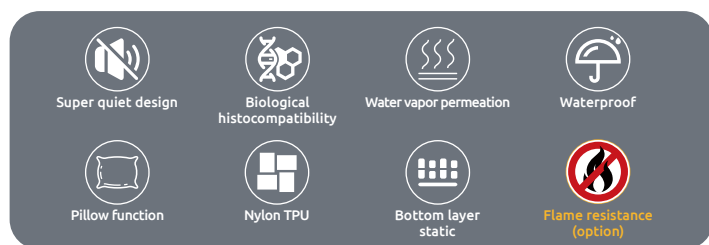

Option:
Top layer and second layer A-B or A-B-C
line alternate, bottom layer and pillow
keep static.



The vents on a cyclic alternating air
mattress are primarily used to regulate
air pressure and ventilation, preventing
overheating or moisture buildup.

MATTRESS T08

Size:	80"L*32"~36"W*9"H(20cells)
Bariatric size:	42"/48"(width)
Mattress material:	Nylon TPU or Pure TPU
Cover material:	PU
Bottom material:	Dacron PVC
Net weight:	8.35Kg
Weight capacity:	210Kg

*Three layer mattress;

Alternating mode: A-B or A-B-C or A-B-C~S or A-B~S;
Pillow Function

PUMP P06C

Power:	110V/60Hz or 220V/50Hz
Air output:	8-10 L/min
Cycle time:	6/9/12/25 mins
Size:	345*115*185mm(L*W*H)
Net Weight:	2.55Kg

*Touch screen

P02



3 or 4 air connector
Size: 258*85*250mm(L*W*H)
Net Weight: 2.5Kg

P06B



Pressure display
Size: 345*115*185mm(L*W*H)
Net Weight: 2.55Kg


More combinations can be
made according to custom-
er needs.

ALTERNATING PRESSURE MATTRESS WITH PUMP



T09+P06B

ISO13485 FDA CE MDR MDMA

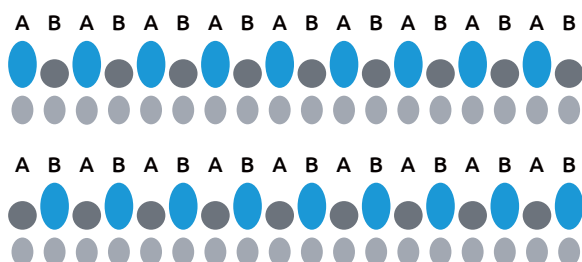

A/B cycle alternating mode

- ▼ Stage IV Intensive care
- ▼ Bottom layer can do static and upper layer A-B or A-B-C or A-B-C~S or A-B~S cycle alternating
- ▼ 2-way stretch PU cover with zipper closure, waterproof, fire retardant(option)
- ▼ CPR for quick air release in emergency
- ▼ Ankle groove design to release pressure
- ▼ Warranty: 1 year



Static bottom air column provides extra support and comfort, and can enhance overall sleep quality.



The ankle groove design reduces local pressure, promotes blood circulation, enhances comfort, reduces friction and compression in the ankle area, and is suitable for users of different body types.

MATTRESS T09

Size:	80"L*32"~36"W*8"H(18cells)
Bariatric size:	42"/48"(width)
Mattress material:	Nylon TPU or Pure TPU, Nylon PVC
Cover material:	PU
Net weight:	8.35Kg
Weight capacity:	210Kg
*Alternating mode: A-B or A-B-C or A-B-C~S or A-B~S; Ankle groove design	

PUMP P06B

Power:	110V/60Hz or 220V/50Hz
Air output:	8-10 L/min
Cycle time:	6/9/12/25 mins
Size:	345*115*185mm(L*W*H)
Net Weight:	2.55Kg
*Pressure display	

P06A



Knob adjust
Size: 345*115*185mm(L*W*H)
Net Weight: 2.55Kg

P06C



Touch screen
Size: 345*115*185mm(L*W*H)
Net Weight: 2.55Kg


More combinations can be made according to customer needs.

ALTERNATING TURNING MATTRESS WITH PUMP



T11+P11

ISO13485 FDA CE MDR MDMA

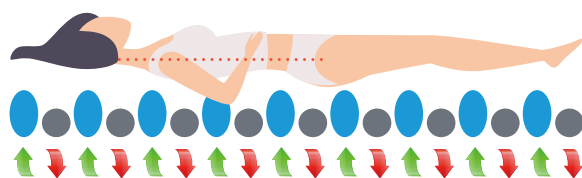

Aircells Cyclic Fluctuation, Alternating Pressure Release



Air cell A

Air cell B



Super quiet design



Biological
histocompatibility



Nylon TPU



Waterproof



Automatic turning



Water
vapor permeation





Multiple
customized colors



Flame resistance
(option)

- ▼ Stage IV Intensive care
- ▼ A-B lines air tubes alternating for anti bedsore/decubitus
- ▼ Automatic turning function (turning degree 30°)
- ▼ CPR for quick air release in emergency
- ▼ Warranty: 1 year

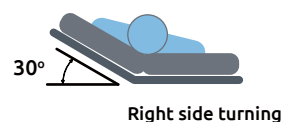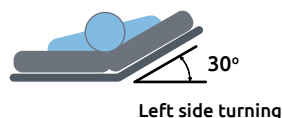

Automatic turning function

Achieve left and right side circulation decompression while lying down. It can significantly reduce the risk of pressure sores from prolonged lying down, enhance blood circulation, and improve overall comfort.

MATTRESS T11

Size:	79"L*35"W*4.5"H(18cells)
Materials:	Nylon TPU or nylon PVC
Max. User Weight:	200 kg
Cover material:	2-ways stretch PU cover, waterproof, easy clean
Working mode:	AB alternating
Function:	Auto turn function(Left/ Right 30°)
*Multiple color customization schemes	

PUMP P11

Power:	110V/60Hz or 220V/50Hz
Air output:	7 L/min
Size:	258*85*250mm(L*W*H)
Net Weight:	2.5Kg
Characteristic:	Ultra silence, noise ≤35dB
	Six grades of pressure adjustable
	Auto turn, Left/ Right turn, No turn
*4 air hose connection	



CPR straps

CPR for quick air release in emergency.

Applicable places

Widely applicable to hospital wards, nursing homes, nursing institutions, etc.

The benefits of alternating turning mattresses

Nursing staff do not need to manually flip patients, reducing workload and improving nursing efficiency.


HIGH-LEVEL+

ALTERNATING TURNING MATTRESS WITH PUMP



T13+P13

ISO13485 FDA CE MDR MDMA

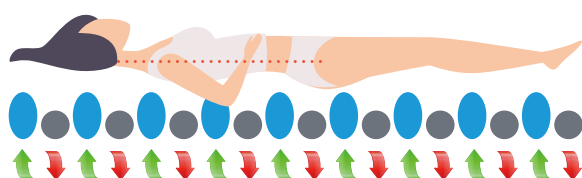

Aircells Cyclic Fluctuation, Alternating Pressure Release



Air cell A



Air cell B



Super quiet design



Biological histocompatibility



Water vapor permeation



Waterproof



Automatic turning



Head up & foot up



Multiple customized colors



Flame resistance (option)

- ▼ Stage IV Intensive care
- ▼ Head up and foot up function
- ▼ Automatic turning function
- ▼ CPR for quick air release in emergency
- ▼ Warranty: 1 year

MATTRESS T13


Size:	79"L*35"W*4.5"H(18cells)
Materials:	Top grade TPU with Nylon material
Max. User Weight:	200 Kg
Cover material:	2-ways stretch PU cover, water-proof, easy clean
Working mode:	AB Alternating
Function:	Auto Turn function(Left/ Right 30°) The function of head up and foot up

*Multiple color customization schemes; one or two layers are optional

PUMP P13

Power:	110V/60Hz or 220V/50Hz
Air output:	7 L/min
Size:	360*210*120mm(L*W*H)
Net Weight:	2.7Kg
Characteristic:	Ultra silence, noise ≤35dB Eight grades of pressure adjustable Auto turn, Left/ Right turn, No turn

*6 air hose connection

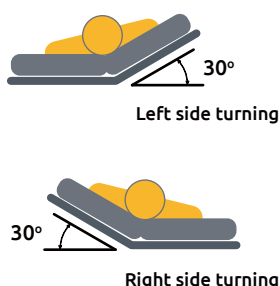

The benefits of alternating turning mattresses

Nursing staff do not need to manually flip patients, reducing workload and improving nursing efficiency.

Automatic turning function

Achieve left and right side circulation decompression while lying down. It can significantly reduce the risk of pressure sores from prolonged lying down, enhance blood circulation, and improve overall comfort.

Head up and foot up function

Improved circulation.


Pressure Relief, providing better overall comfort.

Enhanced breathing.

By improving body alignment and reducing pressure may enhance overall sleep quality.

ANTI-DECUBITUS FOAM MATTRESS

Static Pressure Ulcer
Prevention Series

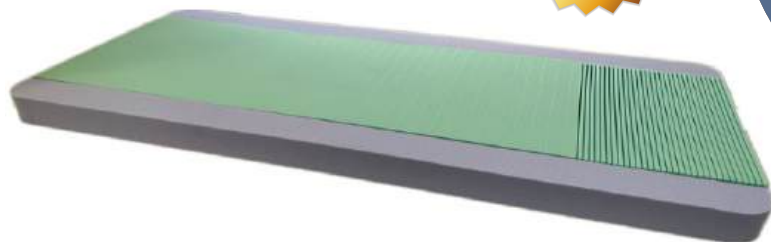

ISO13485 FDA CE MDR MDMA



Anti-decubitus



Biological histocompatibility



Water vapor permeation



Waterproof





Antibacterial (option)



Flame resistance (option)

FOAM MATTRESS

Size:	75"L*35"W*3.9"H
Materials:	Foam
Max. User Weight:	≥ 100Kg
Cover material:	Medical PU material
Warranty:	1 year

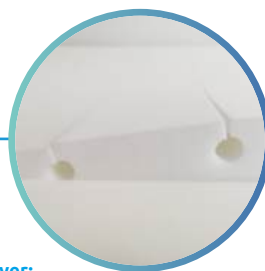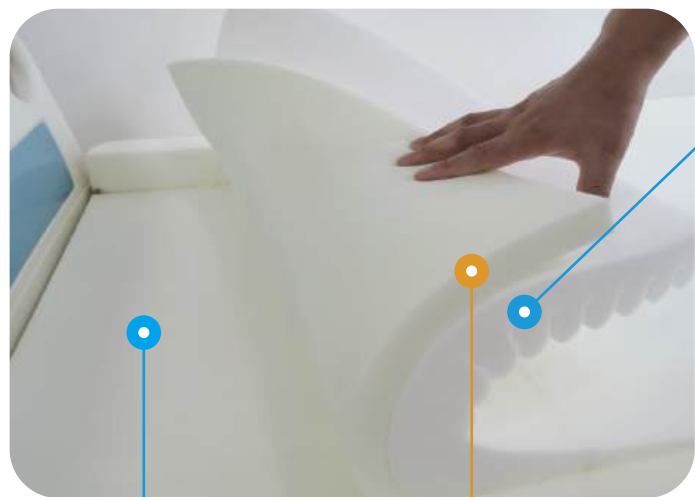

Top Layer:

High elasticity foam

The 2.5cm independent T-shaped columns from the back to the hip area effectively support different body parts of varying shapes and weights, promoting blood circulation and comfort. The leg and foot areas feature a finely cut design that provides a massage effect and adapts to the ankle.

Third Layer:
High hardness foam
Rigid support.




Second Layer:
Slow rebound memory foam
Pressure buffering layer offers weight redistribution and full pressure dispersion.



Designed with a combination of three functionally distinct foams

The upper layer features a special cut groove design for pressure relief. The middle layer is made of slow-rebounding memory foam, which maximizes pressure distribution. The lower layer consists of a high-density foam with an integrated U-shaped groove design. The combination of these three functional foams redistributes the weight, effectively dispersing pressure and preventing the formation of bedsores.

Static Pressure Ulcer Prevention Foam Mattress has right firmness. It provides appropriate and comfortable support for the human body, ensuring that the deformation of the spine will not be caused by discomfort in the softness and hardness of the mattress.

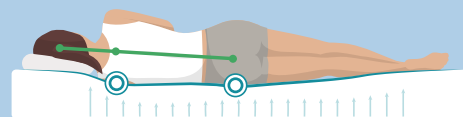

Inflatable & deflated air valves

Can quickly inflate and deflate,
easy to use.

**Hand
Pump**

STATIC BUBBLE BEDSORE PREVENTION TRANSFER MATTRESS

BS01**Pull hole**

Can conveniently and quickly
lift and transfer patients to
designated locations.



Rapid deflation

Biological
histocompatibilityWater
vapor permeation

Waterproof

ISO13485 FDA CE MDR MDMA



POWERED SELF-INFLATION TRANSFER PAD



Super quiet design

Biological
histocompatibilityWater
vapor permeation

Waterproof



Applied to first aid

Pad size (before inflation): 80"x38"

Pad size (after inflation): 75"x32"x5.9"

Material: nylon PVC or nylon TPU

Net weight: 2.5 Kg

Weight capacity: 120 Kg

Full inflation within 60 seconds.

**TM01****Automatic air pump**

Quick automatic inflation
and deflation.

Safety rope

Ensure that patients
do not fall during
the transfer process.



PORTABLE ALTERNATING PRESSURE SEAT CUSHION



C02/C03+P08M

ISO13485 FDA CE MDR MDMA



CUSHION C02

Size:	17.2"L*17.2"W*3.5"H(5cells)
Cell material:	Nylon TPU
Cover material:	PU
Weight capacity:	100Kg
Warranty:	1 year

PUMP P08M

Adapt power supply:	AC100-240Vac 50/60Hz
Fully charged power:	use 8 hours
Cycle time:	12 mins
Pump size:	275*140*95mm(L*W*H)
Pressure range(mmHg):	20 (low), 40 (medium), 60 (high)
Air output:	6~7 L/min
Static:	10, 20 mins
Net weight:	1.8Kg
Lithium battery:	22.2Vdc, 2000mAh
Warranty:	2 years

- ▼ Alternating pressure design
- ▼ Lightweight and portable
- ▼ 4-way stretch PU cover with snap fixture, fire retardant, water-proof
- ▼ Easy to clean and maintain



CUSHION C03


Size:	43"L*19.3"W*3.5"H (11cells)
Cell material:	Nylon TPU
Cover material:	PU
Weight capacity:	100Kg
*Multiple color customization schemes	



ANTI-DECUBITUS SEAT CUSHION (PVC)



ISO13485 FDA CE MDR MDMA

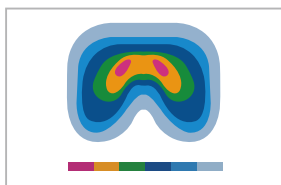

Ordinary seat cushion
Bigger sitting pressure



PVC seat cushion
Sitting area increased 30%
Pressure reduce 1/3



Medical grade materials



Exquisite workmanship and durable



Pressure dispersion



Skin friendly & breathable



Easy to clean



Excellent hip & legs support

- ▶ Pressure dispersion design reduces the risk of pressure ulcers
- ▶ Lightweight and Portable
- ▶ Users can adjust the amount of air to modify firmness and support
- ▶ PVC material is breathable, wear-resistant, and easy to clean, suitable for long-term use

Size	Height	Cell No.
400x400 mm	6 ☒ 8 ☒ 10 ☒	64 ☒ 72 ☐ 81 ☐ 90 ☐ 100 ☐
450x400 mm	6 ☒ 8 ☒ 10 ☒	64 ☐ 72 ☒ 81 ☐ 90 ☐ 100 ☐
450x450 mm	6 ☒ 8 ☒ 10 ☒	64 ☐ 72 ☐ 81 ☒ 90 ☐ 100 ☐
500x450 mm	6 ☒ 8 ☒ 10 ☒	64 ☐ 72 ☐ 81 ☐ 90 ☒ 100 ☐
500x500 mm	6 ☒ 8 ☒ 10 ☒	64 ☐ 72 ☐ 81 ☐ 90 ☐ 100 ☒



NC-802

Size: 450Lx450Wx80H mm
Bubbles: 64 pcs
Cushion material: PVC
Colors: Blue / Black
Net weight: 2.48Kg
Weight capacity: 150Kg
Warranty: 1 year
*More size can option



Pump


Cover

ANTI-DECUBITUS SEAT CUSHION (TPU)



Capacity: 100Kg

ISO13485 FDA CE MDR MDMA

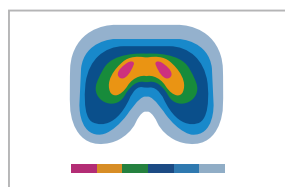

Ordinary seat cushion
Bigger sitting pressure



TPU seat cushion
Sitting area increased 30%
Pressure reduce 1/3



Medical grade materials



Exquisite workmanship and durable



Pressure dispersion



Skin friendly & breathable



Easy to clean



Excellent hip & legs support

- ▼ Pressure dispersion design reduces the risk of pressure ulcers
- ▼ Lightweight and Portable
- ▼ Users can adjust the comfort
- ▼ TPU material is environmentally friendly

Item No.	Size
DC005	410x460x50 mm (25 Bubbles)
DC006	420x420x50 mm (25 Bubbles)
DC007	460x460x50 mm (25 Bubbles)
DC008	420x420x50 mm (36 Bubbles)



DC008

DC-008

Size:	420Lx420Wx50H mm
Bubbles:	36 pcs
Cushion material:	TPU
Net weight:	0.6Kg
Weight capacity:	100Kg
Warranty:	1 year
*More size can option	



Pump



Cover

AIR INFLATABLE CUSHION



ISO13485 FDA CE MDR MDMA




460 x 430 x 33 mm

BL-1

Material: Lycra/ PU/ Suede fabric/ TPU/ Nonwovens/
Anti-slip fabric
Function: Individual inflation/ Deflation/ Adjustment
Weight: 370g/ Bearing: 200kg
Temp: -30-70°C
Usage scenario: Office/ Home/ Car

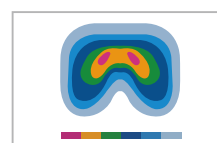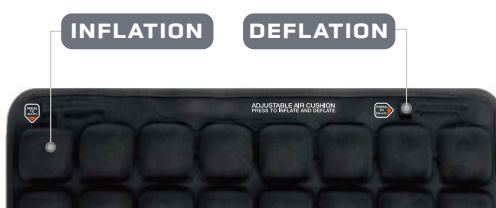

Ordinary seat cushion
Bigger sitting pressure

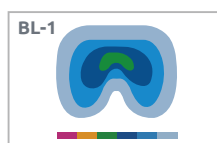

Sitting area increased 30%
Pressure reduce 1/3



400 x 450 x 50 mm

BL-2

Material: Lycra/ TPU/ Nonwovens/ Anti-slip fabric
Function: Individual inflation/ Deflation/ Adjustment
Weight: 370g/ Bearing: 200kg
Temp: -30-70°C
Usage scenario: Office/ Home/ Car

PRESS TO
INFLATE




PRESS TO
DEFLATE



Particular technique makes the cushion flat & elegant.



The bottom of covered with an anti slip layer made
of anti-static sponge.

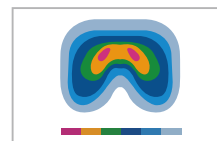

Ordinary seat cushion
Bigger sitting pressure



Sitting area increased 30%
Pressure reduce 1/3



400 x 450 x 50 mm

BL-2S


Material: Lycra/ TPU/ Nonwovens/ Anti-slip fabric
Function: Individual inflation/ Deflation/ Adjustment
Weight: 370g/ Bearing: 200kg
Temp: -30-70°C
Usage scenario: Office/ Home/ Car

PRESS TO
INFLATE
THE LEFT
AREA



PRESS TO
INFLATE
THE RIGHT
AREA


PRESS TO
DEFLATE

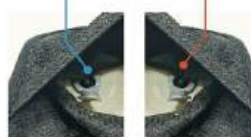

Particular technique makes the cushion flat & elegant.



The bottom of covered with an anti slip layer made
of anti-static sponge.

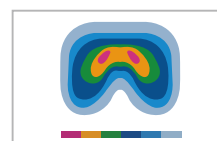

Ordinary seat cushion
Bigger sitting pressure



Sitting area increased 30%
Pressure reduce 1/3

GEL PAD FOAM CUSHION

ISO13485 FDA CE MDR MDMA



WG-ZKZD0004

SPECIFICATIONS

Material: Mesh cloth + Gel + Memory foam 45D
Size: 420L*340W*70H mm
Weight: 800g
Package: 10pcs/Ctn
Ctn: 66x46x35.5 cm

RECOMMENDED APPLICATIONS

Help relieve symptoms of hemorrhoids, prostatitis, postpartum, and lower back pain. It comes with a cooling gel layer at the top to keep cool and comfortable.



WG-ZKZD0005

SPECIFICATIONS

Material: Black flannel + Memory foam 45D
Size: 445L*360W*70H mm
Weight: 1200g
Package: 10pcs/Ctn
Ctn: 63x45.5x37 cm

RECOMMENDED APPLICATIONS

Pain and tension relieve in the lower back, coccyx, buttocks, spine and sciatic nerve areas, with a cooling gel layer on top to keep cool and comfortable.

MEMORY FOAM CUSHION



WG-ZKZD0006

SPECIFICATIONS

Material: Polyester + Memory foam 45D
Size: 445L*360W*70H mm
Weight: 800g
Package: 10pcs/Ctn
Ctn: 63x45.5x37 cm

RECOMMENDED APPLICATIONS

Relieve the pain for hip and leg caused by chronic poor sitting posture. The shape supports your tailbone, waist, spine and buttocks, reducing stress.



WG-ZKZD0008

SPECIFICATIONS

Material: Black flannel + Memory foam 45D
Size: 450L*365W*95H mm
Weight: 675g
Package: 10pcs/Ctn
Ctn: 66x46x38 cm

RECOMMENDED APPLICATIONS

Pain relief, suitable for hemorrhoids, cracks, joint pain, surgery, postpartum recovery, etc.

DVT COMPRESSION THERAPY GARMENT

ISO13485 FDA CE MDR MDMA



Senyang Tri/Uniform Pulse Sequential Garments

Application	Item No.	Size
Foot	DVT-F10	Uniform
Calf	DVT-C20	S/M/L
Boot	DVT-B30	S/M/L
Thigh	DVT-T40	S/M/L



All-round
massage



Air pressure
massage



Breathable
fabric



Time
controllable



Suitable for various types
of air pumps



Convenient &
wearable

DVT sleeves are used to relieve leg pain, promote muscle recovery, and treat shinitis and varicose veins.

DVT Manufacturer Over 20 years / OEM & ODM

Why use DVT sleeves?



Thrombosis

Damaged open and warped valves
Blood flowing in both direction
Blood accumulating in vein

Embolus

When blood moving slowly, a clot or thrombus can form.
A part of the clot can break off and travel in the
bloodstream (embolus).



Use DVT compression therapy garment
for treatment



Normal vein

Blood moving toward
the heart

DVT compression therapy garments play a crucial role in the prevention and treatment of DVT by improving blood circulation, reducing blood stasis, and alleviating swelling and discomfort. For high-risk individuals or patients who have already developed blood clots, compression therapy serves as an effective adjunctive treatment that can significantly reduce the risk of disease progression.



Calf sleeve



Foot sleeve



Calf thigh sleeve

AIR PRESSURE THERAPY SYSTEM



PT1002

ISO13485 FDA CE MDR MDMA



PUMP PT1002

Massage mode: A.B mode available

Pressure range: 20~200mmHg($\pm 10\sim 15\%$)

Timer: 10.20.30 mins

Power consumption: 40VA

Voltage: AC110V/60Hz or AC220V/50Hz

Machine size: 260*170*130mm (L*W*H)

Machine weight: 2.5Kg

Certificate: FDA / 510K(K161907) / CE

*Multiple color customization solutions for air pumps and sleeves.



The Air Pressure Therapy System is a medical device designed to improve blood and lymph circulation, pain edema relief, prevent DVT through intermittent pneumatic compression.

The key feature: Multi-setting, easy application.



Air pressure
massage



All-round
massage



Two massage
mode



Time
controllable



Breathable
fabric



Loop mode

AIR PRESSURE THERAPY SYSTEM



PT1003

ISO13485 FDA CE MDR MDMA

PUMP PT1003

- Massage mode:** A.B mode available
- Pressure range:** 20~200mmHg($\pm 10\sim 15\%$)
- Timer:** 10.20.30 mins
- Power consumption:** 40VA
- Voltage:** AC110V/60Hz or AC220V/50Hz
- Machine size:** 280*175*130mm (L*W*H)
- Machine weight:** 2.5Kg
- Certificate:** FDA / 510K (K181409-prescription) / 510K (K192358-OTC)
- *Multiple color customization solutions for air pumps and sleeves.



After sports or other high-intensive exercises, the pressure therapy system can relief muscle sore and getting body recovery soon.

The key features: Multi-setting, easy carry with handle.



Air pressure
massage



All-round
massage



Two massage
mode



Time
controllable



Breathable
fabric



Loop mode

AIR PRESSURE THERAPY SYSTEM



ISO13485 FDA CE MDR MDMA



PT1005



PUMP PT1005

Massage mode: 6 mode

Pressure range: 20~200mmHg

Timer: Max 30mins

Power adaptor input:

100-240Vac 50/60Hz 0.8A max

Machine size: 285*90*90mm (L*W*H)

Machine weight: 1.6Kg

*Multiple color customization solutions for air pumps and sleeves.



After completing yoga, running or other high-intensity exercises, the pressure therapy system stimulates blood and lymphatic circulation in the arms, legs and other parts, helping the body quickly recover from fatigue.



Air pressure
massage



All-round
massage



Two massage
mode



Time
controllable



Breathable
fabric



Loop mode

AIR PRESSURE THERAPY SYSTEM



PT1006

ISO13485 FDA CE MDR MDMA

6 usage modes that can be freely switched to easily meet different needs.

Mode 1	
Mode 2	
Mode 3	
Mode 4	
Mode 5	
Mode 6	Mode 1 + Mode 2 + Mode 3 + Mode 4 + Mode 5



Size: 206x100x70mm
Pressure range: 20-200mmHg
Massage mode: 6 Mode
Massage cycle time: Max 30min
Power Adaptor input:
100-240Vac 50/60Hz 0.8A
Power Adaptor output: 12Vdc 4.0A



The compact size and sleek appearance embody excellent industrial design language, and the human-machine interaction experience is integrated into ergonomic design. Lightweight and durable, easy to carry during travel.

2 battery capacities are available for selection: 7.4V/5000mAh and 7.4V/10000mAh.



Air pressure
massage



All-round
massage



Six massage
mode



Time
controllable



Breathable
fabric



High capacity
lithium battery



Cuff Name	S	M	L	XL	XXL	XXXL
Vest Cuff	/	/	/	/	/	/
Arm Cuff	700*250mm	/	900*250mm	/	/	/
Waist Cuff	1200*400mm	/	1450*400mm	/	/	/
Hip Cuff	500*600mm	/	600*600mm	/	/	/
Leg Cuff	900*300mm	/	1100*300mm	1100*350mm [double zippers]	1100*400mm [double zippers]	1100*500mm [double zippers]
Calf Cuff	/	/	/	/	/	/
Leg Sleeve Extension Piece: 1100x100mm				Cuff color can be customized		

WHY US?

Providing you with high-quality physical therapy and relaxation experience

Features:

- *Over lap
- *Promote wound healing
- *Promote wound healing
- *User friendly home care
- *Preventing DVT
- *Increase lymphatic and vascular flow
- *Provide convenient sports injury & fatigue recovery
- *Easy to use, improving patient compliance
- *Compact and foldable design, easy to carry



BMW01

POWERED & MANUAL NURSING TRANSFER LIFT



ISO13485 FDA CE MDR MDMA

BMW01 Powered Nursing Transfer Lift

Maximum outline size:

545*725*980mm (L*W*H)

Lifting range: 250mm

Machine weight: 23Kg

Charging requirements:

Dedicated adapter charging

Rated power: 100W/2A

Load weight: ≤150Kg

Electric lift is mainly used as nursing equipment to help disable or elder people move without barriers, used for short-distance displacement and rehabilitation care. Through the control handle, nursing staff can easily transfer the patients. When using, they can push forward or backward or lock the rear wheels, so that the elderly or the disabled can be barrier-free in hospital beds and toilets, the living room, outdoors and other places. It greatly reduces the work intensity of nursing staff and family members, improves nursing efficiency, reduces nursing risks and prevents patients from encountering secondary injuries during transfer.



Painted
steel pipe



Free lift



Strong
bearing capacity



Easy opening
and closing



Easy
installation



Silent
universal wheel



POWERED & MANUAL NURSING TRANSFER LIFT



ISO13485 FDA CE MDR MDMA



BMW02



BMW02 Powered Nursing Transfer Lift

Maximum outline size:

545*725*980mm (L*W*H)

Lifting range: 250mm

Machine weight: 28Kg

Charging requirements:

Dedicated adapter charging

Rated power: 100W/2A

Load weight: ≤150Kg

Features

- I Enhanced Safety:** It reduces the risk of falls and injuries during transfers for both patients and caregivers.
- II Improved Patient Comfort:** The device typically provides a smoother and more comfortable transfer experience, reducing anxiety for patients.
- III Versatile Use:** It can be used in various settings, including hospitals, nursing homes, and at home, offering flexible transfer solutions.
- IV User-Friendly Operation:** Many models feature simple controls, making them easy for caregivers to operate with minimal training.
- V Promotes Independence:** In some cases, patients can use the lift with guidance, enhancing their independence and confidence.



Painted steel pipe



Free lift



Strong bearing capacity



Easy opening and closing



Easy installation



Silent universal wheel



BMW03

POWERED & MANUAL NURSING TRANSFER LIFT



ISO13485 FDA CE MDR MDMA



BMW03 Powered Nursing Transfer Lift

Maximum outline size:

545*725*980mm / 480*280*710mm (L*W*H)

Lifting range: 250+100mm

Machine weight: 33Kg

Charging requirements:

Dedicated adapter charging

Rated power: 100W/2A

Load weight: ≤150Kg

Features

- I Enhanced Safety:** It reduces the risk of falls and injuries during transfers for both patients and caregivers.
- II Improved Patient Comfort:** The device typically provides a smoother and more comfortable transfer experience, reducing anxiety for patients.
- III Versatile Use:** It can be used in various settings, including hospitals, nursing homes, and at home, offering flexible transfer solutions.
- IV User-Friendly Operation:** Many models feature simple controls, making them easy for caregivers to operate with minimal training.
- V Body suspension lifting system:** You can suspend your body onto the displacement machine to compensate for situations beyond the suitable height of the lifting platform.



Painted
steel pipe



Free lift



Strong
bearing capacity




Easy opening
and closing



Easy
installation



Silent
universal wheel

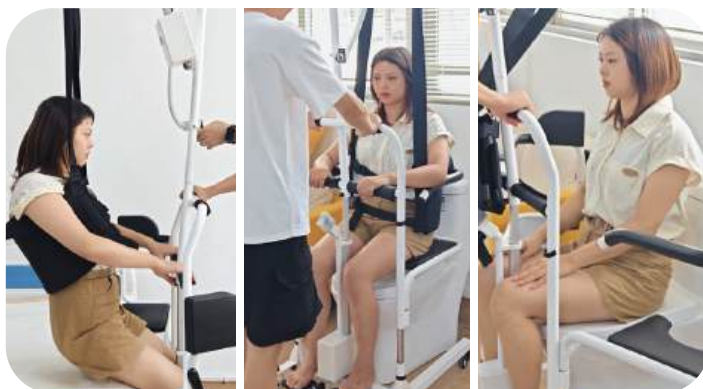

POWERED & MANUAL NURSING TRANSFER LIFT



ISO13485 FDA CE MDR MDMA



BMW05



BMW05 Powered Nursing Transfer Lift

Maximum outline size:

880*610*320mm (L*W*H)

Lifting range: 300mm

Machine weight: 31.5Kg

Charging requirements:

Dedicated adapter charging

Rated power: 72W/1.5A

Load weight: ≤150Kg

Features

- I Enhanced Safety:** It reduces the risk of falls and injuries during transfers for both patients and caregivers.
- II Improved Patient Comfort:** The device typically provides a smoother and more comfortable transfer experience, reducing anxiety for patients.
- III Versatile Use:** It can be used in various settings, including hospitals, nursing homes, and at home, offering flexible transfer solutions.
- IV User-Friendly Operation:** Many models feature simple controls, making them easy for caregivers to operate with minimal training.
- V Remote control:** This function can improve operational flexibility and safety, reduce the physical burden on operators, and allow for precise control from a greater distance, avoiding potential hazards.



Painted steel pipe



Free lift



Strong bearing capacity



Easy opening and closing



Remote control



Easy installation



Silent universal wheel



Mobile toilet

POWERED & MANUAL NURSING TRANSFER LIFT



BMW06

ISO13485 FDA CE MDR MDMA

BMW06 Powered Nursing Transfer Lift

Maximum outline size:

880*610*320mm (L*W*H)

Lifting range: 300mm

Machine weight: 31.5Kg

Charging requirements:

Special adapter charging

Rated power: 72W/1.5A

Load weight: ≤150Kg

Features

- I Enhanced Safety:** It reduces the risk of falls and injuries during transfers for both patients and caregivers.
- II Improved Patient Comfort:** The device typically provides a smoother and more comfortable transfer experience, reducing anxiety for patients.
- III Versatile Use:** It can be used in various settings, including hospitals, nursing homes, and at home, offering flexible transfer solutions.
- IV User-Friendly Operation:** Many models feature simple controls, making them easy for caregivers to operate with minimal training.
- V Remote control:** This function can improve operational flexibility and safety, reduce the physical burden on operators, and allow for precise control from a greater distance, avoiding potential hazards.



Painted
steel pipe



Free lift



Strong
bearing capacity



Easy opening
and closing



Remote
control



Easy
installation



Silent
universal wheel



Mobile
toilet



POWERED & MANUAL NURSING TRANSFER LIFT



ISO13485 FDA CE MDR MDMA



BMW07



Folding state

BMW07

Powered Nursing Transfer Lift

Maximum outline size:

840*595*960mm (L*W*H)

Lifting range:

875~1220mm

Machine weight:

32Kg

Base inside width:

495mm

Base max expandable:

870mm

Folding size:

600*390*920mm (L*W*H)

Load weight:

≤150Kg

Features

- I Enhanced Safety:** It reduces the risk of falls and injuries during transfers for both patients and caregivers.
- II Improved Patient Comfort:** The device typically provides a smoother and more comfortable transfer experience, reducing anxiety for patients.
- III Versatile Use:** It can be used in various settings, including hospitals, nursing homes, and at home, offering flexible transfer solutions.
- IV Folding design:** More convenient for storage and carrying, suitable for long-distance travel.
- V Adjustable legs:** Enhanced adaptability and compatibility to the environment, resulting in a better user experience.



Painted steel pipe



Free lift



Strong bearing capacity



Folding design



Adjustable legs



Easy installation



Silent universal wheel



Mobile toilet



BMW08

POWERED & MANUAL NURSING TRANSFER LIFT



ISO13485 FDA CE MDR MDMA

BMW08 Hydraulic Lift Transfer Chair

Dimension: 518*725*980mm (L*W*H)

Lifting range: 200mm

Machine weight: 30Kg

Load weight: ≤150Kg

Features

- I Enhanced Safety:** It reduces the risk of falls and injuries during transfers for both patients and caregivers.
- II Improved Patient Comfort:** The device typically provides a smoother and more comfortable transfer experience, reducing anxiety for patients.
- III Versatile Use:** It can be used in various settings, including hospitals, nursing homes, and at home, offering flexible transfer solutions.
- IV User-Friendly Operation:** Many models feature simple controls, making them easy for caregivers to operate with minimal training.
- V Hydraulic lifting system:** Adopting hydraulic lifting, the structure is simple, the operation is convenient and fast, and it does not rely on external power sources.



Painted
steel pipe



Free lift



Strong
bearing capacity



Easy opening
and closing



Hydraulic
lifting



Easy
installation



Silent
universal wheel



Mobile
toilet



POWERED NURSING TRANSFER LIFT



ZNHLY001

ISO13485 FDA CE MDR MDMA



ZNHLY001 Powered Nursing Transfer Lift

Dimensions: 891x676x1166mm (L*W*H)

Battery: Medical lithium iron phosphate

Waterproof rating: IP22

Load capacity: 150Kg

Charging: Special adapter

Driving speed: 2KM/H

Control: Joystick+buttons



IP22 waterproof
can shower



Folding design
compact storage



Joystick
operation



Patient
autonomous operation



Silent
operation



Anti impact &
anti self ignition



150 kg
load-bearing capacity



2 km
per hour



2 PANEL CONTROL OR MANUAL OPERATION MODES OF OPERATION

Suspended Seat System

Load weight: $\geq 150\text{Kg}$,
instantaneous load-bearing capacity $\geq 270\text{Kg}$.
Foam cushion, safe and comfortable.

Multi-function control panel

Easy and convenient operation, and
patients can operate on their own.

Fully enclosed joint design

Anti collision, good reliability,
no pollution.



Leg widening mechanism

The leg spacing can be widened to 76 cm
to adapt to a wider range of places.

Legs extension
range 52-76cm

IP22 waterproof

Can be used in the shower,
fully aluminum body.
Lithium iron phosphate batteries
for automobiles, anti impact and
self ignition.

Silent tires

Use low noise.
The speed is 2 kilometers
per hour and can be set.

Folding design

Space saving, convenient
storage, suitable for long-distance
transportation.



MEDICAL AUXILIARY PRODUCT

ISO13485 FDA CE MDR MDMA



MPJ-01

Transfer Pad

Pad size: 1050x720 mm

Material: Nylon, 6mm PE cotton, imitation nylon ribbon, anti-skid leather

- Easy to fold and store



MPJ-02

Auxiliary Tension Band

Material: Nylon+3mm PE cotton

Expand Size: 1200x180x5mm

Folded size: 310x180x20mm

Load-bearing: 100 Kg

Weight: 0.285 Kg

- Waterproof nylon fabric, easy to clean



MPJ-03

Transfer Waist Belt Sling For Stand Assist

Size: 900x120x6 mm

Material: 600D polyester cloth, 4mm EVA, 3mm sponge

Loading: 225 Kg

- Ergonomic compliance
- High Weight Capacity
- Well-made and durable



MPJ-04

Auxiliary Turning Pillow

Material: EVA foam+PU leather

Size: 460x400x90 mm

Weight: 0.65 Kg



MPJ-05

Triangle Pad R Turning Pillow

Pad size: 600x250x200 mm

Material: EVA, cotton cloth



A Leading Manufacturer Of Home Medical Equipment

www.medicalxm.com

ISO13485 FDA CE MDR MDMA

 vip@medicalxm.com

 No.189-1, Tongji South Road, Jimei District, Xiamen,
361000, Fujian, P.R. China.

 00-86-0592-5026095

